

# COMPLEX

# APPROACH *with*



Deta  
Elis  
Holding  
Health  
Innovation  
Human





## COMPLEX "WITHOUT COLD"

The action of the complex is aimed at enhancing the action of the immune system, increasing resistance and reducing recovery time.



## COMPLEX "STOP STRESS"

The complex is aimed at harmonizing the nervous system, relieving chronic fatigue and improving overall well-being.



## COMPLEX "STOP INSOMNIA"

The action of the complex is aimed at improving resistance to mental and psychological stress, reducing nervous excitement.



## COMPLEX "WOMEN'S HEALTH"

The action of the complex is aimed at restoring the functions of the excretory system and reducing psychoemotional stress.



## COMPLEX "MOVEMENT WITHOUT PAIN"

The action of the complex is aimed at relieving joint inflammation, recovering from injuries and improving the functioning of the musculoskeletal system.



## COMPLEX "SPORT PERFORMERS"

The action of the complex is aimed at increasing strength indicators, normalizing muscle tone and relieving pain after training.



## COMPLEX "BODY FIGURE"

The action of the complex is aimed at normalizing metabolic processes, reducing body weight, improving the psychoemotional state and quality of life.



## COMPLEX "STOP ALLERGY"

The action of the complex is aimed at eliminating biological allergens, increasing immunity, cleansing of toxins.



## COMPLEX "HEART AND CIRCULATION PACK"

The action of the complex is aimed at improving blood microcirculation and harmonizing heart rate.



## COMPLEX "DIGESTION WITHOUT PROBLEMS"

The action of the complex is aimed at improving the secretory and motor functions of the digestive system and restoring the functions of the excretory system.



## COMPLEX "MEN'S HEALTH"

The action of the complex is aimed at protecting against infections, regulating sexual activity and restoring physiological parameters.



## COMPLEX "CLEAN SKIN PACK"

The action of the complex is aimed at getting rid of peeling and inflammation, increasing the immunity of the skin and improving the external condition of the skin, hair and nails.



## COMPLEX "STOP PARASITES"

The action of the complex is aimed at detoxifying the body, increasing immunity, harmonizing the work of internal organs and recovering from cleansing from parasites.



## COMPLEX "HEALTHY CELL"

The action of the complex is aimed at reducing the negative effects of the environment and toxic load, preventing the development of undesirable consequences for the body.



## COMPLEX "HEALTHY TEETH"

The action of the complex is aimed at reducing bleeding of the gums, improving the health of the mucous membranes, strengthening local immunity, and early recovery from dental interventions.



## COMPLEX "LYME DISEASE"

The action of the complex is aimed at reducing parasitic and toxic loads, maintaining harmonious functioning of the nervous system, increasing the energy state and improving overall well-being.



## COMPLEX "STOP CORONAVIRUS"

The action of the complex is aimed at enhancing the defences of the organism, increasing resistance, reducing the likelihood of complications and shortening the recovery time.



## COMPLEX "STOP INFECTION"

The action of the programs aims to destroy pathogenic viruses and bacteria, which are the causative agents of the most common infections, as well as to remove toxins from the body.



## COMPLEX "HEALTHY BACK"

The action of the programs is aimed at improving the condition of the spine and back muscles, reducing pain, inflammation and swelling. Facilitates freedom of movement.



## COMPLEX "GASTRITIS+"

The action of the programs is aimed at improving the functioning of the gastrointestinal organs: reduction of inflammatory reactions, restoration of mucous membranes, elimination of unpleasant sensations. It also contributes to the normal functioning of the stomach, pancreas, liver and all parts of the intestine, including the irritable one.



## COMPLEX "HEALTHY SKIN"

The action of the programs is aimed at fighting pathogenic bacteria and fungi of the skin, reducing inflammatory reactions, increasing skin immunity, improving the appearance of the skin, hair and nails, soothing the symptoms of psoriasis and acne.



## COMPLEX "FIRST AID"

The action of the programs is aimed at the fastest elimination of unpleasant sensations: nausea, different types of pain, inflammation, tendency to vomit.



## COMPLEX "HEALTH 60+"

The action of the programs aims to support the health of the elderly, taking into account the main problems with sugar, blood microcirculation, injuries, pressure, digestive, inflammation, dementia, parkinson's etc. that characterize this category of the population. It contributes to the achievement of active longevity.



## COMPLEX "HEALTHY CHILD"

The programs aim to maintain the health of children from 0 to 12 years, taking into account their age characteristics. They contribute to the strengthening of the immune system, better adaptation to the new conditions in kindergarten and school. They protect against childhood diseases and participate in quick recovery. They significantly enhance the improvement of symptoms in children who are on the autism spectrum.

# CONTENTS

## BASIC COMPLEXES

Digestion without problems.....	7
Stop parasites.....	11

## THEMED COMPLEXES

Heart and circulation.....	17
Stop stress.....	21
Stop insomnia.....	25
Women's health.....	29
Men's health.....	33
Movement without pain.....	37
Clean skin.....	41
Sport performers.....	45
Stop allergy.....	49

## PARTICULAR COMPLEXES

Active immune system.....	55
Stop Coronavirus.....	59
Stop Infection.....	63

## SPECIAL COMPLEXES

(Used only after one of the basic complexes)

Body figure.....	69
Healthy Cell.....	73
Lyme Disease.....	77
Healthy teeth.....	81
Gastritis+.....	85
Healthy Back.....	89
Healthy Skin.....	93
First Aid.....	97
Health 60+.....	98
Healthy Child.....	99

## PROGRAMS

The list of manual programs DeVita AP.....	100
Description of automatic programs DeVita AP.....	118
The list of manual programs DeVita Ritm.....	122
Description of automatic programs DeVita Ritm.....	131

# CONTENTS

## REFERENCE INFORMATION

Explanations for the use of programs DeVita Ritm.....	134
Explanations for the use of programs DeVita AP.....	141
The biological role of DeLixir.....	156
Recommended DeVita programs for use with various functional problems.....	158
DeMatrix Applications.....	165

## INTRODUCTION

The complex approach to health by the DEHolding Company is based on the homotoxicology theory by Dr. H.H. Reckeweg. Dr. Reckeweg introduced the concept of "homotoxin", meaning any substance toxic to humans, regardless of their nature. They can be endotoxins (metabolic products or unsplit/not removed in sufficient quantities) or exotoxins (harmful environmental factors, pathogens of infectious diseases, parasites, etc.). Due to the accumulation of such substances in the human organism, certain diseases can appear. Homotoxicology considers health as "A state of freedom of an organism from homotoxins and the absence of functional or organic lesions caused by them".

The complex approach makes it possible to gently, without stress for the organism, normalize the internal environment of the body, increase energy, and improve the quality of life. To achieve this goal, before using any of the special complexes, it is necessary to take a course of a gastrointestinal tract restoration using the basic complexes "Digestion without problems" or "Stop parasites". They are fundamental for releasing the organism from toxic load and normalizing the activity of all excretory organs.

# BASIC COMPLEXES





## COMPLEX DIGESTION WITHOUT PROBLEMS

### WHY IS IT IMPORTANT MAINTAIN HEALTH OF THE GASTROINTESTINAL TRACT?



The gastrointestinal tract is a complex system on which the vital activity of the whole organism depends. In this system the food is processed into substances necessary for the functioning of organs. Also, the processes of obtaining energy and building material for new cells begin in the gastrointestinal tract.



More than half of the world's population suffers from digestive disorders. Their main causes are an unbalanced diet, poor ecology, poor eating habits (overeating, diets, nutrition on the go, convenience foods, etc.). One in four people has a history of a disease of the gastrointestinal tract.



This leads to serious consequences - from problems with body weight to diseases that require surgery. In addition, problems with the gastrointestinal tract also affect overall well-being, performance and appearance. Improper metabolism causes skin problems - acne, dryness or greasiness, dull complexion, wrinkles.



Common digestive problems are pain after eating and on an empty stomach, heartburn, stool disorders, nausea and flatulence. A person is forced to follow a special diet. Discomfort associated with the stomach and intestines affect the lifestyle, making it less active, and puts a limit on travelling.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites disrupt the secretory and motor functions of the stomach, pancreas, liver and intestines. Parasites can block the ducts of the liver and pancreas with their bodies, causing their inflammation. Helminths are also conductors for other microbes. DeVita AP basic antiparasitic programs affect most known pathogens, quickly, effectively and without side effects, removing the parasitic load from the organism.

The special program “No schistosomes (blood-flukes) & trematods” relieves from the parasitic burden of the liver and blood, restoring their functions and reducing intoxication.

Special programs “Healthy Gastrointestinal”, “Help against poisoning”, “Drainage Antitox” improve the condition of intoxication or intestinal flu, neutralize and remove toxins.

“Drainage -General” program is used to cleanse the body after applying basic antiparasitic programs, reduces parasitic intoxication and improves overall well-being.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy day” and “Health” automatic programs normalize the work of the stomach and intestines, which is especially important for people with a violation of the diet, helps to increase the body’s defenses and harmonize the work of all systems.

Special programs “Digestion normalization” and “Healthy pancreas” are aimed at improving the functioning of the pancreas, normalizing the functions of the intestines, which is important not only to restore overall digestion, but also to improve absorption, reduce severity and bloating, restore immunity, metabolism normalization.



**DeVita Energy** restores physiological processes in the digestive system, starts the mechanism of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- pain, bloating;
- heartburn, a feeling of heaviness;
- irregular bowel movements.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements.



# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Energy	Energy or Energy Pro	
2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Energy	Energy or Energy Pro	
2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
5		Special programs	No schistosomes (blood-flukes) & trematods	Daily in the evening after DeVita Ritm programs for 1 <sup>st</sup> -2 <sup>nd</sup> weeks.
6			Drainage Antitox	Daily in the evening after DeVita Ritm programs for 3 <sup>rd</sup> -4 <sup>th</sup> weeks.

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Energy	Energy or Energy Pro	
2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Energy	Energy or Energy Pro	
2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
3		Special programs	Digestion normalization and / or Healthy pancreas	Daily in the afternoon, if needed.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites-Automatic	At night, alternating every other day for 2 weeks

IN CASE OF INTOXICATION	DEVICE	PROGRAM		USAGE
	DeVita Ritm	Special programs	No poisoning - automatic	
DeVita AP	Special programs	Healthy Gastrointestinal, Help against poisoning, Drainage Antitox		Daily for 10 days.

# PREVENTIVE MEASURES FOR DISEASES OF THE GASTROINTESTINAL TRACT



Drink fluids 15-30 minutes before or 2-4 hours after eating, but not with meals. Fluids dilute the gastric acid and interfere with digestion.



Eat fruits separately from other foods: either 30 minutes before or 3-4 hours after eating.



Obtain moderation in eating. Get up from the table with a feeling of slight hunger. The signal of the stomach about saturation does not immediately go to the brain. By eating up to satiety, you risk overeating.



Eat freshly cooked foods. It is advisable to cook food immediately before eating it.



Eat small portions 4-5 times a day. You should not experience severe hunger in between meals, but you should not eat if you don't feel hungry.



When eating, do not get distracted by reading, TV, long conversations. Also, do not eat in a bad mood. All this interferes with the proper digestion of food.



Eat slowly, chew for a long time and carefully so that the food is saturated with saliva as much as possible.



In the evening, it is advisable not to eat 2 hours before bedtime. A feeling of hunger can be quenched with something light - for example, low-fat cottage cheese or kefir.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Pain, discomfort in the abdomen, heartburn, heaviness, bloating.

After the first use  
of DeVita



Reduced pain, relieve from heartburn and heaviness, improved mood.

After a month of using DeVita



Normalization of digestion, elimination of exacerbation of chronic diseases, excellent health.



Use the scheme of the "DIGESTION WITHOUT PROBLEMS" complex together with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX STOP PARASITES

### WHY PARASITE INFESTATIONS ARE DANGEROUS?



Parasites enter the body through skin, food, water, dust, from pets or other people. They can settle in all organs and tissues, and also live under the skin or on the skin. Parasites absorb vital elements, secrete toxins, destroy blood cells.



Parasites trigger many chronic diseases of the digestive tract, bronchopulmonary system, and skin. Any health-improving measures (proper nutrition, physical exercises, hardening procedures) without ridding the body of parasites will not give a positive effect.



Parasitic diseases are dangerous because they are perfectly disguised as other diseases. Symptoms of parasite infection include chronic fatigue, irritability and anxiety, hyperactivity in children, anemia, brittle nails and hair, headaches, appetite disorders, decreased immunity.



Treatment for parasites, especially in advanced cases, is very aggressive towards the body. Antiparasitic drugs have a negative effect on the liver and gastrointestinal tract, killing beneficial microflora. Therefore, it is much easier to prevent the entry and multiplication of parasites in the body.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Bacteria, viruses, fungi, parasites or their toxins are the main cause of poor health and one of the leading causes of death. Parasites weaken the immune system, increase the degree of intoxication, affect organs, change the pH of the internal environment, and create an imbalance in the body. The basic anti-parasitic programs on the DeVita AP device affect the majority of known pathogens, quickly, efficiently and without side effects, removing the parasitic load and signs of inflammation.

“Cleansing” automatic program is used to cleanse the body after applying basic antiparasitic programs, reduces parasitic intoxication and improves overall well-being.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” and “Health” automatic program improve blood circulation in organs, relieve from stress, which helps to increase immunity.

“Complete Cleansing” automatic program affects organs taking into account their daily activity, gently cleansing the body of accumulated toxins and harmonizing the work of all organs and systems.



**DeVita Energy** activates self-regulation mechanisms, helps in the acceleration of the recovery processes and the distribution of energy balance. Allows the body to quickly get in shape after cleansing from parasites, maintains vigor and good mood.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- inflammatory processes, the appearance of allergies, the deterioration of the skin, hair and nails;
- a change in appetite, increased fatigue, impaired bowel movements;
- in the presence of unexplained pain;
- to maintain health in high-risk groups: pet owners, travelers.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Intestine without problems	After "Healthy Day" program.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Intestine without problems	Daily in the afternoon.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
	5		Special programs	No bacteria	Daily in the evening after DeVita Ritm programs, 14 days.
6	No herpes simplex, No herpes, No hepatitis A			If necessary, according to guidelines.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Intestine without problems	After "Healthy Day" program.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites - Automatic	At night, alternating every other day for 2 weeks.	

# PREVENTIVE MEASURES FOR PARASITIC INFECTIONS



Thoroughly heat food and water before use. Boiling is preferable to frying. Try not to eat meat of wild animals, river fish, as well as dried and slightly salted fish.



Before you travel to tropical countries, take the necessary vaccinations. Be wary of local food.



Try not to walk barefoot on the beach - nematodes live in the sand.



After contact with stray animals, wash your hands with alkaline soap (household) - a soft toilet soap may not completely remove the parasite larvae.



Wear gloves when working in the garden.



If you travel often, walk in nature, work in the garden, have contact with pets, carry out drug prophylaxis of the parasites once a year.



Wash fruits and vegetables with soap and brush for at least 1 minute.



If you have problems with the gastrointestinal tract, muscle pain, general malaise, consult a doctor.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Headaches, nausea, lingering cough, feeling tired, weakness, impaired bowel movements.

After the first use  
of DeVita



Improved general condition of the organism.

After a month of  
using DeVita




Energy boost, disappearance of pain, normalization of all systems.



Use the scheme of the "STOP PARASITES" complex together with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements, and your results will become noticeable much faster.



# THEMED COMPLEXES

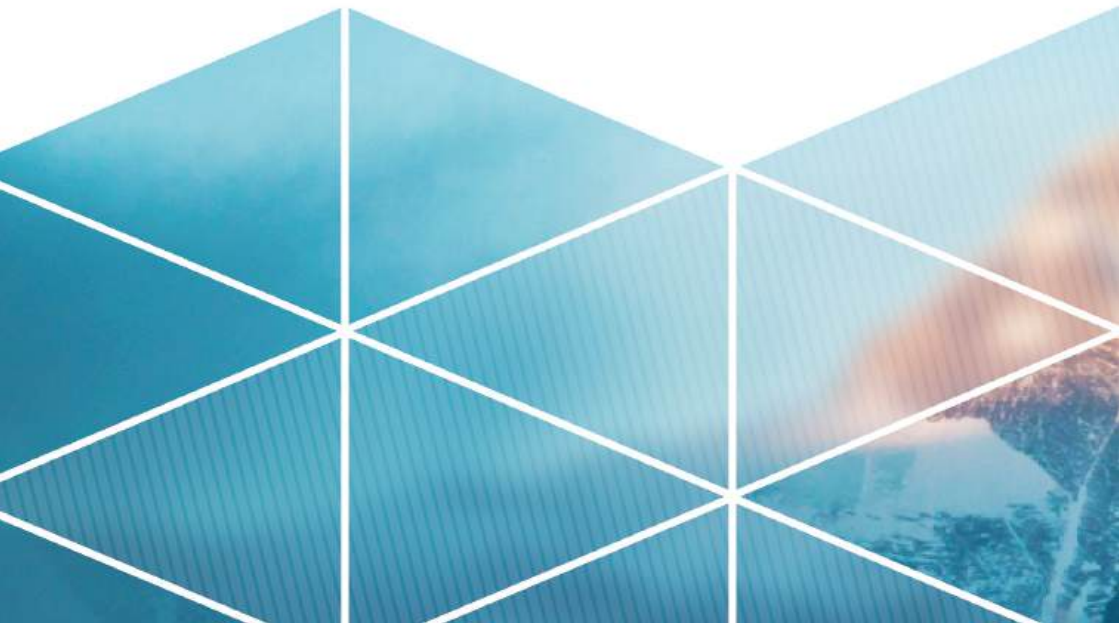





## COMPLEX HEART AND CIRCULATION

### HOW DANGEROUS ARE HEART AND VASCULAR DISEASES?



The main function of the cardiovascular system (CVS) is to provide organs with nutrients, oxygen and energy. Also, decay products and toxins “leave” the blood from organs. The cardiovascular system works constantly, without interruptions, and the slightest malfunctions in its work lead to problems throughout the body.



CVS diseases occupy a leading position among other diseases. Risk factors are stress, negative habits, inactivity, which are seen in the life of almost every person. CVS diseases are also caused by birth defects, injuries, inflammations, metabolic disorders.



The most common diseases are coronary artery disease, hypertension, arrhythmia, and atherosclerosis. They can lead to heart attacks, strokes and hypertensive crises - the main causes of death. In developed countries, these problems appear not only in older age groups - cases of heart attack at 30-35 are known.



Problems of the cardiovascular system are dangerous because they pass asymptotically for a long time. And their first signs are similar to ordinary fatigue - shortness of breath, dizziness, nausea, migraines. Often a person seeks help when the problems with the heart and blood vessels are already in a developed stage.



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites, entering the human body, exert a toxic load on the cardiovascular system and impair its functioning. Microbes destroy the heart and blood vessels, migrating along the vascular system and settling in the heart muscle and its surrounding membranes. DeVita AP basic antiparasitic programs the parasitic most known pathogens, quickly, effectively and without side effects, removing the parasitic load.

A special program "No schistosomes (blood-flukes) & trematods" cleans the blood of parasites that affect the condition of blood vessels and the human heart.

The "Cleansing" program is used to cleanse the body after applying basic antiparasitic programs, removing toxins and improving overall well-being.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

"Healthy Day" automatic program has a beneficial effect on the functioning of the cardiovascular system, helping to open even small capillaries. Improves blood microcirculation, nutrition of tissues and cells. It activates the process of body self-regulation and increases the efficiency of the immune system. The automatic program "Complete Cleansing" helps in the restoration of the excretory system and in harmonizes the work of the body.

Special programs "Healthy heart" and "Cardiac rhythm normalization" restore the heart rhythm, reduce the influence of negative factors, provide normal blood supply to the heart. The "Without headache" program improves cerebral circulation. The "Normal blood pressure" autoprogram helps to normalize blood pressure.



**DeVita Energy** relieves from psycho-emotional stress, which is one of the causes of cardiovascular diseases.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating "junk food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- with a violation of the heart rhythm (arrhythmia)
- with problems with veins and arteries
- with circulatory failure
- with general problems of the cardiovascular system.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Healthy heart, Cardiac rhythm normalization	Daily afternoon, if needed.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Healthy Heart, Without headache, Cardiac rhythm normalization	Daily afternoon, if needed.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
	5		Special programs	No schistosomes (blood-flukes) & trematods	Daily in the evening after DeVita Ritm programs for 1-2 weeks .
6	No chlamydia			Daily in the evening after DeVita Ritm programs for 3-4 weeks.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed for the night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	One of the manual programs (Healthy Heart, Without headache, Cardiac rhythm normalization) or Normal blood pressure	Daily in the afternoon or at night, if needed.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites - Automatic	Daily in the evening after the DeVita Ritm programs, alternating 1 program for 1-2 weeks.	

# PREVENTION MEASURES FOR CARDIOVASCULAR SYSTEM PROBLEMS



Reduce your salt intake. Salt retains water, which increases blood volume and leads to increased pressure.



Watch your blood pressure. Normal indicators are 120/80 mm Hg.



Once a year do a cardiogram.



Watch your blood glucose and cholesterol levels. Excess sugar and cholesterol cause the creation of deposits on the walls of blood vessels and clog them.



Treat tonsillitis and other bacterial infections in time. Chronic inflammation puts a strain on the heart and causes intoxication.



Eat boiled or steamed fish and poultry, avocados, linseed oil, nuts, pumpkin, broccoli, pomegranate, garlic and onions, ginger, turmeric, oat and corn bran, decoctions of hawthorn and rose hips, dried apricots and figs.



Limit refined oil, fatty and fried, sugar and salt.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Fatigue when standing up and walking fast. General fatigue, memory and efficiency loss.

After the first  
use of DeVita



Increased overall body tone,  
heart rate alignment.

After a month of  
using DeVita



Increased endurance, improved  
vascular and memory status,  
lack of shortness of breath.



Use the "HEART AND CIRCULATION" Complex scheme together with DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.



# COMPLEX STOP STRESS

## WHY IS STRESS DANGEROUS?



Stress is the organism's response to tension, negative emotions or excitement. During periods of severe or prolonged stress, the immune, endocrine, autonomic nervous systems suffer. Stressful situations are the cause of many physical and mental illnesses. The danger of stress is that its effects can accumulate gradually.



With severe or prolonged stress: The concentration of attention worsens, migraines appear; Problems with the cardiovascular system arise; Peptic ulcer and gastritis worsen; Immunity decreases In especially severe cases, degeneration of the cells of the spinal cord and brain occurs.



It is hard to work under stress and engage in habitual activities. Irritability, insomnia, and neurosis appear. Physical stamina is reduced. It becomes difficult to concentrate on the usual things – whether it is at work or at home. Relations with loved ones also suffer. An addiction to alcohol and/or drugs may appear.



Severe stressful situations require serious treatment – from psychotherapy sessions to taking medications that can be addictive. The way out of severe stress can last several years. It is much more effective to notice the first signs of stress in time or to prevent stressful situations.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Stress reduces the body's defenses, which is a favorable factor for infestation with parasites. They negatively affect the nervous system, causing fatigue, irritability and impair sleep quality. DeVita AP basic antiparasitic programs work against most known pathogens, quickly, efficiently and without side effects, removing the parasitic load.

"Cleansing" program is used to cleanse the body after using basic antiparasitic programs, removing toxins and improving overall well-being.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

The "Healthy Day" and "Health" automatic programs increase immunity and harmonize the work of all body systems. Help with acclimatization and stressful situations. The automatic "Complete Cleansing" helps in the restoration of the excretory system and in the harmonization of the body.

Special programs "Healthy sleep", "Good mood" and "Balanced nervous system" restore the nervous system, relieve nervous tension, improve the process of falling asleep and the quality of sleep.



**DeVita Energy** relieves stress and chronic fatigue. Restores biorhythms. Launches the functions of regulating the body, helping a person to be in good shape, maintain vigor and good mood throughout the day.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating "junk food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- in case of problems with sleep - insomnia, worsening sleep quality;
- when you feel tired in the morning, lethargy and weakness during the day;
- with nervous strain, severe excitement, frequent manifestations of irritability, apathy.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
			Special programs	Good mood, Balanced nervous system, Healthy sleep	If necessary, 2-3 hours before bedtime.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
				Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
				Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
				Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
			Special programs	Good mood, Balanced nervous system, Healthy sleep	If necessary, 2-3 hours before bedtime.
4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed for the night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
			Special programs	Good mood, Balanced nervous system, Healthy sleep	If necessary, 2-3 hours before bedtime.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites - Automatic	At night, alternating every other day for 2 weeks	

# STRESS PREVENTION MEASURES



Do not suppress spontaneous emotions. Suppressed anger or anxiety leads to stress.



Eat right. Add oranges and bananas, oatmeal, sea fish, herbal tea, avocado, nuts, fresh vegetables, dark chocolate to your diet.



Go to bed no later than 23:00 p.m. Sleep for at least 7-8 hours.



Do meditation or special breathing exercises.



Reduce your coffee, tea, and other caffeinated drinks.



Exercise. Exercising causes the production of endorphins - hormones of good mood, and also the adrenaline levels get higher.



Non-stressful monotonous work or creativity helps with stress - knitting, coloring pictures, etc.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Aggressiveness,  
fatigue, apathy,  
insomnia.

After first use  
of DeVita



Tranquility,  
mood improvement,  
burst of energy.

After a month of  
using DeVita



Positive mood,  
energy,  
poise.



Use the recovery scheme of the complex "STOP STRESS" together with the DeLixir MULTI-ENERGY and DeLixir pH BALANCE smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX STOP INSOMNIA

### WHAT DANGER DOES INSOMNIA HAVE?



The mechanism of sleep is not fully understood, but its primary importance has been proven. During sleep, recovery processes, the structuring of information received during the day, and physical relaxation occur. It is during this period that information is assimilated, tissues of the body are updated and restoration of immunity occurs.



Without a full night sleep, the body does not have time to recover and malfunctions. First of all, the nervous system suffers. Then irritability, aggressiveness, lethargy, apathy and reduced performance occur.

The mood drops sharply, a person becomes uncommunicative. This leads to numerous conflicts at work and at home.



Then physiological problems appear. A sharp decrease in immunity occurs, since immune cells -

T-lymphocytes - are synthesized during sleep. Body weight increases due to the production of the stress hormone - cortisol. Chronic lack of sleep leads to frequent illnesses, fatigue and poor health.



The lack of a full and healthy sleep significantly affects the appearance. Collagen production decreases in the skin of the face, wrinkles appear, and elasticity decreases. The complexion acquires an earthy hue, swelling, dark circles and wrinkles appear under the eyes. Irritations or rashes may occur.



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Microbes and larvae of large parasites can inhabit the organs and tissues of the body, including the nervous system. Their presence and the toxins secreted by them disrupt the system, leading to rapid fatigue, a state of nervous excitement and, subsequently, sleep disturbances. DeVita AP basic antiparasitic programs affect most known pathogenic organisms, quickly, efficiently and without side effects, removing parasitic load.

“Cleansing” program is used to cleanse the body after applying basic antiparasitic programs, removing toxins and, as a result, reducing nervous system excitement and improving overall well-being.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” and “Health” automatic programs improve the functioning of the excretory system, activate the circulatory and immunity systems. Help in coping with physical, mental and psychological stress, create favourable conditions for healthy sleep. The automatic program “Complete Cleansing” helps in the restoration of the functions of the excretory system and harmonizes the body.

Automatic program “Healthy sleep - Automatic” reduces nervous excitement, promotes fast falling asleep and good quality sleep.



**DeVita Energy** relieves from psycho-emotional stress, stabilizes the nervous system, helps in restoring energy. It positively affects the quality of sleep, gives vigor when waking up.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water for at least 1-1.5 l / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- problems with falling asleep;
- sleep disorders (light sleep, waking up in the middle of the night, etc.);
- state of nervous excitement.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
			Special programs	Healthy sleep - Automatic	If necessary, 2 hours before bedtime.
	3	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
				Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
				Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
				Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating every other day.
			Special programs	Healthy sleep - Automatic	If necessary, 2 hours before bedtime.
3	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed for the night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
			Special programs	Healthy sleep - Automatic	If necessary, 2 hours before bedtime.
3	DeVita AP	Special programs	No nematodes and ascarides, No parasites - automatic	At night, alternating every other day for 2 weeks	

## PREVENTIVE MEASURES FOR SLEEP DISORDERS



Maintain a comfortable temperature for sleeping in the bedroom - 19-20 ° C. Walk before bedtime and ventilate the room, and in the warm period, sleep with the window open.



Create a ritual of going to bed: take a hot bath, listen to music, relax. Try not to do work or daily routine in the bedroom.



Try not to sleep during the day. Or at least reduce your daytime sleep to 20-30 minutes.



The key to a good sleep is complete silence and darkness. Use blackout curtains and for protection against extraneous sounds - ear plugs.



A way of falling asleep easier is by using "white noise" sounds - monotonous sound of the radio on an "empty" frequency. Other examples are the sounds of waterfall, fan, wind, rain etc.



Exercise and eat no later than 2 hours before bedtime.



Remove flowers and other plants from the bedroom: they produce oxygen only in the daytime, and at night, on the contrary, they absorb it.



If possible, use blue and purple tones in the design of the bedroom.



Include in your evening menu products with a calming effect: honey, milk, bananas, pumpkin, date fruits.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before  
use



Difficulty in falling asleep, drowsiness, restlessness, apathy.

After the first use  
of DeVita



Emotional and energetic recovery, improvement of general condition.

After a month of  
using DeVita



Cheerfulness, improved efficiency, easy awakening, improved communication.



Use the scheme of the "STOP INSOMNIA" complex together with the DeLixir MULTI-ENERGY and DeLixir DETOX smart cell nutrition supplements, and your results will become noticeable much faster.



## COMPLEX WOMEN'S HEALTH

### WHY ARE WOMEN'S HEALTH PROBLEMS DANGEROUS?



Women are more adapted to long-term stress than men.

But the female body has its own resources, and in case of serious prolonged loads, female health is weakened. Ecology has a huge impact on this matter, as every 5th female resident of any city is faced with "female" problems.



The harbingers of the failure of the female body are: fatigue, bad mood, irritability, constant prolonged fatigue, unwillingness to get up in the morning and loss of interest in intimacy with a man.

If you ignore these changes, they can lead to serious health problems.



Failures in the female genitourinary system lead not only to disruption of the cycle and discomfort, but also to problems with pregnancy, inflammation and even tumors. Common female problems include an irregular cycle, cervical erosion, uterine inflammation and inflammation of fallopian tubes, cystitis, bacterial and fungal infections.



Problems with female health lead to negative changes in appearance.

The skin becomes more oily and body hair growth is increased. Inflammation and acne appear.

Failure in the production of female hormones can lead to rapid weight gain and hair loss.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites and microbes cause inflammation of the organs of the female genitourinary system and even lead to infertility and chronic kidney diseases. DeVita AP basic antiparasitic programs affect most known pathogens, quickly, effectively and without side effects, removing the parasitic load from the organism.

“Cleansing” automatic program is used to cleanse the body after applying basic antiparasitic programs, remove toxins, and has an antiseptic effect.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

The “Healthy Day” and “Health” automatic programs harmonize the work of the body, the energy of the chakras, activate the circulatory system, digestion processes, the excretory system, normalize the functions of the genitourinary system and improve the quality of intimate life. The automatic program “Complete Cleansing” contributes to the restoration of the main excretory system and the harmonization of the body.

Special programs “Female sexuality”, “Women’s health”, harmonize sexual desire, regulate the reproductive system, facilitate the condition of women during menstruation and menopause.



**DeVita Energy** relieves from psycho-emotional stress, restores and strengthens physiological indicators. It gives energy and vitality, fights against sexually transmitted diseases, removes psychological burden.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- violation of the menstrual cycle;
- severe menstrual pain;
- to harmonize the work of the body during menopause;
- decrease in sexual desire.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.









# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Women's health	Daily in the afternoon.
	4			Female sexuality	If necessary, according to guidelines.
	5	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	6			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	7			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	8			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Women's health	Daily in the afternoon.
	4			Female sexuality	If necessary, according to guidelines.
	5	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
6	Special programs		No chlamydia	Daily in the evening after DeVita Ritm programs.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Women's health	Daily in the afternoon.
	4			Female sexuality	If necessary, according to guidelines.
5	DeVita AP	Special programs	No nematodes and ascarides, No parasites - automatic	At night, alternating every other day for 2 weeks	

-  Visit your gynecologist and urologist at least once per year.
-  Observe hygiene. Use special cleansers - regular soap disrupts the microflora of the genitourinary system.
-  Take hormonal drugs and contraceptives only if they are prescribed by your doctor. Do not take medicines, especially antibiotics, for longer than 7 days without special instructions.
-  Have a regular sex life, be sure to use means of protection against sexually transmitted diseases.
-  Protect the genitourinary system from excessive heat or cold.
-  Do not lift weights and control physical activity to avoid prolapse of the uterus and other problems.
-  Walk at least 20 minutes a day to avoid stagnation in your pelvic organs.
-  Exclude sports and sexual intercourse during menstruation.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Bad mood, decreased libido, discomfort and cycle disruptions.

After the first use  
of DeVita



Burst of energy, decrease of pain, elimination of irritability and apathy.

After a month of  
using DeVita



Improved mood, normalization of the reproductive system, increased libido.



Use the scheme of the "WOMEN'S HEALTH" complex scheme together with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX MEN'S HEALTH

### WHY ARE MALE HEALTH PROBLEMS DANGEROUS?



Even young men may have problems with erection due to unhealthy lifestyle and urological diseases, the most common of which is prostatitis. Another dangerous disease is prostate adenoma, a benign tumor of the prostate gland.

Every second male over age of 50, suffers from prostate adenoma.



Acute prostatitis is difficult, with increased body temperature up to 39-40 ° C. Urination is difficult and very painful. If untreated, prostatitis becomes chronic. It can affect the bladder and kidneys, cause inflammation of the prostate and sepsis. And the source of infection can not always be determined.



Prostate adenoma occurs due to a combination of reasons - age, hormonal disorders, obesity, stress, etc. Prostate adenoma is dangerous and complications cause: acute urinary retention, urinary incontinence, problems in sexual life, and infections of the genitourinary system. Adenoma often requires surgical treatment.



Due to prostatitis and prostate adenoma, not only the urination process is disturbed, but also the potency. Also, the complications of these diseases can lead to infertility. Often men are faced with psychological disorders - this leads to depression, neurosis, chronic stress and other neurological diseases.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites cause inflammation of the genitourinary system. Diseases proceed unnoticed by humans, having a serious impact on the function of the urinary and reproductive systems. DeVita AP basic antiparasitic programs affect most known pathogens, quickly, effectively and without side effects, removing the parasitic load from the organism.

“Cleansing” automatic program is used to cleanse the body after applying basic antiparasitic programs, remove toxins, and has an antiseptic effect.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” and “Health” automatic programs harmonize the work of the organism, normalize the functions of the genitourinary system and improve the quality of intimate life. The automatic program “Complete Cleansing” helps in main excretory systems restoration and harmonizes the organism.

Special programs “Male sexuality”, “Men’s Health” and “Stop prostatitis” increase self-confidence and sexual desire, regulate sexual activity, ensure the normal functioning of the genitourinary system.



**DeVita Energy** relieves from psycho-emotional stress, restores and strengthens physiological indicators. Gives energy, reduces fatigue accumulated throughout the day.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the device. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- violation of potency;
- psychological problems associated with the “neurosis of expectation of failure”;
- problems with the bladder and prostate.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Men's health or Men's health - automatic	Daily in the afternoon.
	4	Male sexuality		If necessary, according to guidelines.	
	5	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	6			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	7			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	8			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Men's health or Men's health - automatic	Daily in the afternoon.
	4	Male sexuality		If necessary, according to guidelines.	
	5	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
6	Special programs		No chlamydia	Daily in the evening after DeVita Ritm programs.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Men's health or Men's health - automatic	Daily in the afternoon.
				Male sexuality	If necessary, according to guidelines.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites - automatic	At night, alternating every other day for 2 weeks	

## PREVENTIVE MEASURES FOR UROGENITAL PROBLEMS



Maintain a regular and stable sex life. Excessive activity and prolonged abstinence lead to disruption of the genitourinary system.



Exercise at least 20 minutes a day.



Avoid hypothermia, especially in the groin area.



When changing your regular partner, be sure to use condoms. Sexually transmitted infections can cause prostatitis.



Keep a track of your bowel movements. Frequent constipation leads to prostatitis.



4 times a year, have a massage of the lower back area.



Have a routine examination by a urologist 2 times a year.



Diversify your diet with fresh fruits, vegetables, herbs, seafood, nuts and seeds.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Pain in the lower abdomen,  
frequent urges to urinate,  
sexual problems.

After the first use  
of DeVita



Pain relief,  
mood improvement,  
surge of energy.

After a month of  
using DeVita



Stable genitourinary system  
functioning, stable erection,  
self-confidence.



Use the scheme of the "MEN'S HEALTH" complex in conjunction with DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX MOVEMENT WITHOUT PAIN

### WHY SPINE AND JOINT PROBLEMS ARE DANGEROUS?



Problems with spine and joints are a fee that's humans have to pay for upright posture. Having only two points of support, the human body is forced to constantly maintain balance. Strictly speaking, each person is fated with problems with the spine. Because of this, the joints experience increased stress.



We exacerbate these problems by leading a sedentary lifestyle, lifting weights, or not paying attention to posture. A huge role is played by shoes, which are designed to soften the load on the spine, but do not always fulfill this function. Meanwhile, over 70% of people have spinal curvature.



The most common problems are impaired mobility, severe pain and inflammation. And this indirectly affects the appearance, causing obesity, varicose veins, problems with gastrointestinal tract and general fatigue. In especially severe cases, problems with the spine and joints lead to restricted movement and disability.



Diseases of the back and joints can go unnoticed for a long time, but at the most unexpected moment declare about themselves with severe pain. Treatment can, basically, only relieve symptoms - reduce pain, improve mobility and blood circulation. Therefore, it is important to take care of the health of the spine and joints on time.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

The parasite larvae that enter the body settle on the membrane lining the joint cavity, or in the periarticular tissues, causing inflammation. Parasites also affect the spine. DeVita AP basic antiparasitic programs affect most known pathogens, quickly, effectively and without side effects, removing the parasitic load from the organism.

“Cleansing” automatic program is used to cleanse the body after applying basic antiparasitic programs, reduce the parasitic intoxication and improve overall well-being.

Special programs “No bacteria” and “Drainage Antitox” help to reduce toxic load in the presence of specific groups of microorganisms and salts of heavy metals in the body that have a negative effect on the musculoskeletal system.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” and “Health” automatic programs improve blood circulation, gently reduce the toxic load, which enhances immunity, improve the outflow of metabolic products from damaged areas. “Complete Cleansing” automatic program helps in main excretory system restoration and harmonizes the organism.

The programs “Recovery from injuries”, “Healthy joints - automatic” and “Joints without pain - automatic” relieve from inflammation, pain, swelling, activate metabolic processes in damaged tissues, normalize muscle tone, improve the functioning of the ligamentous apparatus.



**DeVita Energy** triggers self-regulation mechanisms, restores physiological processes in organs and tissues. Gives vigor and energy for movement.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- pain in the joints and spine;
- injuries and inflammation of the muscles and joints;
- pinched nerves.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Joints without pain - automatic or Healthy joints and / or Recovery from injuries	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Joints without pain - automatic or Healthy joints and / or Recovery from injuries	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night
5	Special programs		No bacteria, Drainage Antitox	Daily in the evening after DeVita Ritm programs for 14 days.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Joints without pain - automatic or Healthy joints and / or Recovery from injuries	Daily in the afternoon, if necessary.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites - automatic	At night, alternating every other day for 2 weeks	

## PREVENTIVE MEASURES FOR SPINAL AND JOINT PROBLEMS



Use vitamins and minerals in spring and autumn.



Watch your posture. Lean your back against a flat wall - your nape, shoulder blades, buttocks and feet should touch the wall. If necessary, consult a doctor and purchase a corrective corset.



Sleep on an orthopedic mattress and pillow. In extreme cases - on a dense hard surface. Fall asleep lying on your back.



Exercise and strengthen your back muscles and abs. Developed abdominal muscles relieve pressure on the lower back.



When doing sports and lifting weights, follow the correct technique - do not round your back, strain your abs.



Eat sour-milk products rich in calcium, as well as products with gelatin content - meat with gelatin and jelly.



Consult a doctor if you want to try a new physical activity - yoga, step aerobics, strength training, etc.



Be wary of osteopathy, manual therapy, mechanical massage.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Pain in the back and joints,  
fatigue, headache.

After the first use  
of DeVita



Decreased pain  
and tiredness.

After a month of  
using DeVita



Fast recovery,  
ease of motion,  
no pain.



Use the scheme of the "MOVEMENT WITHOUT PAIN" complex together with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX CLEAN SKIN

### WHAT IS CAUSING SKIN PROBLEMS?



The condition of the skin is an indicator of human health. Failures in the work of some body systems: digestive, endocrine, immune, circulatory, are reflected on the skin of the face and body. Also, the condition of the skin and hair is affected by ecology, stress, the use of low-quality cosmetics.



The most common facial skin problems are acne and blackheads, age spots, and rosacea. Basically, these problems have age-related causes associated with hormone imbalance. However, facial skin conditions may worsen at an earlier age. The reason for this may be improper skin care or the selection of poor-quality cosmetics.



Common scalp problems are dandruff and seborrhea, which cause excessive oiliness or dryness of the skin. This affects the health of the hair - they become thin, brittle and weak, begin to fall out. Problems with the scalp can be dermatological, they are also often caused by microorganisms - bacteria and fungi.



Body skin problems are similar to scalp problems. Those are dryness in the winter season, acne and blackheads, especially in the back, low skin tone and sagging. These problems are caused by both external pathogens and internal malfunctions in the body, as well as improper care, unbalanced diet and daily routine.



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Viral, bacterial, fungal and other types of infections cause problems with the skin, hair, and nails. Dandruff, the formation of crusts, scales, pustules, allergic skin reactions often indicate the presence of pathogenic organisms. The DeVita AP basic antiparasitic programs affect the majority of known pathogens, quickly, efficiently and without side effects, removing the parasitic load and signs of inflammation.

"Cleansing" automatic program is used to cleanse the body after applying basic antiparasitic programs, reduces parasitic intoxication and improves both overall well-being and skin condition.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

"Healthy Day" and "Health" automatic program improve blood circulation, normalize the functioning of the gastrointestinal tract. Gently reduce toxic load, improve the nutrition of skin cells and hair follicles, normalize cellular metabolism and activate regeneration processes. The automatic program "Complete Cleansing" contributes to the restoration and harmonization of the body.

The programs "Beautiful Skin", "Beautiful Hair and Nails" relieve from inflammation, activate metabolic processes, improve skin tone, strengthen hair and nails, slow down the aging process.



**DeVita Energy** launches self-regulation mechanisms, restores physiological processes in organs and tissues, enhancing the action of the DeVita AP and DeVita Ritm programs.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating "junk food" (chips, crackers, etc.) – this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- inflammatory processes on the skin, acne;
- deterioration of hair and nails;
- hair loss, dandruff.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Beautiful skin and / or Beautiful hair and nails	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Beautiful skin and / or Beautiful hair and nails	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
5	Special programs		No fungi, No herpes zoster, No herpes simplex	If necessary, according to guidelines.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Beautiful skin and / or Beautiful hair and nails	Daily in the afternoon, if necessary.
4	DeVita AP	Special programs	No nematodes and ascarides, No parasites - automatic	At night, alternating every other day for 2 weeks	

# PREVENTIVE MEASURES FOR PROBLEMS WITH FACE, BODY AND HEAD SKIN



With pronounced problems that arise for no apparent reason, conduct an examination of the body for internal pathologies.



Do not use other people's cosmetic and hygiene items, do not try on other people's hats.



Try to avoid prolonged sun exposure, limit your visits to the solarium. Use cosmetics with UV filters.



Monitor the balance of hormones, with a predisposition to hormonal disruptions, undergo an examination at least once a year.



Use only high-quality cosmetics, monitor their expiration date. Thoroughly clean your face with a special cleanser in the morning and evening.



Do not perform scrubs and mechanical cleanings frequently. Frequent exfoliation causes irritation and dryness.



In the cold period, be sure to moisturize the skin after a shower, use thermal water for the face every few hours.



Red fish, eggs, greens, nuts and legumes are useful for good skin and hair condition.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before  
use



Acne, dandruff,  
age spots,  
sagging skin.

After the first use  
of DeVita



Improved skin tone,  
reduced rashes,  
regulated oil balance.

After a month of  
using DeVita



Dark spots and inflammation  
becomes less noticeable,  
dandruff disappears.



Use the scheme of the "CLEAN SKIN" complex with the DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.



## COMPLEX SPORT PERFORMERS

### HOW DO TRAINING AND SPORTS AFFECT THE ORGANISM?



Physical activity seriously changes the metabolism and vital functions of the body. Athletes need more nutrients and vitamins, relaxation and full sleep. Bones and joints experience increased stress. Metabolism is accelerated, the need for fluid intake increases.



When doing fitness or sports, we want to get high results in a short time. However, the body's recourses are limited. Due to active metabolism during training, a large expenditure of energy occurs, the accumulation of decay products, which are removed only after 24 hours, increases the toxic load on the liver and kidneys.



Intense and active workouts cause post-workout muscle pain and a feeling of fatigue. In addition, injuries are common in sports. For best results, you need to resume training only after a full recovery, but this is fraught with a significant loss of time. And in sports, every day training is important.



Each training or competition is a great stress for the body. The accumulated stress negatively affects the general state of health and can reduce working capacity, and in the future - cause heart problems. The success of training is influenced by psychological factors - mood, personal relationships.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites and microorganisms can significantly undermine health. This affects the quality of training and recovery. DeVita AP basic antiparasitic programs affect the majority of known pathogens, quickly, efficiently and without side effects, removing the parasitic load and signs of inflammation.

“Cleansing” automatic program is used to cleanse the body after applying basic antiparasitic programs, reduces parasitic intoxication and improves overall well-being.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” and “Health” automatic programs improve blood circulation, facilitating quick recovery after training. The work of the digestive tract is normalized, due to which the load on the liver is reduced, the elimination of toxins is accelerated. “Complete Cleansing” automatic program helps to restore and harmonize the body.

The “Warm-up” and “Workout” programs help to increase strength indicators, allow you to train longer and more intensively, give energy and vigor. The programs “Recovery from injuries”, “Joints without inflammation”, “Joints without pain” relieve from inflammation and pain, normalize muscle tone, improve the functioning of the ligamentous apparatus, and help recover after an injury.



**DeVita Energy** activates the mechanisms of self-regulation, restores the physiological processes in organs and tissues. Due to this, recovery is much faster.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- intense physical exertion or sports training;
- pain in the joints, muscles and spine after exertion;
- injuries and inflammation of the organs of the musculoskeletal system;
- for the best recovery after training.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Warm-up, Workout	Accordingly, before and during training.
				Joints without pain and / or Joints without inflammation and / or Recovery from injuries	Daily in the afternoon, if necessary.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
7	Mode 4 automatic			Daily at night during the 4 <sup>th</sup> week.	

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic	Daily in the morning.
	3		Special programs	Warm-up, Workout	Accordingly, before and during training.
4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Warm-up, Workout	Accordingly, before and during training.
	4			Joints without pain and / or Joints without inflammation and / or Recovery from injuries	Daily in the afternoon, if necessary.
5	DeVita AP	Special programs	No nematodes and ascarides, No parasites - automatic	At night, alternating every other day for 2 weeks.	



Start your workout with more complex exercises for large muscle groups, and end with isolated or cardio exercises.



Watch your pulse. The maximum heart rate is calculated from the formula  $214 - \text{age} \times 0.8$ . The pulse during cardio training should be 50-59% of the maximum, with light fitness - 60-69%, during endurance training - 70-79%, during anaerobic training - 80-89%.



Eat right. In terms of quality composition, proper nutrition for an athlete should be close to the formula: 30% - proteins, 60% - carbohydrates, 10% - fats.



Use natural energy sources - ginseng extract, guarana, succinic acid.



Variety of exercises and training methods - muscles quickly get used to uniformity.



If possible, train at the same time. Choose it depending on your daily routine and activity level.



If you have severe muscle pain, take a hot bath and do self-massage.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before  
use



Fatigue, muscle tension  
and pain,  
long recovery time.

After the first use  
of DeVita



Burst of energy,  
reduction of pain.

After a month of  
using DeVita



Increased efficiency,  
quick recovery, improved  
stamina and strength.



Use the scheme of the "SPORT PERFORMERS" complex together with the DeLixir MULTI-ENERGY, DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX STOP ALLERGY

### HOW DANGEROUS ARE ALLERGY MANIFESTATIONS?



Allergy is a disrupted mechanism for the implementation of the body's immune response. The body mistakenly takes a substance as a threat to life and forces the immune system to produce antibodies and various chemicals, one of which is histamine, which causes allergy symptoms.



In addition to food, wool and pollen, the role of an allergen can be played by chemical substances, some metals, alcohol, cold and even sunlight. There is also an allergy to substances formed inside the body itself. Often the symptoms are similar to the common cold or usual fatigue, making diagnosis difficult.



Allergic reactions may be the cause of development of bronchial asthma, hay fever, allergic dermatitis. The most serious complications are convulsions, acute vascular insufficiency, stomach cramps. The most severe allergic reactions - anaphylactic shock and Quincke's edema - lead to instant death.



According to statistics from the World Health Organization, about 40% of the world's population suffer from allergies. Diagnosing and identifying allergens is a complex and time-consuming process, and treatment is reduced to reduce symptoms and alleviate the general condition of the body, but does not eliminate the allergy in general.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Bacteria, fungi, viruses, helminths are biological allergens. The basic antiparasitic programs of DeVita AP affect most microorganisms, quickly, efficiently and without side effects, removing the pathogenic load.

Parasitic burdening of the liver and blood reduces their functions, enhances intoxication and allergic reactions. Special programs “No schistosomes (blood-flukes) & trematods”, “No allergies”, “No herpes” relieve the pathogenic load.

The “Cleansing” program is used to cleanse the body after applying antiparasitic programs, the removal of toxins, which can also be allergens.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” and “Health” automatic programs have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. Strengthens the immune system, increases vitality. Special programs “Without allergens” and “Anti-cold” relieve from itching, reduce lacrimation, swelling of the mucous membranes, and activate immune defense in case of allergies. The automatic program “Complete Cleansing” helps in the restoration of the basic functions of the excretory system and in the harmonization of the body.



**DeVita Energy** harmonizes the work of the body, energizes, activates the process of elimination of toxins. Improves the general condition when allergic reactions manifest, reduces lethargy and drowsiness arising when taking antihistamines.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

Allergy is a consequence of intoxication of the body. To minimize its manifestations, it is recommended to carry out cleansing measures within 6 months. In the first three months it is necessary to use the devices according to schemes 1-3, in the next 3 months - alternate schemes 2 and 3. When cleansing and restoring the body, minor exacerbations are possible, in this case, a break in the use of devices is necessary.

In acute and severe conditions, special programs can be included several times a day.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices together with DeLixir DETOX and DeLixir pH BALANCE smart cell nutrition supplements (if there is no food allergy).

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy or Energy Pro		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Without allergens	Daily in the morning, alternating with each other every other day.
	4	Anti-Cold, Free breath		If necessary, according to guidelines.	
	5	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	6			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	7			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	8			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.
	9			No allergies	In the evenings after the DeVita Ritm programs.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Without allergens	Daily afternoon.
	4	Anti-Cold, Free breath		If necessary, according to guidelines.	
	5	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
	6		Special programs	No allergies	In the evenings after the DeVita Ritm programs.
7	No schistosomes (blood-flukes) & trematods			In the evenings after the DeVita Ritm programs, 14 days.	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Complete Cleansing-automatic	Daily from 8 a.m.
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If necessary, at night.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Mode program		1-2 times in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3		Special programs	Without allergens	Daily afternoon.
	4	Anti-Cold, Free breath		If necessary, according to guidelines.	
	5	DeVita AP	Special programs	No allergies	In the evenings after the DeVita Ritm programs.
6	No nematodes and ascarides, No parasites - automatic			At night, alternating every other day for 2 weeks	

## PREVENTION MEASURES FOR ALLERGIC REACTIONS



If you already have an allergy to any substance, go through a full allergen test. Often an allergy-prone organism cannot tolerate several foods.



Check the gastrointestinal tract. The presence of helminths, dysbiosis and diarrhea problems exacerbate the course of allergic reactions.



If you are allergic to dust, store bedding, pillows, and books in tightly closed cabinets. Wash clothes weekly at a temperature of at least 60 °C.



When hay fever starts, perform frequent wet cleaning of your personal and living space every day. Take off your outdoor shoes and clothes at the doorstep - pollen and plant seeds accumulate on them.



For food allergies, avoid foods with strong allergenic properties: honey, peanuts, fish and seafood, spices, chocolate, citrus fruits, tomatoes, strawberries.



Use the medicine only as directed by your doctor. Keep track of the expiration date of drugs - expired drugs emit volatile allergenic substances.



If you are allergic to mold, do not close the doors to the bathroom and the kitchen. Throw away moldy foods. Check the flower pots - mold may form in them.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Runny nose, itching,  
swelling, skin rashes,  
lacrimation.

After first use  
of DeVita



Reduced itching, relieve  
from swelling, better general  
condition.

After a month of  
using DeVita




Reduced allergy manifestation,  
reduced rashes,  
breathing relief.



Use the recovery scheme of the "STOP ALLERGY" complex together with DeLixir DETOX and DeLixir pH BALANCE smart cell nutrition supplements (if there is no food allergy), and your results will become noticeable much faster.



# PARTICULAR COMPLEXES

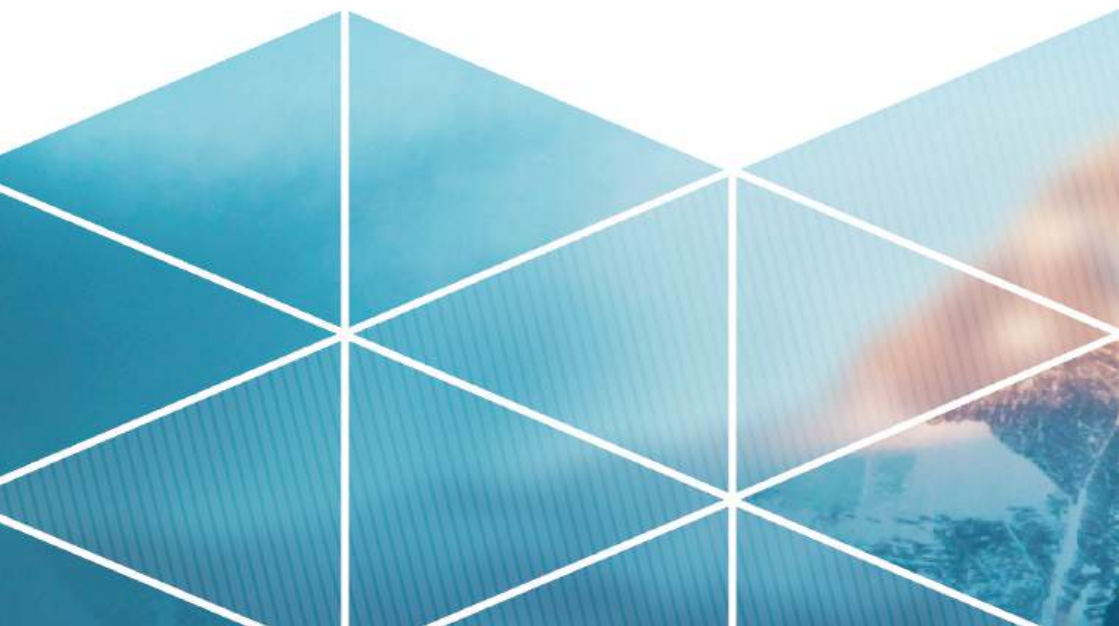



## COMPLEX WITHOUT COLD

### WHY ARE DIGESTIVE DISEASES DANGEROUS?



The immune system protects the body, participates in its life activity and in adaptation processes to changing conditions.

Immune cells should always be on guard, because the body is attacked by pathogenic microbes every second.

To some of their species, long-term immunity is developed, but not to colds and flu viruses.



Frequent catarrhal diseases lead to disorders

in the work of organs and an even greater decrease in immunity. Microbes freely enter the body and destroy it.

As a result, the common cold results in serious complications. They require outpatient treatment - from 2 weeks to several months of bed rest.



The most common complications are inflammations:

sinuses of the nose (sinusitis, frontal sinusitis); middle ear (otitis media, labyrinthitis); respiratory organs (bronchitis, pneumonia, pleurisy).

Every complication affects major organ systems and causes severe intoxication.



Medicines against flu and colds also cause negative reactions.

Antipyretic drugs reduce immunity and contribute to inflammation. Vasoconstrictor nasal drops affect the mucous membrane and small vessels, causing chronic edema. Antitussive drugs produce sputum and respiratory infections.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

DeVita AP's basic antiparasitic programs affect most known pathogens. The use of the programs helps in coping with the disease quickly, efficiently and without side effects. In case of a long and severe course of the disease, combine DeVita AP programs from schemes 2 and 3.

The "Cleansing" program is used to cleanse the body after basic antiparasitic programs, remove toxins, strengthen immunity and recover from the disease. As a result, overall well-being improves.



**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

Automatic programs "Healthy Day" and "Health" improve local and general blood circulation, digestive organs work, which is important for non-ideal diet. Help to boost immunity and harmonize the work of all systems. Help with acclimatization and stressful situations. Automatic program "Without cold - automatic" normalizes the bronchopulmonary system, relieving inflammation of the upper respiratory tract, relieving coughing. The "Throat without problems" program relieves throat discomfort.



**DeVita Energy** relieves stress and chronic fatigue. It neutralizes the effects of negative environmental factors and domestic radiation, helps with weather sensitivity. Launches the body's regulation functions, helping a person to be in good shape, maintain vigor and good mood throughout the day.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - the impact of the devices. Remove or minimize coffee in your diet. Avoid eating "junk food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY and DeLixir pH BALANCE smart cell nutrition supplements.





## USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

SCHEME №1: in the season of colds and flu, in close contact with sick people, with overcooling

	DEVICE	PROGRAM	USAGE
SCHEME №1	1 DeVita Energy	Energy Pro, Health	Once in the morning.
	2 DeVita Ritm	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
	3 DeVita AP	Support immune system - automatic	1-2 times per week.

SCHEME №2 : at the first sign of bacterial and viral infections (SARS, ARI, tonsillitis, rhinitis, etc.), except for flu.

	DEVICE	PROGRAM	USAGE
SCHEME №2	1 DeVita Energy	Energy Pro, Health	Once in the morning.
	2 DeVita Ritm	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
		Anti-Cold, Free breath	Daily afternoon, 5-7 days.
	4 DeVita AP	Antivirus + Without Toxins	First 2 days, after DeVita Ritm programs.
		Mode 2 automatic + Anti-Cold	Daily, alternating with each other for 5-7 days.
		Throat without pain and/or Against cold - General	If needed

SCHEME №3: at the first signs of flu.

	DEVICE	PROGRAM	USAGE
SCHEME №3	1 DeVita Energy	Energy Pro, Health	Once in the morning.
	2 DeVita Ritm	Healthy Day-automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
		Anti-Cold automatic and Throat without problems	Daily afternoon, 5-7 days.
	4 DeVita AP	Support immune system - automatic	First 2 days, after DeVita Ritm programs.
	5 DeVita AP	No flu - automatic	Daily, 5-7 days.
		Mode 2 automatic	A week after influenza programs, daily for 5-7 days.

# PREVENTIVE MEASURES DURING THE EPIDEMIC OF FLU AND COLD



Perform frequent wet cleaning of your personal and living space.



Ventilate the room every day.



While being outside, keep your legs, head and neck warm and dry.



Change your disposable hygiene mask every 3 hours.



Use disposable wipes to cover your nose and mouth when you sneeze or cough, and throw them away immediately after use.



Eat garlic, red bell pepper, cranberries, rose hips, carrots, sea fish and seafood, yogurt without additives, kefir, honey, broccoli.



In addition to juices, tea and other liquid products, be sure to drink plenty of clean water at least 1-1.5 l / day.



Moisten the mucous membranes with saline (1 tsp salt in 1 liter of water). Spray or drip it into your nose every 1-2 hours.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Fatigue in the morning, complications after illness, lethargy and weakness during the day.

After first use  
of DeVita



A burst of energy, improvement in well-being and mood.

After a month  
of using DeVita



Easy rise in the morning, vigor, increased adaptive reserves of the body, reduced morbidity.



Use the recovery scheme of the complex "WITHOUT COLD" together with the DeLixir MULTI-ENERGY and DeLixir pH BALANCE smart cell nutrition supplements, and your results will be noticeable much faster.



# COMPLEX STOP CORONAVIRUS- PREVENTION

## BRIEFLY ABOUT CORONAVIRUS



Coronaviruses are a family of viruses that includes 40 types of RNA-viruses from the group of coronaviruses of animal origin, that infect humans and animals. The name is associated with the structure of the virus, the structures on the virus itself, which resemble a crown, thanks to which the virus enters the cell of a person or animal.



Sources of coronavirus infections can be an infected person or animal. The main transmission routes: airborne, dust, fecal-oral, contact. The infection rate increases during winter and spring.



According to WHO, the virus spreads only within a radius of 1 meter around an infected person. In air, the virus retains its activity for 8-10 hours, in water – up to 9 days. In the most favorable conditions, but outside the human or animal body – 48 hours.



The most effective means to prevent infection are to treat your hands with alcohol-containing antiseptics and thoroughly wash your hands with soap. Do not touch with dirty hands your eyes, nose, and mouth. Masks are ineffective because they get wet from breathing and the virus penetrates through the moisture in the mask.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist discoverer Dr. Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency; it is the increase in amplitude which leads to the death of the microorganism without harming other cells.

DeVita AP antimicrobial programs help the body's efforts to destroy respiratory viruses and bacteria. The programs quickly and effectively help to cope with the infection, as well as help the body recover from the disease.

Special programs "Without cough" and "Without pneumonia" are used to reduce inflammation and facilitate the discharge of sputum from the respiratory tract.



**DeVita Ritm** – here the action is based on the studies of the scientists Dr. Reinhold Voll and A.F. Cramer, which found that each organ has a healthy frequency of oscillation. DeVita Ritm gently affects the body, pointing out the frequency of the vibrations of healthy organs. As a result, the body is capable of self-healing.

Autoprograms "Healthy Day" and "Health" have a beneficial effect on the work of the whole organism, strengthen defenses, increase vitality. They strengthen the defenses, increase vitality. Using the "Anti cold – automatic" automatic program helps to restore immunity and increase the energy potential of the respiratory system.



**DeVita Energy** harmonizes the function of the body, energizes, activates the removal of toxins, starts the processes of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING THE DEVITA DEVICES

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the effectiveness of the devices. Remove or minimize coffee in your diet.

Avoid eating "fast food" (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING CASES:

- if you have to communicate with sick people;
- if your work is connected with a large number of people;
- if you have been in countries with an unfavorable situation for this disease;
- if you often get colds;
- if you want to protect yourself from infection.



### COMBINE THE DEVICES WITH THE DELIXIR SMART CELL NUTRITION SUPPLEMENTS

For best results, it is recommended to use the DeVita devices in conjunction with the smart cell nutrition supplements DeLixir pH Balance, DeLixir Top Form, DeLixir Multi-Energy.

# APPLICATION SCHEMES

SCHEME №1 For infection prevention:

SCHEMA №1	DEVICE		PROGRAM		APPLICATION
	1	DeVita Energy		Energy Pro, Health	Once in the morning.
2	DeVita Ritm	Automatic programs	Healthy Day - automatic	1 <sup>st</sup> day, in the morning.	
3			Health - automatic	2 <sup>nd</sup> day, in the morning.	
4		Manual programs	Immune protection	Use after the automatic program "Healthy Day – automatic" on the DeVita Ritm device.	
5	DeVita AP	Automatic programs	Without cold - automatic	2-3 times a week, at night.	
6		Manual programs	Stop Coronavirus- Prevention- New	2-3 times a week after the use of the DeVita Ritm device.	

SCHEME № 2 During the infection, at its first symptoms:

SCHEMA №2	DEVICE		PROGRAM		APPLICATION
	1	DeVita Energy		Energy Pro, Health	Once in the morning.
2	DeVita Ritm	Automatic programs	Healthy Day - automatic	Alternate every second day for 2 weeks. Use in the morning.	
3			Health - automatic		
4			Anti cold - automatic	Use after the automatic program "Healthy Day – automatic" on the DeVita Ritm device (first week) and after the automatic program "Health - automatic" (second week).	
5	DeVita AP	Automatic programs	Without cold - automatic	Alternate every second day during the first week. Use at night.	
6			Healthy breathing - automatic		
7			Restoration 2 - automatic	Use after the automatic program "Healthy Day – automatic" on the DeVita Ritm device (first week).	
8			Mode 2 - automatic	Alternate every second day during the second week. Use at night.	
9			Healthy breathing - automatic		
10	Manual program	Without cough, Without pneumonia	If necessary, use after the automatic program on the DeVita Ritm device "Health - automatic" (second week).		

SCHEME № 3 Maintaining health after infection (2 weeks from start):

SCHEMA №3	DEVICE		PROGRAM		APPLICATION
	1	DeVita Energy		Energy Pro, Health	Once in the morning.
2	DeVita Ritm	Automatic programs	Healthy Day - automatic	Alternate every second day for 2 weeks. Use in the morning.	
3			Health - automatic		
4			Anti cold - automatic	Use after the automatic program "Healthy Day – automatic" on the DeVita Ritm device.	
5	DeVita AP	Special programs	Healthy heart	Use after the program "Anti cold - automatic" on the DeVita Ritm device.	
6			Healthy kidneys	Use after the automatic program on the DeVita Ritm device "Health - automatic".	
7	DeVita AP	Automatic programs	Restoration 2 - automatic	Alternate every second day for 2 weeks. Use at night.	
8			Mode 4 - automatic		

## DISEASE PREVENTION RECOMMENDATIONS



Avoid visiting public places: shopping centers, sports and entertainment events, public transport during rush hour.



Avoid close contact and staying in the same room with people who have visible signs of SARS (coughing, sneezing, nasal discharge).



Wash your hands with soap and water thoroughly after returning home, or after contact with strangers. Do not touch your eyes, nose, mouth with dirty hands.



Sterilize gadgets, office equipment and the surfaces you touch.



Limit tight hugs and handshakes when greeting.



Use only personal hygiene items (towel, toothbrush).



Regularly carry out wet cleaning with disinfectants and frequent ventilation of the living spaces.



You will get the most effective result by combining preventive measures with the use of DeVita-technologies and DeLixir cell nutritionstrition

Before  
use



Dry cough, fever up to 39°C and weakness throughout the day.

After the first use  
of DeVita



Improvement of well-being, temperature drop and burst of energy.

After a month of  
using DeVita



Increase of adaptive reserves of the organism, disease reduction, vigor.



Use the recovery scheme of the "Stop Coronavirus- Prevention" Complex with the smart cell nutrition supplements DeLixir MULTI-ENERGY, DeLixir pH BALANCE and DeLixir TOP FORM for much faster results.



## COMPLEX STOP INFECTION


### HOW DOES INFECTION OCCUR?



Infectious diseases can equally affect both adults and children. Some progress mildly, while others involve dysfunction of specific organs. Among these infections, there is a group considered “childhood” infections because they are more commonly observed in childhood, and after recovery, they leave a strong immunity.



Viruses enter the human body in various ways and remain “silent” for a certain time. The incubation period can last from several hours to several months. Severe viral infections suppress the immune system, increasing the likelihood of bacterial or fungal flora joining, which further exacerbates the individual’s condition.

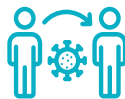

Bacterial infections enter the body through contact with an infected person or animal, consumption of contaminated food and water, skin and mucous membrane damage, and insect bites. After penetration, the symptoms of the disease do not appear immediately. Some people may not become ill but can be carriers for a long time.



Bacterial and viral infections are often confused with each other. It is very important to correctly identify the type of infection at an early stage to choose the appropriate treatment method and necessary protection. Since any of these infections weaken the body, it is necessary to choose the right course of action to restore energy and immunity.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Special DeVita AP programs are effective against most viral and bacterial agents. They not only help effectively fight them but also reduce inflammation, pain, and other undesirable manifestations of corresponding infections.

The “Mode 4” automatic program prevents the likelihood of fungal flora joining. The “Cleansing” automatic program is used to reduce the toxic load on the body after experiencing an infection.



**DeVita Ritm** – here the action is based on the studies of the scientists Dr. Reinhold Voll and A.F. Cramer, which found that each organ has a healthy frequency of oscillation. DeVita Ritm gently affects the body, pointing out the frequency of the vibrations of healthy organs. As a result, the body is capable of self-healing.

The automatic programs “Healthy Day - automatic” and “Health - automatic” have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. They strengthen the defenses, increase vitality. The automatic program “Complete Cleansing” helps in main excretory system restoration and harmonizes the body.

The use of special programs normalizes blood supply and nutrition to the affected organs, triggers regeneration processes, improving their performance.



**DeVita Energy** harmonizes the body’s work, energizes, activates the removal of toxins, activates the processes of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet.

Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- at the first sign of infection;
- to prevent possible contamination;
- to recover from an infection.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir Multi-Energy, DeLixir pH and DeLixir Detox smart cell nutrition supplements.



# APPLICATION SCHEMES

SCHEME №1: At the first signs of infection (use for 2 weeks, if necessary, extend up to 1 month).

SCHEME №1	DEVICE		PROGRAM		APPLICATION
	1	DeVita Energy	Mode program	Energy Pro, Health	Once in the morning.
2	DeVita Ritm	Auto-programs	Healthy Day*, Health*	During the day after DeVita Energy programs or from 8:00 a.m. Alternating with each other every other day.	
3		Special programs	One of the programs*: Free breath, Throat Without Problems, Intestine without problems, Men's/Women's health (depending on the problem)	After the DeVita Ritm programs "Health"	
4	DeVita AP	Auto-programs	Restoration 2	At night.	
5		Special programs	One of the manual programs (depending on the problem) + Without toxins	After the DeVita Ritm program "Healthy Day"	

SCHEME №2: Recovery after infection (to be used within one month)

SCHEME №2		DEVICE		PROGRAM		APPLICATION
		1-2 weeks	1	DeVita Energy	Mode program	Energy Pro, Health
2	DeVita Ritm		Auto-programs	Healthy Day*, Complete Cleansing*	In the daytime from 8:00. Alternating with each other every other day.	
3			Special programs	Recovery after infection**	After the DeVita Ritm programs "Healthy Day"	
4	DeVita AP		Auto-programs	Mode 4 automatic	Daily at night during the 2nd week.	
3-4 weeks	1	DeVita Energy	Mode program	Energy Pro, Health	Once in the morning.	
	2	DeVita Ritm	Auto-programs	Healthy Day*, Active recovery - automatic**	In the daytime from 8:00. Alternating with each other every other day.	
	3		Special programs	One of the programs*: Intestine without problems, Men's/ Women's health, Healthy joints, Without allergens (depending on the problem).	After the DeVita Ritm programs "Healthy Day"	
	4	DeVita AP	Auto-programs	Cleansing*	At night during the 2nd week	

\*Programs from the "General list of programs 1" complex.

\*\* Program from the "First aid" complex

## MEASURES FOR THE PREVENTION OF INFECTIOUS DISEASES



Limit contact with sick people;



Observe personal hygiene;



Perform wet cleaning of your living space.



Fight disease carriers (flies, mosquitoes, fleas, etc.).



Timely perform prophylactic vaccinations.



Take measures to strengthen the immune system (exercise, intake of vitamins, outdoor walks).



Eat a balanced diet, using natural fresh products.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before use



Elevated body temperature, weakness, decreased appetite.

After 2 days of using DeVita



Decrease in temperature, improvement in well-being.

After a month of using DeVita



Cheerfulness, activity, good mood.



Use the scheme of the "STOP INFECTION" complex in conjunction with DeLixir MULTI-ENERGY, DeLixir pH BALANCE and DeLixir Detox smart cell nutrition supplements, and your results will be noticeable much faster.



# SPECIAL COMPLEXES

(ARE USED ONLY AFTER ONE OF THE BASIC COMPLEXES)





## COMPLEX BODY FIGURE

### WHY IS EXCESS WEIGHT DANGEROUS?



According to WHO, more than 3 billion people on Earth are already overweight, and this number is constantly increasing. Children also begin to suffer from problems with body weight - about 40 million children under the age of 15 weigh much more than the norm. Doctors sound the alarm - being overweight leads to many health problems.



With an increase in body weight, the load on the musculoskeletal system increases. Adipose tissue forms between the internal organs, disrupting their work. Vascular problems and shortness of breath appear. Due to the increased load on the legs, varicose veins develop. All this reduces the immune defense, and overweight people are more likely to get sick.



The main factor in gaining excess weight is unhealthy diet and lack of activity. This is due to the accelerating rhythm of life. We travel by car and public transport, walk a little, have a bite on the run. After work, we feel tired, and we have absolutely no strength and desire to play sports.



For most overweight people, losing weight is a very difficult task. Weight loss is a complex process that includes changing the diet, including physical activity in the daily regimen and the use of special supplements or drinks that help balance the metabolism.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasitic forms that live in the human body cause malfunctions in the work of organs and, as a result, metabolic disorders. And this is one of the main factors in weight gain. DeVita AP antiparasitic programs affect most parasites that negatively affect metabolism, destroying them and effectively removing parasitic load.

“Weight Reduction” automatic program reduces the parasitic load and destroys viruses that cause excess weight. In addition, it normalizes the functioning of the autonomic system, maintains the tone of the body and improves the psychoemotional state.

“Reduce Points” is used to cleanse the body of toxins, normalize lymphatic drainage, maintain immunity and increase energy.

**DeVita Ritm** gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“High Tonus” automatic program is necessary to relieve toxic load from the liver, lymphatic system and kidneys, activate their work, and improve blood circulation. This helps to increase energy and improve the tone of the body. In addition, the automatic program launches natural rejuvenation processes.

“Healthy Metabolism” automatic program harmonizes the work of the body, regulates the activity of the autonomic system, stabilizes the appetite, and activates the digestive and excretory organs. This normalizes metabolic processes and the psychoemotional state, helps to reduce body weight and improve the quality of life.

**DeVita Energy** activates self-regulation mechanisms, helps accelerate the recovery process and the distribution of energy balance. Allows the body to maintain shape in the process of cleansing and losing weight, maintain vitality and good mood.



## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the effects of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water for at least 1-1.5 l / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- to reduce your weight;
- to improve the figure for the summer period or an important event;
- for athletes: to “dry” the body before the competition.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir ULTRA SLIM+ smart cell nutrition supplement.

SCHEME	DEVICE		PROGRAM		USAGE
	1	DeVita Ritm	Auto-program	High Tonus- automatic	1 <sup>st</sup> day starting at 15:00
2	Auto-program		Healthy Metabolism- automatic	2 <sup>nd</sup> day starting at 15:00	
3	DeVita AP	Auto-program	Weight Reduction- automatic	1 <sup>st</sup> day after the usage of DeVita Ritm or at night	
4		Auto-program	Reduce points- automatic	2 <sup>nd</sup> day after the usage of DeVita Ritm or at night	



1. It is recommended to work on this scheme after the general cleansing of the body using the complex "Digestion without problems" or "Stop parasites".
2. Further work on the scheme involves the use of the DeVita Energy device before the "High tonus" and "Healthy Metabolism" programs: "Energy Pro" and "Slimness" in the morning, followed by DeVita Ritm automatic programs: "Healthy Day" and "Health" (see scheme No. 1 "Digestion without problems")

Attention: when using this program scheme, increased sweating and urination are possible.

# PREVENTION MEASURES FOR INCREASED BODY WEIGHT



Eat 5 times a day small portions.



The last meal should be 2 hours before bedtime.



Drink plenty of water. Water reduces appetite and helps to involve deferred fats in the metabolism, as well as remove toxins.



Eat fruits, they allow you to effectively lose weight.



Vitamins also accelerate the metabolism in cells.



Anti-cellulite massage of problematic areas accelerates blood circulation, the influx to the site of exposure increases, which significantly speeds up the metabolism.



Exercise and physical activity allow you to effectively lose weight and adjust your figure.



Observe your sleep patterns. This contributes to the renewal of brain cells, which, in turn, leads to an acceleration of the metabolism.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before  
use



Increased weight, bad mood,  
heaviness and stiffness in  
movements.

After 1st week of  
using DeVita



Significant improvement  
in immunity, increased overall  
body tone.

After a month of  
using DeVita



Well-being, good mood,  
lightness and ease of  
movement.



Use the scheme of the "BODY FIGURE" complex together with DeLixir ULTRA SLIM+ smart cell nutrition supplement, and your results will be noticeable much faster.





# COMPLEX HEALTHY CELL

## WHY CARCINOGENS ARE DANGEROUS?



Carcinogens are certain factors under the influence of which a person is more likely to develop tumor processes. Carcinogens include organic (food), inorganic (household chemicals) substances, and ionizing radiation. The development of the process is influenced by bad habits and the ecology of the area.



Smoking is the most significant cause of cancer. Nutrition is another reason for the development of pathological processes, especially in the digestive tract. Civilization costs lead to the fact that food is contaminated with harmful environmental substances and chemicals of agriculture, including carcinogens.



Moderate doses of solar radiation have a beneficial effect on the human body. But the desire to tan better leads to overheating, skin burns and thermal shock. Such abuse of sunbathing can lead to serious disorders in the body.



Infectious agents (some viruses, bacteria, fungi) with prolonged exposure to the body and certain conditions are also risk factors for the development of such diseases. This was proven back in 1910 and confirmed in the 1980s.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Prolonged presence of carcinogenic factors in the body: biological, chemical, physical, can lead to severe intoxication and uncontrolled division and proliferation of atypical cells in any human organs. The complex antiparasitic programs on the DeVita AP device help in toxic load and inflammatory reactions reduction, neutralize infectious pathogens.



“Cleansing” automatic program is used to remove toxins and reduce the toxic load on the excretory system after applying the complex of antiparasitic programs.

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

“Healthy Day” and “Health” automatic programs have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. Strengthens defenses, increases vitality. The use of special programs helps to reduce the load on specific organs and systems, restores their energy potential. The automatic program “Complete Cleansing” helps in main excretory system restoration and harmonizes the body.



**DeVita Energy** harmonizes the body's work, energizes, activates the removal of toxins, activates the processes of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- if you have a bad diet;
- if you have bad habits;
- if you get sick often;
- if you live in an unfavourable ecological region;
- if you are under constant stress.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy Pro, Protection		Once in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
	5			Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
	6			Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
	7			Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy Pro, Protection		Once in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day
	3		Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night
	5		Special programs	No schistosomes (blood-flukes) & trematods	Daily in the evening after DeVita Ritm programs for 1 <sup>st</sup> -2 <sup>nd</sup> weeks
6	Drainage Antitox			Daily in the evening after DeVita Ritm programs during 3 <sup>rd</sup> -4 <sup>th</sup> weeks	

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy Pro, Protection		Once in the morning
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Alternate every other day. Turn on at 8:00a.m.
	3		Special programs	Organ Correction Programs	After the DeVita Ritm program "Healthy Day"
	4	DeVita AP	Auto-programs	Cell Protection 1 - automatic Cell Protection 2 - automatic	Alternate every other day within 2 weeks. Turn on at 5:00p.m.
5	Cell Protection 3 - automatic Cell Protection 4 - automatic			Alternate every other day within 2 weeks. Turn on at 5:00p.m.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Ritm	Auto-programs	Healthy Day - automatic	1 <sup>st</sup> day at 8:00 a.m.
				Complete Cleansing - automatic	2 <sup>nd</sup> day at 8:00a.m.
	2		Special programs	Organ Correction Programs	After the DeVita Ritm program "Healthy Day"
3	DeVita AP	Auto-programs	Cleansing - automatic	1 <sup>st</sup> day at night	
4		Special programs	Without Papilloma - General	2 <sup>nd</sup> day in the evening and at night	

Programs of 3<sup>rd</sup> and 4<sup>th</sup> months to repeat for another 2 cycles!

## PREVENTIVE MEASURES FOR NEOPLASMS DEVELOPMENT



Exclude sweets, white bread, convenience foods, margarine, smoked meats from ration.



Give up bad habits (alcohol, tobacco).



Timely treat inflammatory infections.



When staying in the open sun, use sunscreen, a hat, clothes that cover open areas of the body as much as possible.



Move more. This contributes to the physiological course of all processes in the body.



Strengthen immunity with regular walks in the fresh air and normalization of the psycho-emotional state.



Have regular physical examinations.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Weakness, fatigue,  
decreased immunity.

After a week of  
using DeVita



Increased  
vitality.

After a month of  
using DeVita



Strengthened immunity,  
improved performance.



Use the scheme of the "HEALTHY CELL" complex together with DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements, and your results will become noticeable much faster.



## COMPLEX LYME DISEASE

### WHAT DANGERS DO TICK BITES CARRY?



Ticks are one of the most diverse and ancient groups of arthropods on Earth. As a rule, mites are free-living, but parasitic blood-sucking species are also found. Among them, the most famous are the ixodid ticks (Ixodoidea). They carry the causative agents of human diseases: tick-borne encephalitis, tick-borne borreliosis (Lyme disease), etc.



Ticks are active in the period: late March, early April, increasing numbers by mid-May. The maximum number is preserved until the end of June, and sometimes until September. Ticks prefer moderately shaded and moistened deciduous and mixed forests. Ticks lie in wait for their prey, sitting on the ends of blades of grass, sticks and twigs sticking up.



Tick saliva is able to dissolve the walls of blood vessels and reduce immunity. Infectious agents are also transmitted through saliva, which pose a danger to human health. Agent transfer occurs in the first minutes after being bitten. Tick-borne encephalitis is the most dangerous of tick-borne infections. Tick-borne borreliosis is the most common.



Infectious agents caught by a tick bite infect the nervous system, musculoskeletal system, and the cardiovascular system. They lead to intoxication and further damage to the liver and kidneys, a general decrease in immunity and a state of constant fatigue.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Parasites that enter the body cause severe intoxication with damage to the nervous, cardiovascular systems and musculoskeletal system. The DeVita AP complex of antiparasitic programs acts on specific pathogens and associated parasites, quickly, effectively and without side effects, removing pathogenic toxic load and signs of inflammation.



“Cleansing” automatic program is used to remove toxins and reduce the toxic load on the excretory system after applying a complex of antiparasitic programs.

**DeVita Ritm** – In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

The automatic programs “Healthy Day” and “Health” have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. Strengthen defenses, increase vitality. Special programs “Nervous system- general”, “Healthy joints” and “Healthy heart” help reduce the load on these organs and restore their energy potential. “Complete Cleansing” automatic program helps in main excretory system restoration and harmonizes the body.



**DeVita Energy** harmonizes the body's work, energizes, activates the processes of toxin elimination. Improves the general condition of intoxication resulting from infection after a bite. Reduces lethargy and drowsiness.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- immediately after a tick bite;
- if the tick bite was earlier;
- with problems with the cardiovascular system, joints, kidneys, liver, resulting from a tick bite.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy Pro, Health		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
			Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
	4	DeVita AP	Auto-programs	Mode 1 automatic	Daily at night during the 1 <sup>st</sup> week.
				Mode 2 automatic	Daily at night during the 2 <sup>nd</sup> week.
				Mode 3 automatic	Daily at night during the 3 <sup>rd</sup> week.
				Mode 4 automatic	Daily at night during the 4 <sup>th</sup> week.

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy Pro, Health		Once in the morning.
	2	DeVita Ritm	Auto-programs	Healthy Day automatic, Health-automatic	Daily in the morning, alternating with each other every other day.
			Special programs	Digestion normalization and / or Intestine without problems	Daily in the afternoon, if needed.
	4	DeVita AP	Auto-programs	Cleansing-automatic	Daily at night.
				Special programs	No schistosomes (blood-flukes) & trematods
Special programs			Drainage Antitox		Daily in the evening after DeVita Ritm programs for 3-4 weeks.

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy Pro, Health		Once in the morning.
	2	DeVita Ritm	Auto-programs	Health Day automatic, Health - automatic	Alternate every other day. Turn on time 8:00.
			Special programs	Nervous system - General, Healthy heart	After DeVita Ritm "Healthy Day" program, alternating every other day.
	4	DeVita AP	Auto-programs	Lyme 1 - automatic, Lyme 2 - automatic	Alternate every other day for 2 weeks. Turn on time 17:00.
Lyme 3 - automatic, Lyme 4 - automatic				Alternate every other day for 2 weeks. Turn on time 17:00.	

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Ritm	Auto-programs	Healthy Day - automatic	1 <sup>st</sup> day from 8:00 a.m.
	2		Complete Cleansing - automatic	2 <sup>nd</sup> day from 8:00 a.m.	
	3	Special programs	Nervous system - General, Healthy joints		After DeVita Ritm's "Healthy Day" program, alternating every other day.
	4	DeVita AP	Auto-programs	Cleansing-automatic	1 <sup>st</sup> day at night.
Nerve Regeneration - automatic				2 <sup>nd</sup> day at night.	

Programs of 3rd and 4th months are to repeat for 2 more cycles!

# PREVENTIVE MEASURES FOR TICK BITES



When in tick habitats, avoid dark tones in clothing, as ticks are harder to spot on a dark background.



Pull outerwear into your pants and pants into your socks. If there is no hood, put on a hat.



Every 15 minutes, inspect your clothes, and if possible, inspect the head and body, especially above the waist, ticks are most often found there.



Use chemical protection against ticks. Treat exposed body parts, clothing, sleeping bag, tent and other items.



If you most often spend time in a summer cottage, where ticks live, then this territory should be treated with a special insecticidal agent that kills ticks.



Arriving home, remove clothing outside the living room and carefully inspect it, paying particular attention to folds, seams, pockets.



If you find a sucking tick, do not remove it yourself. Contact the regional clinic or district emergency room.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Headache, weakness,  
lack of appetite.

After a week of  
using DeVita



Reduced pain, increased  
overall tone.

After a month of  
using DeVita



Normalization of the functional  
and psychological state.



Use the scheme of the "LYME DISEASE" complex together with DeLixir pH BALANCE, DeLixir DETOX and DeLixir MULTI-ENERGY smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX HEALTHY TEETH

### WHAT ARE THE DANGERS OF ORAL CAVITY PROBLEMS?



The human mouth is the main gateway to the body. Through it, not only vital substances enter the body, but also pathogenic forms. Therefore, it is necessary to respond in a timely manner to various changes and strange symptoms, because they can indicate not only problems with teeth, but also other organs.



The most common problems are bad breath, tooth decay, sores on the mucous membrane, plaque on the tongue.

These manifestations may indicate the presence of viral (herpes), bacterial (streptococcus) and fungal (candida) infections.



Violation of the blood microcirculation, the metabolism can lead to bleeding and inflammation of the gums, which may cause loosening and loss of teeth. The same reasons often cause inflammation of the lips and salivary glands.



The condition of the oral cavity (teeth, mucous membrane, gums, tongue) is an indicator of the work of many internal organs. It is affected by: prolonged use of medications, immune failure, inflammatory processes, unbalanced nutrition, bad habits, vitamin deficiency and other factors.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Failure to comply with hygiene rules, bad habits, the presence of chronic inflammatory processes in the body creates a favorable environment for the development of pathogenic organisms in the oral cavity, and as a result, the occurrence of problems with teeth, gums and other organs of the oral cavity. DeVita AP antiparasitic programs help reduce pathogenic infectious load, while Special programs reduce inflammatory reactions and help recover faster after dental interventions.

“Cleansing” automatic program is used to remove toxins and reduce the toxic load on the excretory system after applying the complex of antiparasitic programs.



**DeVita Ritm** –In this device, the development of scientists R. Voll and A. F. Kramer is embodied, which found that each organ has a healthy frequency of oscillations. DeVita Ritm gently affects the body, by giving it the frequency of vibrations of healthy organs. As a result, the body works on self-healing.

The automatic programs “Healthy Day” and “Health” have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. Strengthen defenses, increase vitality. The use of special programs improves blood microcirculation, strengthens gums, reduces bleeding, and restores mucous membranes. The automatic program “Complete Cleansing” normalizes the functioning of the main excretory systems and the coordination of the work of the whole organism.

**DeVita Energy** harmonizes the body’s work, energizes, activates the removal of toxins, starts the processes of self-regulation.



## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the impact of the devices. Remove or minimize coffee in your diet. Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water, at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- if you have a bad diet;
- you have bad habits;
- you are prone to frequent colds;
- you have bleeding gums, loose teeth, damage to the oral mucosa;
- you have planned dental procedures (treatment, removal, implantation, etc.).



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir MULTI-ENERGY, DeLixir pH BALANCE smart cell nutrition supplements.

Scheme 1. Use after dental interventions (removal, implantation, treatment, etc.). Duration of use 2 weeks:

SCHEME №1 DURING 1 <sup>st</sup> MONTH	DEVICE	PROGRAM		USAGE	
	1	DeVita Energy	Energy Pro, Health		Once in the morning.
2	DeVita Ritm	1 <sup>st</sup> day	Auto-program	Healthy Day - automatic	1 <sup>st</sup> day. Turn on at 8:00 a.m.
			Special program	Recovery from injuries	
3		2 <sup>nd</sup> day	Auto-program	Health - automatic	2 <sup>nd</sup> day. Turn on at 8:00 a.m.
			Special program	Strong teeth	
4	DeVita AP	Auto-programs	Restoration 1 automatic		1 <sup>st</sup> day in the evening after DeVita Ritm programs.
Mode 2 automatic			1 <sup>st</sup> day. At night.		
Restoration 2 automatic			2 <sup>nd</sup> day. At night.		

In the presence of abscesses, fistulas add the program "Without fistula" daily.

Scheme 2. Use in chronic problems of the oral cavity:

The first month according to the scheme of 1 month of the complex "Digestion without problems".

SCHEME №2 DURING 2 <sup>nd</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Energy	Energy Pro, Health	
2	DeVita Ritm	Auto-programs	Healthy Day - automatic, Health - automatic	In the daytime from 8:00. Alternating every other day.
		Special programs	Strong teeth	
3			After the DeVita Ritm program "Healthy Day".	
4	DeVita AP	Auto-programs	Healthy Gums - automatic Cleansing - automatic	At night. Alternating every other day.

SCHEME №3 DURING 3 <sup>rd</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Ritm	Auto-programs	Healthy Day - automatic, Complete Cleansing - automatic
Special programs			Strong teeth	
2			After the DeVita Ritm program "Healthy Day".	
3	DeVita AP	Auto-programs	Mode 1-4 automatic	If needed, at night.

SCHEME №4 DURING 4 <sup>th</sup> MONTH	DEVICE	PROGRAM		USAGE
	1	DeVita Ritm	Auto-programs	Healthy Day - automatic, Health - automatic
Special programs			Strong teeth	
2			After the DeVita Ritm program "Healthy Day".	
3	DeVita AP	Auto-programs	No parasites automatic	At night, daily for 1 week.

# PREVENTIVE MEASURES FOR DENTAL PROBLEMS



Brush your teeth twice a day: morning and evening.



Use a damp brush to clean the gums, tongue, and inside of the cheeks..



Timely treat inflammatory infections of the ENT and other internal organs.



Use dental floss to clean the interdental spaces.



Give up bad habits.



Strengthen immunity, normalize the psycho-emotional state, eat right.



Check with your dentist regularly.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements.

Before  
use



Soreness and bleeding  
gums, decreased appetite.

After the first use  
of DeVita



Reduced soreness, improved  
overall well-being.


After a month of  
using DeVita



Reduced bleeding gums,  
strengthened immunity.



Use the scheme of the "HEALTHY TEETH" complex together with the DeLixir pH BALANCE and DeLixir MULTI-ENERGY smart cell nutrition supplements, and your results will be noticeable much faster.

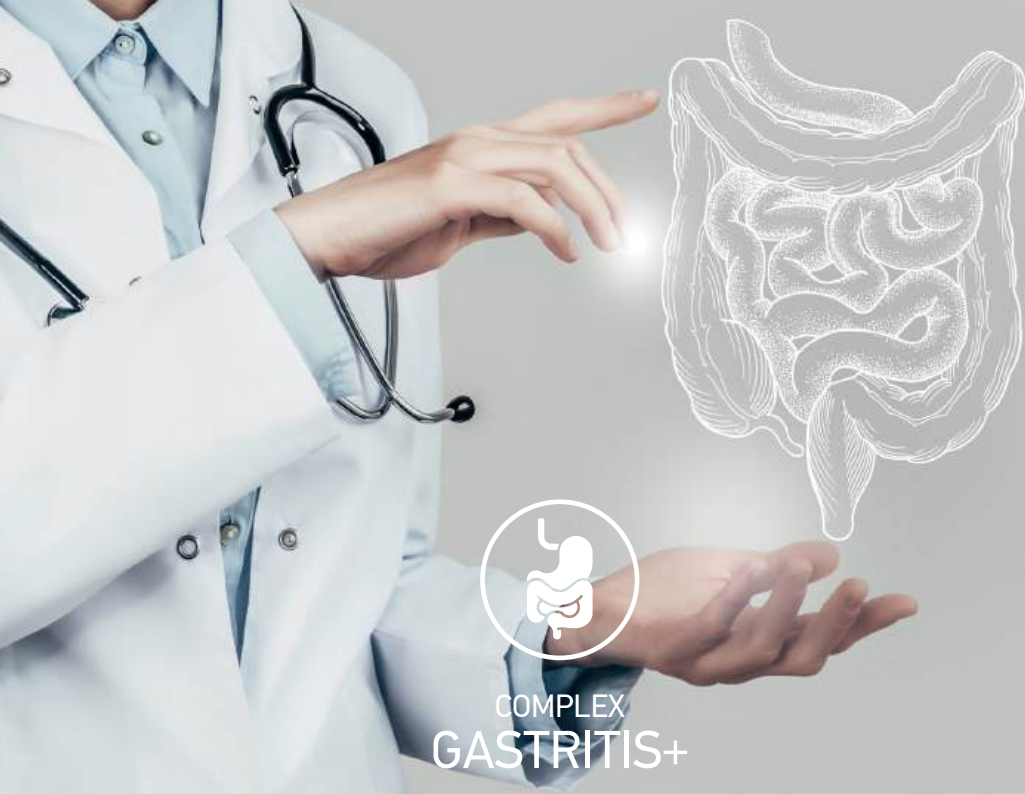

## COMPLEX GASTRITIS+

### WHAT DANGERS DO DIGESTIVE SYSTEM PROBLEMS CAUSE?



The digestive system is one of the most important systems in the body. In it, there is a transformation of nutrients into the necessary building materials for the tissues of the body. If there is a failure at any stage of digestion, then the well-being and quality of life are noticeably impaired.



The most common problems that people see doctors for are stomach pain, heartburn, and belching. According to statistics, 60% of the world's population suffers from gastritis. Harmful habits, improper nutrition, stress, the presence of specific bacteria, and some hereditary factors contribute to this.



Problems with one digestive organ do not pass without a trace for others. For example, stomach diseases can provoke problems with the duodenum. Liver diseases - problems with the pancreas, etc. Therefore, it is important to solve not the problem of one organ of the gastrointestinal tract, but to improve the work of the digestive system as a whole.



A serious problem for the world's population is infection with trematodes (flukes). After pinworms and roundworms, opisthorchis are the third most common infection. These parasites affect the bile ducts, gallbladder, and ducts of the pancreas. They can live in the body for up to 25 years.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Non-compliance with hygiene rules, harmful habits, poor diet, frequent stress, create a favorable environment for the formation of inflammatory processes in the organs of the digestive system. The DeVita AP programs contribute to reducing the pathogenic load, while special programs reduce inflammatory reactions.

“Cleansing” automatic program is used to remove toxins and reduce the parasitic load on the excretory system after applying the complex of antiparasitic programs.



**DeVita Ritm** – here the action is based on the studies of the scientists Dr. Reinhold Voll and A.F. Cramer, which found that each organ has a healthy frequency of oscillation. DeVita Ritm gently affects the body, pointing out the frequency of the vibrations of healthy organs. As a result, the body is capable of self-healing.

The automatic programs “Healthy Day-automatic” and “Health-automatic” have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. They strengthen the defenses, increase vitality. The automatic program “Complete Cleansing” helps restoring and harmonizing the body’s functioning.

Using special programs contributes to improving blood microcirculation, restoring mucous membranes, and normal functioning of the digestive organs.



**DeVita Energy** harmonizes the body’s work, energizes, activates the removal of toxins, activates the processes of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the effectiveness of the devices. Remove or minimize coffee in your diet.

Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- If you have an unhealthy diet;
- you have harmful habits;
- you experience frequent abdominal pain;
- you have heartburn, acid reflux, bloating;
- you suffer from constipation.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir pH and DeLixir Detox smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME NO.1 FIRST MONTH	According to the scheme of 1st month of the "Digestion without problems" complex
----------------------------	--

SCHEME NO.2 SECOND MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Energy	Basic	Energy or Energy Pro	Once in the morning
	2	DeVita Ritm	Automatic programs	Healthy day*, Health*	In the daytime after DeVita Energy programs or from 8:00. Alternating every other day.
	3		Special programs	One of the programs: "Stomach+automatic" or "Healthy intestine automatic" or "Normal blood sugar" or "Healthy pancreas" (depending on a problem)	After DeVita Ritm "Healthy Day" program
	4	DeVita AP	Automatic programs	One of the 3 automatic programs: "Active gastro" or "Bowel without irritation" or "Without constipation" (depending on a problem), Cleansing*	In the nighttime. Alternating one of the 3 automatic programs with the "Cleansing" program every other day
	5		Special programs	"Without fluke", "Without toxins"	If needed. After DeVita Ritm "Health" program

SCHEME NO.3 THIRD MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Ritm	Automatic programs	"Healthy day"* , "Complete cleansing"*	In the daytime from 8:00. Alternating every other day
	2		Special programs	Without disbacteriosis	After DeVita Ritm "Healthy day" program
3	DeVita AP	Automatic programs	Stop parasites*	In the nighttime, daily for 1 week	

SCHEME NO.4 FOURTH MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Ritm	Automatic programs	Healthy day*, Health*	In the daytime from 8:00. Alternating every other day.
	2		Special programs	One of the programs: "Stomach+automatic" or "Healthy intestine automatic"	After DeVita Ritm "Healthy day" program.
3	DeVita AP	Automatic programs	"Healthy pancreas", "Blood sugar under control", "Stop diabetes"	At night for 1 month, if needed	

\* Programs of "General list of complexes 1"

## PREVENTING STOMACH AND INTESTINAL DISEASES



Eat healthily, maintain a regular eating schedule.



Quit smoking, avoid alcohol consumption.



Adhere to hygiene rules.



Avoid stress, think positively.



Timely treat gastrointestinal diseases.



People suffering from gastritis need to adhere to a gentle diet, excluding irritating foods such as coffee, chocolate, spices, smoked and fried foods. Eat fractionally 5-6 times a day in small portions.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before use



Pain, heartburn, bloating, flatulence, and bowel movement problems.

After first use  
of DeVita



Reduction of pain, relief of heartburn, heaviness, and bloating.

After a month of  
using DeVita



Normalization of digestion, elimination of exacerbations of chronic diseases, excellent well-being.



Use the scheme of the "Gastritis +" complex in conjunction with DeLixir pH BALANCE and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster.





## COMPLEX HEALTHY BACK

### WHAT ARE THE RISKS OF BACK PAIN?



In recent times, an increasing number of people are experiencing back pain. Most often, this is associated with a sedentary lifestyle, endless stress at work and home, and internal disturbances in the body. In addition, there is a tendency for this problem to “rejuvenate”, i.e. more and more young people suffer from this ailment.



The back is the support of our body, it consists of the spine, ribs, pelvic bones and large and small muscles located on them. This complex “structure” gives a person the ability to move and protects internal organs from damage. The general condition of the whole organism, the quality and duration of life depend on the health of the back.



The spine is the frame of our body. Its health is the key to excellent health and the prevention of many diseases. Disruption of the functionality of the spine leads to diseases of the spine itself and to chronic diseases of the internal organs and extremities.



Muscles are the active part of the musculoskeletal system, contracting, they act on the bone levers, setting them in motion. A decrease in muscle activity occurs due to fatigue or spasms, which can occur due to overcooling, muscle overload, injury, inflammation, or intoxication.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

A sedentary lifestyle, harmful habits, poor diet, frequent stress, create a favorable environment for the formation of back pain. “DeVita AP programs contribute to reducing inflammatory reactions, while special programs reduce pain sensations.”

The “Cleansing” automatic program is used to reduce the toxic load on the tissues and organs of the musculoskeletal system, improving mobility and overall well-being.



**DeVita Ritm** – here the action is based on the studies of the scientists Dr. Reinhold Voll and A.F. Cramer, which found that each organ has a healthy frequency of oscillation. DeVita Ritm gently affects the body, pointing out the frequency of the vibrations of healthy organs. As a result, the body is capable of self-healing.

The automatic programs “Healthy Day-automatic” and “Health-automatic” have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. They strengthen the defenses, increase vitality. The automatic program “Complete Cleansing” helps restoring and harmonizing the body’s functioning.

Using special programs contributes to improving blood microcirculation, restoring muscle and connective tissue of the musculoskeletal system, and normalizing the functioning of muscles and the spine.



**DeVita Energy** harmonizes the body’s work, energizes, activates the removal of toxins, activates the processes of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the effectiveness of the devices. Remove or minimize coffee in your diet.

Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- when experiencing pain and fatigue in the back area;
- for injuries and inflammation of muscles and the spine;
- Posture disorders
- If your work involves little mobility.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir Multi-Energy, DeLixir Collagen+ and DeLixir Detox smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME NO.1 FIRST MONTH	According to the scheme of 1st month of the "Movement without pain" complex
----------------------------	---

SCHEME NO.2 SECOND MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Energy	Basic	Energy or Energy Pro	Once in the morning
	2	DeVita Ritm	Automatic programs	"Healthy day"*, "Health"*	In the daytime after DeVita Energy programs or from 8:00. Alternating every other day
	3		Special programs	Healthy spine automatic	After DeVita Ritm "Healthy day" program
	4	DeVita AP	Automatic programs	One of 2 automatic programs: Low back without pain, Healthy spine (depending on the problem), Cleansing*	In the nighttime. Alternate one of 2 automatic programs along with Cleansing program every other day.

SCHEME NO.3 THIRD MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Ritm	Automatic programs	"Healthy day"*, "Complete cleansing"*	In the daytime from 8:00. Alternating every other day
	2		Special programs	Healthy spine automatic	In the nighttime, daily.
	3	DeVita AP	Special programs	"Stop bacteria"*, "Drainage antitox"*	After DeVita Ritm "Healthy day" program

SCHEME NO.4 FOURTH MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Ritm	Automatic programs	"Healthy day"*, "Health"*	In the daytime from 8:00. Alternating every other day.
	2		Special programs	Healthy spine automatic	After DeVita Ritm "Healthy day" program
	3	DeVita AP	Automatic programs	Один из 2-х автоматов: Поясница без боли, Здоровая спина (в зависимости от проблемы)	At night for 1 month.

\* Programs of "Movement without pain" complex.

# PREVENTION OF DISEASES OF THE SPINE AND BACK MUSCLES



Engage in sports or actively walk for at least half an hour a day.



Strengthen your back muscles.



Maintain good posture, properly arrange your workspace.



Lift heavy weights correctly, distribute the load evenly when carrying bags.



Choose a good mattress and pillow



Control your body weight



Eat a balanced diet. Optimal foods for the back include: eggs, dairy products, fish, liver, seeds, beans, nuts, corn, buckwheat, walnuts, flaxseed oil, cucumbers, potatoes, carrots, pumpkin, apples, bananas.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before use



Prolonged pain in the lower back, a feeling of pressure in the spine, numbness and sharp pain, inability to move freely.

After 2 days of using DeVita



Reduction of pain sensations, feelings of pressure, and numbness.

After a month of using DeVita



Absence of pain, restoration of free movement.



Use the health improvement scheme of the 'HEALTHY BACK' complex." in conjunction with DeLixir pH COLLAGEN, DeLixir DETOX and DeLixir MULTI-ENERGY, smart cell nutrition supplements your results will be noticeable much faster.



## COMPLEX HEALTHY SKIN

### WHY IS SKIN HEALTH IMPORTANT?



The skin is the outer covering of the body, actively participating in all its life processes, performing a number of important functions: respiratory, protective, secretory, receptor, and others. It is the largest organ of the human body, often one of the first to signal internal organ diseases.



With age, a person's adaptive capabilities of all organs and tissues, including the skin, deteriorate. In elderly and old age, people's skin is more prone to various dermatoses, dryness, and the appearance of spots. During this period, any damage to the skin can have a long period of regeneration and often becomes infected.



Skin damage largely depends on the phototype - the lighter the skin, the more susceptible it is to the negative effects of environmental factors (sun, cold, wind), and the faster it loses elasticity, ceases to perform protective functions, becomes thin, and develops pigmentation.



Failure to observe personal hygiene or weakening of protective functions can allow pathogenic fungi to enter the body. As a rule, this occurs through personal contact with an infected person, their belongings, or by walking on surfaces contaminated with fungi (baths, floors). Pathogens are capable of affecting not only the skin, nails, and hair but also underlying tissues and even organs.

## MECHANISM OF ACTION



**DeVita AP** – the action is based on the phenomenon of resonance, in particular, the studies of the scientist and physician Royal Raymond Rife. Rife discovered that each microorganism has an individual vibration frequency, an increase in the amplitude of which leads to the death of the microorganism without harming other cells.

Viral, bacterial, and fungal infections, harmful habits, poor nutrition and care, frequent stress are causes of skin, nail, and hair problems. The DeVita AP programs contribute to reducing inflammatory reactions, decreasing infectious and toxic load.

The “Cleansing” automatic program is used to cleanse the body after applying a complex of antiparasitic programs.



**DeVita Ritm** – here the action is based on the studies of the scientists Dr. Reinhold Voll and A.F. Cramer, which found that each organ has a healthy frequency of oscillation. DeVita Ritm gently affects the body, pointing out the frequency of the vibrations of healthy organs. As a result, the body is capable of self-healing.

The automatic programs “Healthy Day - automatic” and “Health - automatic” have a beneficial effect on the work of the whole organism, increasing its resistance to foreign agents. They strengthen the defenses, increase vitality. The “Complete Cleansing” automatic program promotes the restoration and harmonization of the body’s functioning.

Using special programs activates metabolic processes, improves blood microcirculation, alleviates inflammatory manifestations, improves the condition of the skin, nails, and hair, and slows down the aging process.



**DeVita Energy** harmonizes the body’s work, energizes, activates the removal of toxins, activates the processes of self-regulation.

## DIRECTIONS FOR USE



### NUTRITIONAL RECOMMENDATIONS FOR USING DEVITA

Do not drink carbonated alcoholic and non-alcoholic drinks - they reduce the effectiveness of the devices. Remove or minimize coffee in your diet.

Avoid eating “junk food” (chips, crackers, etc.) - this can cause digestive problems. Drink clean drinking water at least 1-1.5 L / day to remove toxins from the body.



### USE DEVITA DEVICES IN THE FOLLOWING SITUATIONS:

- for inflammatory conditions on the skin and nails;
- occurrence of pigmentation, warts, cracks;
- skin damage.



### COMBINE THE DEVICES WITH DELIXIR SMART CELL NUTRITION

For best results, it is recommended to use the DeVita devices in conjunction with the DeLixir Collagen+ and DeLixir Detox smart cell nutrition supplements.

# APPLICATION SCHEMES

SCHEME NO.1 FIRST MONTH	According to the scheme of 1st month of the "Clean skin" complex
----------------------------	--

SCHEME NO.2 SECOND MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Energy	Basic	Energy or Energy Pro	Once in the morning
	2	DeVita Ritm	Automatic programs	"Healthy day"*, "Health"*	It the daytime after DeVita Energy programs or from 8:00. Alternating every other day
	3		Special programs	"Healthy skin" or "Smoothing of wrinkles" and "Facial tone" to maintain youthfulness	
	4	DeVita AP	Automatic programs	One of the 3 automatic programs: "Skin without pimples", "Healthy nails", "Stop skin fungus"*	In the nightttime. Alternate one of 3 automatic programs along with Cleansing program every other day.
	5		Special programs	Any of the manual programs depending on the problem. "Without toxins"	

SCHEME NO.3 THIRD MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Ritm	Automatic programs	"Healthy day"*, "Complete cleansing"*	In the daytime from 8:00. Alternating every other day
	2		Special programs	"Healthy skin" or "Smoothing of wrinkles" and "Facial tone" to maintain youthfulness	
3	DeVita AP	Automatic programs	Stop parasites*	In the nighttime daily for 1 week.	

SCHEME NO.4 FOURTH MONTH	DEVICE	PROGRAM		APPLICATION	
	1	DeVita Ritm	Automatic programs	"Healthy day"*, "Health"*	In the daytime from 8:00. Alternating every other day.
	2		Special programs	"Healthy skin" or "Smoothing of wrinkles" and "Facial tone" to maintain youthfulness	
	3	DeVita AP	Automatic programs	One of the 3 automatic programs: "Skin without pimples", "Healthy nails", "Stop skin fungus"	At the nighttime for 1 month.
4	Special programs		Any of the manual programs depending on the problem. "Without toxins"		

\*Programs of "Clean skin" complex.

# RECOMMENDATIONS TO MAINTAIN SKIN HEALTH



Promptly treat acute and chronic illnesses



Eat a rational and balanced diet.



Consume an adequate amount of fluids.



Eliminate unhealthy habits.



Follow a rational work and rest schedule, get enough sleep.



Avoid aggressive exposure of the environment to the skin.



Maintain personal hygiene: take a shower daily, do not use others' hygiene products and cosmetics, use antibacterial agents less frequently.



You will get the most effective results by combining the preventive measures with the use of DeVita technologies and DeLixir smart cell nutrition supplements

Before use



Inflammation, redness, scaling, itching.

After 2 days of using DeVita



Reduction of inflammation, itching.

After a month of using DeVita




Improvement of the skin, hair, and nails condition. Good mood and well-being.



Use the scheme of the "HEALTHY SKIN" complex with the DeLixir COLLAGEN+ and DeLixir DETOX smart cell nutrition supplements, and your results will be noticeable much faster



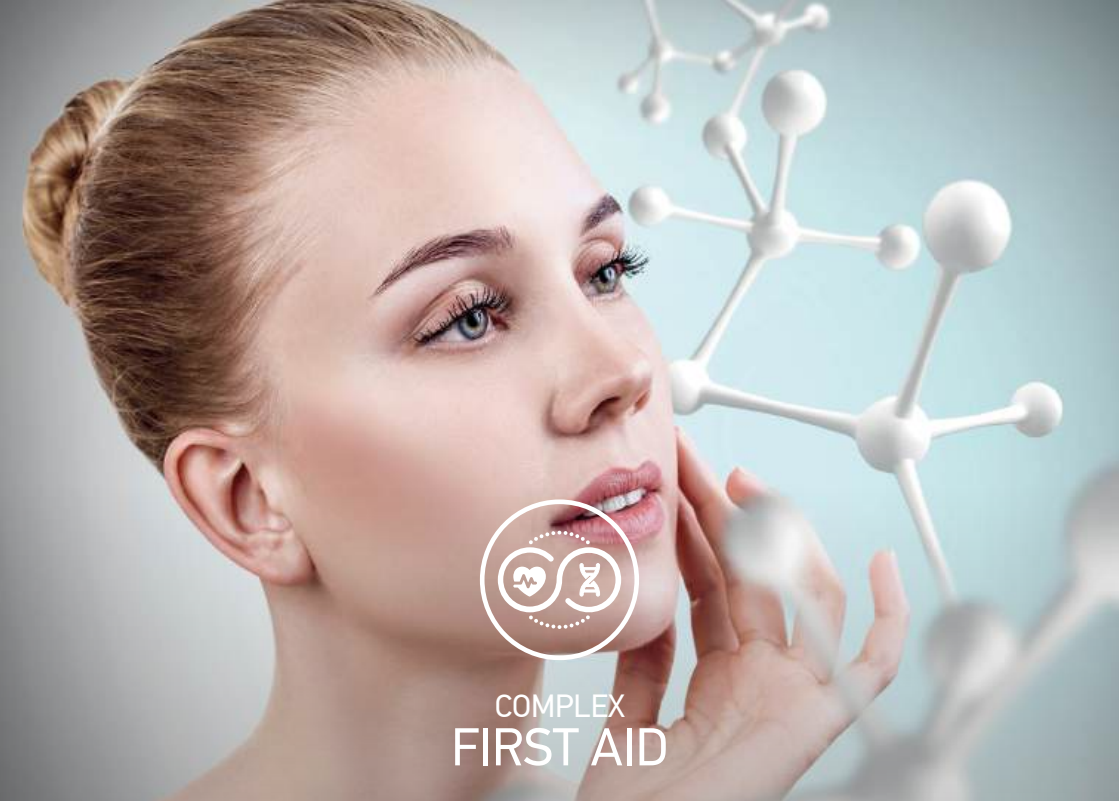

## COMPLEX FIRST AID

THE PROGRAMS OF THE “FIRST AID” COMPLEX ARE USED:

①

as needed;

②

as additional programs to the complexes:

- “Digestion without problems”
- “Stop insomnia”
- “Stop allergy”
- “Movement without pain”
- “Female Health”
- “Male Health”
- “Heart and Blood Vessels”

depending on the problem.



## COMPLEX HEALTH 60+

THE PROGRAMS OF THE “HEALTH 60+” COMPLEX ARE USED:

- ① as needed;
- ② as additional programs to the complexes:
  - “Digestion without problems”
  - “Stop insomnia”
  - “Stop allergy”
  - “Movement without pain”
  - “Healthy back”
  - “Female Health”
  - “Male Health”
  - “Heart and Blood Vessels”
  - “Without cold”
  - “Clear skin”

Depending on the problem, one of the complexes is chosen as a basis and works according to the specified scheme, supplementing it with the necessary programs of the complex “Health 60+”.



## COMPLEX HEALTHY CHILD

THE PROGRAMS OF THE HEALTHY CHILD COMPLEX ARE USED:

- ① as needed;
- ② as additional programs to the complexes:
  - “Digestion without problems”
  - “Stop insomnia”
  - “Stop allergy”
  - “Healthy back”
  - “Without cold”
  - “Clear skin”
  - “Stop parasites”

Depending on the problem, one of the complexes is chosen as a basis and works according to the specified scheme, supplementing it with the necessary programs of the complex “Healthy child”.

**Attention!** When working according to the schemes of the complexes, the “Healthy day -automatic” program for children under 14 is replaced with the “Day (child) - automatic”.

# LIST OF DEVITA AP MANUAL PROGRAMS

Nº	The name of the program	Program description
1	No staphylococcus	<p>According to statistics, 20% of the population are carriers of staphylococci. Staphylococci are common pathogens causing infections of the skin, soft tissue, upper respiratory tract and musculoskeletal system. Common sources of infection are hospitals and doctors who become bacteria carriers, as well as people taking antibiotics. Staphylococci cause secondary infections to influenza and wound infections, post-operative suppuration and food poisoning. The program eliminates the cause of the infection, which can have a harmful effect on the human body. It is also effective as part of a comprehensive treatment for diseases caused by these bacteria.</p> <p><b>Application:</b> Daily for up to 2 weeks, depending on how quickly the symptoms of infection subside.</p>
2	No streptococcus	<p>Streptococci are a widespread group of bacteria that can cause various diseases. They colonize the skin, respiratory tract and intestines of humans and can lead to acute inflammations in the ENT area, urogenital tract, abdominal organs, skin and soft tissue. Streptococcal infections are more severe than staphylococcal infections and can lead to joint, heart and kidney diseases as well as maxillary sinusitis, etc. The program is effective both in a complex and as a single program against otitis, tonsillitis and maxillary sinusitis.</p> <p><b>Application:</b> Daily for up to 2 weeks, depending on how quickly the symptoms of infection subside.</p>
3	No chlamydia	<p>Chlamydia is one of the most common sexually transmitted infections. Chlamydiosis very often occurs without or with mild symptoms. It primarily affects the urogenital tract. Women are more susceptible to this infection, and children can be infected when passing through the birth canal. Chlamydiosis can cause serious complications. Therefore, it is very important to take preventive measures in time, pay attention to personal hygiene and use DeVita technologies.</p> <p><b>Application:</b> Daily for 2-3 months.</p>
4	No helicobacter	<p>Helicobacter pylori is a chronic pathogen living in the human stomach and duodenum that can cause various dangerous diseases such as gastritis, stomach ulcers and hepatitis. According to the WHO, 60% of the world's population is infected with H. pylori. The infection is easily transmitted and is often considered a family disease, as if a family member is infected, there is a 95% chance that it will be passed on within the family. It is usually transmitted through the sharing of personal hygiene items and dishes and through kissing. This bacterium can also cause symptoms such as acid reflux, flatulence, stomach pain and a feeling of fullness.</p> <p><b>Application:</b> Daily for 2 weeks.</p>
5	No bacteria	<p>Many diseases are caused by bacteria. The most important feature of diseases caused by bacteria is the production of bacterial toxins, which cause inflammation and poisoning and attack the internal organs. The "No bacteria" program is used to alleviate the symptoms of intoxication and improve the general condition in case of bacterial infection of the digestive and respiratory tract caused by pneumococci, listeria, klebsiella, etc.</p> <p><b>Use:</b> Daily until the symptoms disappear.</p>
6	No hepatitis A	<p>In developing countries and areas with poor hygiene, there is an increased risk of contracting hepatitis A infection as the virus is transmitted through dirty hands, food and water, including seawater. If you are planning a trip to Asia or Africa, please observe the following safety measures: Pay attention to your personal hygiene, drink only bottled water and use the "No hepatitis A" program if you suspect an infection.</p> <p><b>Application:</b> Daily up to 2-3 times a day for at least one month.</p>
7	Without toxins	<p>The program detoxifies the body, improves drainage function and consequently increases immunity and vitality. It is used as a detox program after the antiparasitic programs but also as a single program for one month. Unlike other detox programs, this program does not use up the body's own energy reserves.</p> <p><b>Application:</b> Once a day after the antiparasitic programs. It is part of most auto-programs.</p>

8	No herpes zoster	<p>Herpes zoster (shingles) is caused by the varicella-zoster virus - the same virus that causes chickenpox. After a person has had chickenpox, the virus remains in a dormant state in certain nerve cells. Years later, it may reactivate as shingles, triggered by hypothermia, a cold or chronic illness: after an initial period of burning, tingling, itching, or stinging in the affected area, the characteristic rash appears in a stripe or band-like pattern along a nerve path, affecting one side of the body. The program aims to inactivate the virus and alleviate the main symptoms.</p> <p><b>Application:</b> Daily up to 3-4 times a day until the undesirable symptoms disappear. Use for at least one month in combination with the "No herpes" program.</p> <p><b>Important:</b> Do not use for chickenpox.</p>
9	No herpes simplex	<p>Herpes simplex has been known since ancient times. It is estimated that every third person has the herpes simplex virus in their body. It appears on the mucous membranes of the lips, mouth and genitals. If the immune system is weakened by hypothermia or overheating, an infectious disease or during menstruation, this can cause a rash on the mucous membranes. The virus is usually transmitted through contact with bodily secretions (blood, saliva, semen, vaginal secretions) of an infected person. The program effectively reduces the stress caused by herpes and fights off the formation of cold sores on the mucous membranes.</p> <p><b>Application:</b> At the first signs of a rash, use every three hours and then 1-2 times a day for one month in combination with the "No herpes" program.</p>
10	No herpes	<p>Over 90% of the world's population are carriers of the herpes virus. The virus is transmitted in various ways, including through saliva, sexual intercourse, hygiene items and droplet infection. In the beginning, the disease can be asymptomatic, then characteristic symptoms appear: headache, general malaise, fever and skin rash. The program helps reduce recovery time and relieve pain and inflammation.</p> <p><b>Application:</b> Use 1-2 times a day for one month together with the "No herpes zoster" or "No herpes simplex" program.</p>
11	Healthy lymph	<p>Lymph is a clear fluid that circulates through our tissues to cleanse them. It plays a major role in maintaining the water balance and homeostasis of the cellular microenvironment and is able to eliminate pathogenic microorganisms and toxins from the tissues. If the lymph flow is disturbed, the lymph cannot fulfill its main function, which leads to pathological processes. The program promotes the regulation of lymphatic drainage, the strengthening of lymphatic vessels and improves the condition of the lymph nodes, which leads to a reduction in swelling and better removal of metabolic waste products.</p> <p>The use of the program strengthens the body's immune system.</p> <p><b>Application:</b> If necessary, use once a day or according to the schemes of the complexes. It is part of the automatic program "Cleansing-automatic".</p>
12	Anti-Cold	<p>The common cold is one of the most widespread infectious diseases in the world. The causative agents are respiratory viruses (adenoviruses, respiratory syncytial virus) that affect the upper respiratory tract (nose, pharynx, trachea, larynx). Diseases caused by this group of viruses are contagious, because they are transmitted by airborne droplets and are not treated with antibiotics. Complications from a cold may include rhinitis, sinusitis, otitis media, etc. Very often, the weakening of the immune system can lead to the development of a bacterial infection. Therefore, it is advisable to use this program in combination with the antibacterial programs "No Staphylococcus", "No Streptococcus". The program is aimed at reducing inflammation of the mucous membranes of the nose, throat, bronchi, and activation of local immunity.</p> <p><b>Application:</b> at the first sign of a cold, the first 3 days - up to 2-3 times a day. Then 1 time per day for 7-14 days.</p>
13	No mould & fungi	<p>Fungi can be found everywhere. In the soil, in the air, in our bodies. In modern apartments and houses with heating systems, moisture can lead to mold growth. Popular ecological niches for mold are bathrooms, libraries and poorly ventilated living spaces. Fungi can enter the body through food, water, air and old dusty things. Molds cause inhalation allergies, as well as cross-reactions that manifest themselves on the skin. Molds can penetrate through the skin barrier wherever you take a steam bath, walk barefoot and there are suitable conditions for mold growth (moisture, high temperatures, organic residues such as dirt or sebum). Mold and other forms of fungi cause diseases of the respiratory tract, digestive tract and skin. To prevent fungal infections and reduce allergic reactions, we recommend using the "No mould &amp; fungi" program. The program is effective against black mould.</p> <p><b>Application:</b> 1-3 times a day for 4 weeks. To treat skin and nail fungal infections use it together with the "No fungi" program.</p>

14	No candida	<p>The program protects you from infection with the common fungus <i>Candida</i>. It colonizes the normal intestinal, vaginal and oral flora. However, if the immune system is weakened, the <i>Candida</i> fungus multiplies actively and damages the body. Diseases can be also caused by pathogenic <i>Candida</i> fungi that are transmitted through food, household items, after taking strong antibiotics or by sexual intercourse with an infected person. The symptoms are an allergic skin rash, vaginal mycosis, itching in the genitals and dysbiosis. You can prevent the proliferation of <i>Candida</i> by using the "No candida" program.</p> <p><b>Application:</b> The program is used 1-3 times a day for 2-4 weeks together with the "No protozoans" program.</p>
15	No fungi	<p>Fungal diseases in humans are infectious diseases caused by pathogenic fungi. The greatest danger of a fungal disease is that the affected areas are often an entry point for many other infectious agents. Fungi most frequently infect the skin and its appendages (hair, nails) and the mucous membranes of the internal organs. Treatment requires a comprehensive approach and personal hygiene. The program actively combats fungi such as <i>Trichophyton</i>, <i>Microsporum</i> and <i>Mucor</i>.</p> <p><b>Application:</b> To prevent fungal infections, it is advisable to apply the program every time after visiting a sauna, swimming pool or beach and after contact with an infected person. In the case of an existing fungal infection, the program is applied 1-2 times a day for one month.</p>
16	No protozoans	<p>The program protects you from the negative impact of most protozoa, which infect one in five people. Protozoa are transmitted through dirty hands, contaminated water or food. An infected person suffers from severe abdominal pain, bloating, constant fatigue and depression. The liver and pancreas can also be affected. Certainly, these symptoms make life difficult for us. There is also a risk of infecting your relatives. The program affects amoebae and giardia and kills them. Follow hygiene rules and the golden rule, which is "don't drink tap water", and use the "No protozoans" program regularly as a preventive measure to protect yourself from these microorganisms.</p> <p><b>Application:</b> 1-2 times a day for 10-14 days. In case of re-infection, the sources of infection should be determined.</p>
17	Intimate health	<p>There are many pathogens that cause diseases of the urogenital tract. These diseases are often caused by not just one, but several pathogenic microorganisms. The infections can be acute or chronic. Over time, a chronic infection can cause narrowing of the urethra, joint inflammation and infertility, and increase the risk of an ectopic pregnancy. Sometimes an exacerbation can occur during a cold or during psychological stress. Infections are usually transmitted via sexual intercourse and sometimes during childbirth from an infected mother. Bacteria such as ureaplasma can colonize the body asymptotically for years and go unnoticed. A key protective factor of intimate health is optimal microflora, and today's DeVita technologies work to kill pathogenic microflora and support the body's immune defenses. The program effectively combats parasites and pathogens such as <i>Ureaplasma</i>, <i>Toxoplasma</i> and <i>Mycoplasma</i>.</p> <p><b>Application:</b> Daily for a month.</p>
18	No schistosomes (blood-flukes) & trematods	<p>Trematodes (flukes) are flatworms that can cause serious illness in humans. They infect the blood (schistosomes), liver, intestines and lungs. The infection is transmitted by chewing on blades of grass, especially near bodies of water or when swimming in them, by eating raw salad or vegetables and by eating undercooked fish and red meat. The program aims to remove most types of flukes, including schistosomes, regardless of their developmental stage and location in the body.</p> <p><b>Application:</b> Daily for 14 days.</p>



19	No nematodes and ascarides	<p>The worm eggs, which are invisible to the naked eye, enter the human mouth through undercooked meat, unwashed fruit and dirty hands. The worms living in the host's body feed on its tissue, deprive it of nutrients (causing excessive hunger) and release toxins. An infected person loses weight, may suffer from hair loss and dandruff. The most common helminths are roundworms (ascarids), which can cause injury to the wall of the small intestine, intestinal obstruction, peritonitis, bronchitis, bronchial pneumonia, coughing up blood, as well as liver and pancreatic diseases. Ascariasis triggers an allergic reaction that manifests itself as a skin rash and bronchial asthma. The application of the program prevents the roundworms from feeding, which leads to their death and natural elimination from the digestive tract. The program has no toxic side effects. It is recommended that pet owners, parents of infants, animal keepers and anyone whose job involves dealing with animals follow the preventive measures once a quarter.</p> <p><b>Application:</b> The program is used daily for 7-10 days. At least 10 sessions in total. The cycle should be repeated after 14 days.</p>
20	No parasites, General	<p>Parasites can enter the body through the skin as well as through food, water and dust. Even with a healthy lifestyle, parasites cause health problems. A healthy diet, exercise and hardening of the body do not bring significant positive results without ridding the body of parasites. Parasites are the cause of many chronic diseases, such as gallbladder inflammation, gallstone disease, pancreatitis, colitis, diabetes, bronchial asthma and dermatitis. Symptoms such as irritability, anxiety, chronic fatigue, appetite disorders, anaemia, hyperactivity in children, headaches, problem skin, brittle nails and hair and immunodeficiency can be signs of parasite infections. Parasitic diseases are dangerous because they disguise themselves as other diseases. The prevention of parasitic diseases consists primarily of following hygiene rules and applying the "No parasites, General" program. These measures serve to maintain your health and the health of those around you.</p> <p><b>Application:</b> use the program daily for 7-10 days. At least 10 sessions in total. Repeat the course of treatment after 14 days.</p>
21	No parasites, Digestive tract	<p>Parasites in the gastrointestinal tract and liver cause biliary dyskinesia, cholangitis, gallbladder inflammation and hepatitis. In the intestine, they cause inflammation of the small intestinal wall and a reduction in nutrient absorption. Toxins released by parasites cause joint and muscle pain, affect the central nervous system and trigger depressive states, restlessness, nervousness and sleep disorders. The program aims to rid the body of parasites in the gastrointestinal tract, liver and pancreas.</p> <p><b>Application:</b> The program is to be used daily for 7-10 days. The course should be repeated after 14 days.</p>
22	Anti-Flu	<p>Influenza is an acute viral infection that affects people of any age. Up to 500 million people get the flu every year. Severe influenza can cause irreversible damage to the cardiovascular system, respiratory tract, and central nervous system. The virus is transmitted very easily by airborne droplets (by sneezing, coughing, talking). The most common complication of influenza is pneumonia. The program is aimed at eliminating various strains of influenza viruses and alleviating inflammatory manifestations.</p> <p><b>Application:</b> to prevent influenza during an epidemic, the program is to be used once every 2 days. If necessary - daily up to 2-3 times a day, alternating it with the "Against flu-General" program. After that it is necessary to use the "Drainage-General" program. It can be used together with the programs "No staphylococcus", "No streptococcus", "No candida" in order to avoid bacterial and fungal burden, as well as the programs "Antivirus", "Anti-Cold".</p>
23	Against flu - General	<p>Influenza is a disease that affects the upper and lower respiratory tract, causes significant intoxication and can lead to serious complications. Influenza always starts acutely and quickly (the first symptoms appear within 48 hours of infection). The main symptom is fever (40 °C and higher). It can last up to 5 days and be accompanied by shortness of breath, breathing difficulties, impaired consciousness and seizures. In people with a weakened immune system, the elderly and children, such temperatures can cause toxic shock syndrome (a condition causing vital organs to fail). If you use the program "Against flu-General", you can improve your condition (in case of illness), as well as restore your health quite quickly and without complications when using several programs against influenza viruses in combination.</p> <p><b>Application:</b> At the first signs of flu and elevated body temperature, use 1-3 times in a row or alternating with the "Anti-Flu" program. Be sure to drink plenty of fluids and maintain good personal hygiene.</p>

24	Healthy gastrointestinal	<p>Gastrointestinal infections can be caused by pathogens like bacteria, viruses, or parasites entering the body and causing inflammation. People of all ages can be affected, most often children and those who take antibiotics. Gastrointestinal infections are accompanied by nausea, vomiting, diarrhea, abdominal pain, fever. Symptoms of viral infections also include: runny nose and redness of the throat, as well as pain by swallowing, lack of appetite, loss of strength. The program is aimed at eliminating rotavirus, which causes gastroenteritis, and pathogenic Escherichia coli strains, normalizes well-being and prevents dehydration, which is harmful to our body.</p> <p><b>Application:</b> On the first day, use 4-5 times in a row to achieve a lasting result. Then use 1-2 times a day for 3-5 days, alternating with the "Drainage-General" and "Antivirus" programs. Make sure to drink plenty of fluids.</p>
25	Help against poisoning	<p>When enterobacteria or enteroviruses enter the body with food, water, or through the mucous membranes of the respiratory tract, an infection develops, which is commonly referred to as poisoning. This infection can be latent or show symptoms of its presence: headache, muscle pain, abdominal pain, fever, vomiting, stool disorder. Eventually it can lead to serious organ damage. The program helps the body cope with infection caused by enteroviruses and some types of salmonella, accelerating the recovery process, and reducing the likelihood of complications.</p> <p><b>Application:</b> Use it 4-5 times in a row on the first day to achieve sustainable results. Then use it 1-2 times daily for 3-5 days, alternating with the program "Drainage-General" and "Antivirus". Be sure to drink as much water as possible.</p>
26	Antivirus	<p>Viral infections are a large group of diseases caused by viruses. Viruses are connecting links between living and nonliving organisms. Since viruses have genetic material, but do not have their own cell, they are only intracellular parasites. They are difficult to destroy because they replicate inside host cells. Viruses can be transmitted by airborne droplets, through various objects, kissing and sexual intercourse. Our immunity plays a large role in destroying viruses in the body. The Antivirus program is aimed at combating various viruses, and in particular the cytomegalovirus and Epstein-Barr virus. It is also used to attenuate viruses and enhance the body's immunity.</p> <p><b>Application:</b> daily for a month.</p>
27	Drainage, General	<p>It is one of the basic detoxification programs. The lymphatic system is responsible for detoxification and lymphatic drainage. During an illness, the level of toxins in the body rises significantly. The lymphatic system is unable to eliminate them all, which can lead to intoxication. In addition, the natural drainage of wound fluid can be disrupted as a result of an illness or injury. The accumulation of fluid containing pathogenic microbes can cause severe inflammation and worsen the patient's condition. By using the "Drainage - General" program, you protect your body from the toxic load that often occurs as a result of various diseases and injuries. The program is aimed at neutralizing toxins and detoxification without stressing the lymphatic system.</p> <p><b>Application:</b> Once after using one of the DeVita AP programs. It is part of the automatic program "Cleansing-automatic".</p>
28	Drainage Antitox	<p>The extracellular matrix is the main area where toxins accumulate, leading to metabolic disorders, to the development of many serious chronic diseases and to premature ageing. Tissue drainage is a means of detoxifying the tissues. It helps remove salts of heavy metals and radionuclides from the body, improves metabolic processes and the general condition of the body.</p> <p><b>Application:</b> In case of severe poisoning use daily together with the program "Drainage-General".</p>
29	Mexican detoxification	<p>A method for deep and complete detoxification of the human body proposed by Dr. R. Rife. The program helps eliminate metabolic waste and toxins, improve the drainage function of the lymphatic system, increase the detoxification capacity of the liver, stimulate the excretory function of the kidneys and strengthen the immune system.</p> <p><b>Application:</b> Daily for one month, provided there are no serious diseases of the liver, kidneys or lymphatic system. If you do have such diseases, it is recommended to use the automatic program "Complete cleansing-automatic" on DeVita Ritm to perform a gentle detoxification of the body.</p>



30	No allergies	<p>An allergy is a hypersensitivity reaction of the immune system to certain substances: food, wool, dust, pollen, chemical substances and infectious agents. Today, over 85% of the world's population suffers from various allergies. The most common causes of allergies are a poor diet and an unhealthy lifestyle. The symptoms of allergies are a sudden runny nose, sneezing or tearing, redness and itching of the skin. The program alleviates the main symptoms of an allergy, improves general well-being and helps stabilize the immune system.</p> <p><b>Application:</b> Daily, if indicated and to alleviate the condition.</p>
31	Weight Control	<p>The program aims to combat adenovirus 36, which often causes obesity. It also regulates the autonomic nervous system and accelerates weight loss. The program also supports the body's energy reserves.</p> <p><b>Application:</b> It is recommended to use the program according to the "Body figure" complex.</p>
32	Without Salmonella	<p>The genus Salmonella is a member of the family Enterobacteriaceae. Salmonellae are resistant to the environment: they can survive in household dust, in open water, in frozen meat and sausage products, poultry eggs, dairy products from 3 to 6 months. They are killed at temperatures of 70 degrees Celsius and above. Infection is accompanied by: high fever, abdominal pain, vomiting, liquid stools. However, carriers may have no symptoms.</p> <p><b>Application:</b> Daily for 2 weeks.</p>
33	Without Yersinia	<p>Yersinia is an enterobacterium that causes inflammation in the small intestine, but can also affect other organs such as the liver, kidneys, heart and joints. Yersinia species are found in soil, water and food. They can tolerate low temperatures, allowing them to survive and multiply for a long time in the refrigerator. They can remain in the human body until the immune system is weakened, without triggering any inflammation. After that, an exacerbation begins.</p> <p><b>Application:</b> Daily for two weeks.</p>
34	Without Proteus	<p>Bacteria of the genus Proteus are part of the normal human intestinal flora and facultative human pathogens. As they multiply, they cause damage to the digestive tract organs and the urinary tract, birth defects and wound infections. The pathogenicity of this bacterium manifests itself in case of uncontrolled intake of antibiotics, or when bacteria enter from the outside through food, water or unwashed hands.</p> <p><b>Application:</b> Daily for two weeks.</p>
35	Without Klebsiella	<p>Klebsiella is part of the natural intestinal flora and a facultative pathogenic bacterium. In elderly and immunocompromised people, it can attack the organs of the digestive tract, the respiratory tract and the urogenital tract, causing e.g. purulent pneumonia or a urinary tract infection. This can lead to severe intoxication and consequently to functional disorders of the affected organs. Klebsiella can be transmitted through food, water and contact with infected people.</p> <p><b>Application:</b> Daily for two weeks.</p>
36	Without Campylobacter	<p>Campylobacter bacteria are a common cause of intestinal infections. They penetrate the mucous membranes and lymphatic vessels of the intestine and other organs. The infection is transmitted through livestock, poultry, consumption of contaminated meat, milk or water, less frequently through damaged skin and bites from infected animals, also transplacental or during childbirth. People at risk are rural residents, livestock and poultry farm workers, as well as tourists visiting developing countries.</p> <p><b>Application:</b> Daily for two weeks.</p>
37	Without Clostridia	<p>Clostridium difficile is part of the normal intestinal microbiota in 1–3% of healthy adults and 15–20% of infants. In cases of dysbiosis, it can cause damage to the walls of the large intestine. Clostridium botulinum is widespread in the soil and seawater; it produces one of the most potent bacterial toxins known, botulinum toxin.</p> <p><b>Application:</b> Daily for 2 weeks.</p>
38	Without Shigella	<p>Shigella is the causative agent of bacillary dysentery and a leading cause of bacterial diarrhea worldwide. Infection can occur through fecal-oral route, contaminated surfaces, water and food. It can also be transmitted by flies and cockroaches.</p> <p><b>Application:</b> daily for 2 weeks.</p>
39	Without Rotavirus	<p>Rotavirus (often called "stomach flu") multiplies in the cells of the gastrointestinal mucosa, causing severe diarrhea and dehydration. It is usually transmitted through food.</p> <p><b>Application:</b> Daily for 2 weeks.</p>

40	Anti-virus Gastrointestinal	Enterovirus and Coxsackie virus are enteroviruses that affect the central nervous system, muscles, heart, skin, respiratory system and the urogenital tract. They are transmitted by droplet and contact infection. Waves of infection rise in warm seasons. <b>Application:</b> Daily for two weeks.
41	Without Trematodes - General 1	Trematodes are parasitic flatworms commonly known as flukes. They infest the liver, pancreas, intestines, lungs, blood vessels and brain. The paths of infection are: swimming in natural waters in Asian, African and South American countries and eating raw or undercooked fish, shellfish, mussels and meat. <b>Application:</b> Use the two programs (Without Trematodes - General 1 and 2) in combination daily for two weeks.
42	Without Trematodes - General 2	
43	Stop Amoeba	Entamoeba histolytica is a protozoan parasite. It is the causative agent of amoebic dysentery and can infect and damage the large intestine, liver, lungs and other organs. It enters the body through contaminated food or water. <b>Application:</b> Daily for two weeks.
44	Without Enterococcus	Enterococci colonize the human intestine. When taking antibiotics, when the protective function of the skin and mucous membranes is impaired, in serious illnesses or in older people, enterococci can cause various disorders in the small intestine, urinary tract, heart, vagina, liver and skin. <b>Application:</b> Daily for 2 weeks.
45	Without nematodes and ascarididae - New	The program is aimed at fighting parasitic roundworms. It includes most nematode species that infect the human and animal body. <b>Application:</b> Daily for 2 weeks.
46	Stop Adenovirus - General	Adenoviruses are a group of viruses that cause acute infectious diseases affecting respiratory organs, eyes, and lymph nodes. These viruses are usually transmitted by droplet infection during the cold season. <b>Application:</b> Daily for 2 weeks.
47	Stop Coronavirus	The virus is widespread in Europe, Asia, and North America. It is transmitted by airborne droplets and causes damage to the respiratory tract, even to the point of destruction of pulmonary alveoli. <b>Application:</b> Daily for 2 weeks.
48	Stop Penicillium	Penicillium is a genus of mold that infects the organs of the urogenital and respiratory systems, the lymph nodes and the skin and causes severe intoxication and allergies. <b>Application:</b> Daily for two weeks.
49	Without Rickettsia	Rickettsia is an intracellular parasite. It can infect any organ of the body, causing significant CNS disorders and circulatory problems. Rickettsia causes diseases like typhus and Rocky Mountain spotted fever. It can be transmitted through a tick bite, food and direct contact, and is widespread in developing countries. <b>Application:</b> daily for 2 weeks.
50	Without Papilloma - General	Human papilloma viruses infect the skin and mucous membranes and cause warts. Some of them can lead to malignant tumors, with women being particularly at risk. HPV can be transmitted through direct contact, during sexual intercourse and at birth from mother to baby. <b>Application:</b> Daily for 2 weeks.
51	Without Papilloma	When HPV infects the mucous membranes of the urinary tract, formations occur in the renal pelvis and bladder. With prolonged exposure to unfavorable factors, these formations can not only reduce the functionality of the organs, but are also able to degenerate. <b>Application:</b> Daily for 2 weeks.

52	Without Trichomonas	Trichomonas is a protozoan parasite that infects the urogenital tract and causes trichomoniasis, a common sexually transmitted infection. The infection is almost always transmitted by unprotected sexual contact with an infected person, less frequently through contact with the toilet seat, in the swimming pool, by sharing towels, washcloths, hygiene items or by putting on other people's clothing (underwear). Once the pathogen enters the bloodstream, it can disguise itself as blood cells. The greatest danger does not usually come from the protozoa themselves, but from the pathogens they transmit: chlamydia, candida, gonococci, cytomegalovirus and herpes viruses. <b>Application:</b> Daily for two weeks.
53	Without Leptospira	Leptospira is a bacterium that can cause damage to capillaries, liver, kidneys, and muscles. It is spread worldwide. It can be transmitted through domestic and wild animals, through damaged mucous membranes and skin, and from natural water sources. <b>Application:</b> Daily for 2 weeks.
54	Without Cryptococcus	Cryptococcus is a yeast-like fungus. It can infect the central nervous system, lungs, skin and mucous membranes. The pathogen is often found in soil contaminated with bird droppings. It is transmitted by airborne infection. <b>Application:</b> Daily for 2 weeks.
55	Without Toxoplasma	The worldwide spread of Toxoplasma infection is very high, at about 90%. The pathogen usually remains in the body without any symptoms. It can be transmitted via the placenta or acquired in other ways. Toxoplasma infects the nervous system, eyes, heart, muscles, lungs, intestines and skin. <b>Application:</b> Daily for 2 weeks.
56	Without Ureaplasma	Ureaplasma are facultative pathogenic bacteria that cause diseases of the urogenital tract and joints as well as pathologies of pregnancy and the foetus. They can be transmitted sexually or during pregnancy or childbirth. <b>Application:</b> Daily for 2 weeks.
57	Without Echinococcus	A major problem are the larvae of the tapeworm, which usually form cysts. These cysts press on the surrounding tissue, leading to atrophy and necrosis. The larvae invade the liver, lungs, kidneys, bones and brain. Infection occurs through contact with infected animals, picking berries and herbs and drinking contaminated water. Butchers, shepherds and tanners are also at particular risk. <b>Application:</b> Daily for 2 weeks.
58	Without Pseudomonas aeruginosa	Pseudomonas aeruginosa is a multidrug-resistant, opportunistic pathogen. It infects the digestive tract, the respiratory system, the central nervous system, and other organs and systems. Sources of infection are humans and animals. Transmission routes: droplet infection, food infection, contact infection. <b>Application:</b> Daily for 2 weeks.
59	Against Lyme Disease 1	Lyme disease is the most common disease that can be transmitted through the bite of an infected blacklegged (Ixodes) tick. It affects the central nervous system, the heart and the joints. <b>Application:</b> After a bite, use the program daily for 2 weeks in combination with the program "Against Lyme Disease 2". For symptoms caused by previous bites, see the "Lyme Disease" complex.
60	Against Lyme Disease 2	The tick-borne encephalitis virus is transmitted through the bite of an infected tick and, more rarely, through the consumption of raw milk from an infected cow or goat. The endemic areas are: Siberia, Russian Far East, the Urals. The central nervous system is most affected. <b>Application:</b> After a bite, use the program daily for 2 weeks in combination with the program "Against Lyme Disease 1". For symptoms caused by previous bites, see the "Lyme Disease" complex.
61	Healthy Cell - 1	Numerous pathogenic factors lead to the accumulation of toxins in the body, which can cause cell degeneration. The programs "Healthy Cell - 1" and "Healthy Cell - 2" reduce the risk of cancer, help regulate tissue changes, inhibit inflammation and strengthen the immune system. <b>Application:</b> The programs "Healthy Cell - 1" and "Healthy Cell - 2" can be combined; see the complex "Healthy Cell".
62	Healthy Cell - 2	

63	Healthy Gums	<p>Gum disease (gingivitis) is caused by viruses, bacteria, and fungi. At risk may be people wearing dentures and braces. The program helps to reduce the pathogenic load, reduce the inflammatory process, and improve trophic tissue Gum inflammation (gingivitis) is caused by viruses, bacteria and fungi. People with dentures and braces are particularly at risk. The program helps reduce the pathogenic load, inhibit inflammation and improve tissue trophicity.</p> <p><b>Application:</b> Daily for 1-2 months.</p>
64	Mucous Membrane Regeneration	<p>The formation of aphthous ulcers on the oral mucosa indicates inflammation, which can be triggered by viruses, bacteria, injuries or mineral deficiencies. The program inhibits inflammation, serves to restore the oral mucosa and strengthens local immunity.</p> <p><b>Application:</b> This program is to be used daily for 2 weeks together with the "Healthy Mouth" program.</p>
65	Healthy Mouth	<p>Inflammation of the oral mucosa (stomatitis) is caused by pathogens, injuries to the oral cavity, dental procedures and an unbalanced diet. Stomatitis can involve the tongue, the palate and the mucous membrane of the cheeks and lips. The program inhibits inflammation, restores the oral mucosa and strengthens local immunity.</p> <p><b>Application:</b> Daily for 2 weeks. It can be combined with the program "Mucous Membrane Regeneration" and should be used after any dental procedures that involve damage to the oral mucosa.</p>
66	Without Fistula	<p>Pus in the gums is very often accompanied by the formation of a fistula. This can have serious consequences, from tooth loss to inflammation of the deep tissue. The program helps prevent the formation of pus and reduce inflammation.</p> <p><b>Application:</b> Daily for 2-4 weeks. It should be used immediately after dental procedures such as extraction, implantation or treatment.</p>
67	Stop Injury (teeth)	<p>The program relieves pain, accelerates the healing process and reduces swelling.</p> <p><b>Application:</b> Daily for 1-2 weeks. It should be used immediately after dental procedures such as extraction, implantation or treatment.</p>
68	CNS Without Problems	<p>The program helps improve the function of the nervous system after exposure to various stresses (infections, trauma, etc.). It promotes the restoration of nerve connections between the central nervous system and the periphery, relieves swelling and pain.</p> <p><b>Application:</b> Daily for 1-6 months. Preferably in combination with the "Healthy Nerves" program (see the "Lyme Disease" complex).</p>
69	Healthy Nerves	<p>Injury to the nerve fibers can lead to functional disorders of the nervous system. The program prevents the destruction of the myelin sheath of nerve fibers, thereby preserving their conductivity.</p> <p><b>Application:</b> Daily for 1-6 months. Preferably in combination with the program "CNS without problems" (see the complex "Lyme Disease").</p>
70	Clean Sinuses	<p>Sinusitis is an inflammation of the mucosal lining of one or more paranasal sinuses, which can occur as a result of infectious diseases of the nose and respiratory tract as well as after injuries to the facial area. The program helps reduce inflammation in the sinus area, improve mucus drainage and relieve pain in acute conditions.</p> <p><b>Application:</b> Daily for 2 weeks.</p>
71	Clean Sinuses - General	<p>The biggest problem for many people is chronic sinusitis, which may have no obvious symptoms but persists over a long period of time. The causes can be upper respiratory tract infections, polyps, deviated septum or decayed teeth. The program helps inhibit inflammation in the sinus area, reduce swelling and prevent exacerbations.</p> <p><b>Application:</b> Daily for one month if the problem is already present or to prevent an exacerbation during the cold season</p>
72	Ears Without Inflammation	<p>Ear inflammation can be caused by pathological processes in the nasopharynx, a deviated nasal septum, or ear injuries. The program helps reduce swelling, inflammation and pain.</p> <p><b>Application:</b> Daily for 2 weeks. Alternate with the DeVita Ritm program "Healthy Ears".</p>

73	Throat Without Pain	Inflammation of the pharyngeal tonsils (Tonsillitis) often occurs after a sore throat or other infectious diseases of the pharynx. It can be acute or chronic. The infection can be transmitted externally by airborne droplets or internally from a chronic internal infection. The program helps reduce tonsil inflammation, swelling, pain and toxic load. <b>Application:</b> Daily for 2 weeks. Alternate with the DeVita Ritm program "Throat Without Problems".
74	Against Cold - General	Inflammation of the nasal and pharyngeal mucosa (rhinopharyngitis) is usually caused by viruses and bacteria. The highest incidence of illness is observed between seasons. The program relieves inflammation and swelling, and prevents exacerbations. <b>Application:</b> Daily for 2 weeks.
75	Anti-inflammation - General	Inflammation is the body's response to injury or exposure to a pathogenic agent. The program helps the body to eliminate the source of inflammation, regardless of its location, to reduce the influence of pathogenic factors, and to restore homeostasis. <b>Application:</b> DDaily for 1-4 weeks. It is advisable to combine this program with any other programs in order to reduce inflammation.
76	Stop Migraine	Migraine is one of the most common types of headache. Symptoms include severe, often pulsating and throbbing pain, which can be limited to one side of the head and is aggravated by movement, light or noise. The program relieves vascular spasms, alleviates pain and promotes blood circulation in the head. <b>Application:</b> As needed.
77	Antiseptic	Skin injuries such as scratches, cuts and burns are entry points for pathogens. As a result, the wounds may take a long time to heal, become inflamed and cause pain. The program helps disinfect the wound, reduce inflammation and accelerate healing. <b>Application:</b> Daily for 1-4 weeks. It is advisable to combine it with any other programs to reduce inflammation.
78	Healthy Bile	Poor and irregular diet (prolonged fasting, alcohol, snack and spice consumption, too spicy, fatty and fried food, avoidance of vegetable oil), stress, low fluid intake and infectious diseases of the liver are the causes of gallbladder disorders. The program promotes the flow of bile, helps restore gallbladder function and relieves pain. <b>Application:</b> Daily for 1-6 months, preferably in the evening. Alternatively, every other day.
79	Without cough	Cough is a protective reflex reaction of the body that occurs due to irritation of the respiratory tract mucous membrane. It can be a sign of inflammation of the respiratory system, which can be caused by pathogenic viruses, bacteria, fungi and other agents, as well as allergens. The program alleviates the condition in the presence of acute or chronic cough, reduces inflammation of the mucous membrane of the respiratory tract. <b>Application:</b> daily for 2-4 weeks.
80	Without bacterial pneumonia	According to scientists, the first place among the causes of lower respiratory tract inflammation is occupied by microbial infections. Most bacteria that affect the respiratory system are part of the natural micro flora of the organism, but under adverse conditions they become pathogenic. The program fights against such microorganisms as: Streptococcus pneumoniae, Chlamydomphila pneumoniae, Mycoplasma pneumoniae, Klebsiella pneumoniae, etc. <b>Application:</b> daily for 2 weeks.
81	No Coronavirus - New	A new improved and expanded program that allows you to fight with all types of coronavirus, including new mutated forms. <b>Application:</b> daily for 2 weeks. In order to prevent the disease 2-3 times a week.
82	Without pneumonia	Inflammation of the lung tissue can be caused by pathogenic flora of various forms: viral, bacterial, fungal, etc. It arises not only as an independent form, but also as a result of the common cold, decreased immunity, and lack of vitamins. The program relieves from swelling and inflammation of the lower parts of the respiratory system, improves the functionality and energy of these organs, increases local immunity. <b>Application:</b> daily for 2-4 weeks. It is advised to combine with the program "Without cough".

83	Low back without inflammation	<p>As a result of pinched and inflamed spinal nerve roots in the lumbar and/or sacral spine, many people may experience sharp pain, often spreading into the leg. The most susceptible to the disease are persons who, due to their occupation, are forced to lift heavy weights, be in a bending position or sit for long periods of time. The program relieves pain and reduces inflammation.</p> <p><b>Application:</b> daily for 2-4weeks.</p>
84	Stop acne bacteria	<p>Propionibacterium is a conditionally pathogenic flora of human skin. Under the right conditions, as a result of blockage of the sebaceous gland ducts and increased sebum production, a favorable environment is created for their reproduction inside the hair follicles with the formation of acne (inflammation of the sebaceous glands). The activation of the pathogenetic properties of Propionibacteria may be linked to sex, age, and genetics. The program fights pathogenic microflora, and reduces skin inflammation.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without toxins" and "Skin without acne" programs on DeVita AP and the corresponding programs on DeVita Ritm (see Clear skin complex)</p> <p>In case of skin predisposition to acne formation, in order to prevent inflammation, it is desirable to use the program every quarter according to the above scheme.</p>
85	Without microsporium	<p>A fungus of the genus Microsporium primarily affects the skin, scalp, and hairy areas of the skin (beard, mustache), causing an infectious disease in animals and humans (ringworm). It is mainly transmitted from animals, less frequently from person to person or through household items. The program destroys pathogenic fungi, and it improves the condition of human skin.</p> <p>Application: daily for a month, in conjunction with the "Without toxins" program on DeVita AP and the corresponding programs on DeVita Ritm (see Clear skin complex).</p>
86	No onychomycosis	<p>Onychomycosis is a fungal lesion of the nails, usually caused by fungi of the genus Trichophyton, Epidermophyton floccosum and microsporidia. Men and the elderly are most likely to be affected. The development of onychomycosis provoke trauma, flat feet, varicose veins, hyperhidrosis, endocrine and immunological disorders, blood diseases, vascular pathology, prolonged use of hormonal drugs, antibiotics and cytostatics, wearing tight shoes and refusal to comply with common norms of personal hygiene. The program causes the death of pathogenic fungi, reduces inflammation, and helps restore the integrity of the skin and nails.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Nails without fungus", "Without trichophyton", "Without toxins" programs on DeVita AP and the corresponding programs on DeVita Ritm (see "Clean skin" complex). It is mandatory to maintain personal hygiene. To consolidate the result, it is desirable to repeat the cycle in a month. In order to prevent the occurrence of the disease - once or twice a year (fall-spring).</p>
87	Zygomycetes, stop!	<p>People with low immune defenses are most likely to get infected with these fungi. Persons at risk include long-term users of immunosuppressive medications, persons with diabetes mellitus, cancer or acquired immunodeficiency syndrome. The fungus affects internal organs, more commonly the CNS, lungs and GI tract. The program inhibits the development of fungal flora, and it improves the person's condition.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without toxins" program on DeVita AP and the corresponding programs on DeVita Ritm depending on the concerned organ. To prevent infection - 2 times a year for people with the above-mentioned health problems.</p>
88	Nails without fungus	<p>Mycosis of feet and nails is an infectious disease caused by pathogenic or conditionally pathogenic fungi of the species Trichophyton, Epidermophyton, Candida, Scytalidium. Risk factors are: violation of personal hygiene; sharing shoes; visiting public baths, swimming pools, and beaches; climatic conditions; injuries; wearing tight shoes; chronic diseases, etc.</p> <p>The program is aimed at destroying pathogenic flora, reducing inflammation, restoring the condition of the skin of the feet and nails.</p> <p><b>Application:</b> daily for a month, in conjunction with the "No onychomycosis", "Without trichophyton", "Without toxins" program on DeVita AP and the corresponding programs on DeVita Ritm (see Clear skin complex). It is mandatory to maintain personal hygiene. To consolidate the result, it is desirable to repeat the cycle in a month. In order to prevent the occurrence of the disease - once or twice a year (fall-spring), in case of visiting a swimming pool, bath, beach - every time after the visit.</p>

89	Without trichophyton	<p>Fungi of the genus Trichophyton affect the scalp, smooth skin and nails. Infection occurs through contact with a sick person or their belongings, as well as through contact with animals and their care items. The disease can occur in both acute and chronic forms. The program inhibits the development of fungi and inflammation.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without trichophyton", "Nails without fungus", "Without toxins" programs on DeVita AP and the corresponding programs on DeVita Ritm (see "Clean skin" complex). It is mandatory to maintain personal hygiene. To consolidate the result, it is desirable to repeat the cycle in a month. In order to prevent the occurrence of the disease - once or twice a year (fall-spring).</p>
90	Mycosis fungoides	<p>Mycosis fungoides affects primarily the lymphoid-reticular system of the skin, and subsequently the internal organs. The program helps to reduce the spread of fungus into the tissues and reduce the development of inflammation.</p> <p><b>Application:</b> daily, in conjunction with the "Without toxins" program on the DeVita AP and the "Healthy day" and "Health" programs on the DeVita Ritm.</p>
91	Solar keratosis	<p>The main cause of formation is due to prolonged exposure of the skin to direct sun rays, but manifests itself only with age when the immune forces are weakened. Therefore, it is most common in the elderly. Factors such as heredity, wearing tight clothing, taking certain medications, etc. can be triggers of the disease. The program helps to reduce the chances of keratosis manifestation and progression.</p> <p><b>Application:</b> daily, in conjunction with DeVita Ritm programs (see "Clean skin" complex).</p>
92	Stop warts	<p>Warts are benign skin neoplasms caused by human papillomaviruses. The transmission of papillomavirus can occur through household and sexual contact. The program fights HPV, preventing them from spreading in the body and forming neoplasms. It is necessary to consult a dermatologist regarding the removal of formations from the skin.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without toxins" program on DeVita AP and the programs DeVita Ritm (see Clear skin complex). Family members and sexual partners should use the program at the same time to prevent re-infection.</p>
93	Stop pimples	<p>Inflammation of hair follicles, provoked by the spread of staphylo-streptococcal infection with poor personal hygiene and/or the presence of chronic intoxication, infection, hypovitaminosis, leads to the formation of furuncles. The program inhibits the activity of bacterial flora, reduces inflammatory manifestations.</p> <p><b>Application:</b> daily for a month, in conjunction with the program "Mode 2 - automatic".</p>
94	Fingers without inflammation	<p>In the presence of microtraumas and injuries to the perionychium area of the fingers or toes and the development of bacterial, most often - staphylococcal, less often - streptococcal or mixed flora in the tissues of the fingers can lead to acute purulent inflammation - paronychia. The program reduces tissue inflammation, swelling, improves overall health, and reduces the risk of complications.</p> <p><b>Application:</b> daily for a month, in conjunction with the program "Mode 2 - automatic".</p>
95	No scleroderma	<p>Scleroderma is a connective tissue disease associated with impaired blood supply and subsequent thickening of tissues. The program improves blood supply to the tissues, slowing down the thickening processes.</p> <p><b>Application:</b> daily for a month, in conjunction with the program DeVita Ritm. The cycle is repeated 2 times a year.</p>
96	Erysipelas, stop	<p>Erysipelas is an infectious disease affecting the skin and mucous membranes. It is transmitted by airborne droplets, sometimes by contact. Injuries and microtraumas of the skin and mucous membranes serve as entry gates for infection. The program prevents the spread of the disease, helps to reduce inflammation and intoxication.</p> <p><b>Application:</b> daily, in conjunction with the "Without toxins", "No streptococcus" programs on the DeVita AP and the "Healthy day" and "Health" programs on the DeVita Ritm. The programs should be used in cycles 2-3 times a year for a month.</p>
97	Relieve psoriasis	<p>Psoriasis is a chronic non-infectious skin disease. Occurrence is associated with hereditary predisposition, immunity disorders, endocrine and nervous system disorders. The program helps reduce inflammation, prevent complications, and improve the quality of life.</p> <p><b>Application:</b> daily, in conjunction with DeVita Ritm programs (see "Clean skin" complex).</p>

98	Leprosy	<p>Lepra or leprosy is a chronic infectious disease affecting mainly the skin, peripheral nerves, upper respiratory mucosa and visual organs. It is characterized by a prolonged cycle (from a year to decades) with periodic exacerbations and can be disguised as other diseases, most often skin diseases. The program is aimed at destroying bacteria - the cause of the disease, reducing the inflammatory process, and preventing complications.</p> <p><b>Application:</b> in conjunction with the "Without toxins", "Restoration 2- Automatic," programs on the DeVita AP and the "Healthy day" and "Health" programs on the DeVita Ritm.</p>
99	Without Epidermophyton	<p>A genus of fungi that cause lesions (mycoses) of the skin, hair, and nails. It is the most common group of mycoses and occurs in people of all ages, sexes, occupations and health conditions. More commonly settles in the folds of the skin. The program inhibits the development of fungal flora, reduces inflammation, and improves the quality of life.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without toxins" program on DeVita AP and the corresponding programs on DeVita Ritm (see Clear skin complex). It is mandatory to maintain personal hygiene.</p>
100	Restoration of blood circulation	<p>A circulatory disorder is a condition in which blood does not flow in sufficient quantities to the brain, internal organs, and extremities. The problem can lead to severe complications, hypoxia, provoke deterioration of memory, hearing, etc. Such disorders develop amid some infectious or chronic diseases. The program helps to improve blood supply to organs and tissues especially after infectious diseases.</p> <p><b>Application:</b> daily for a month, in conjunction with the programs on DeVita Ritm (see Heart and blood vessels complex).</p>
101	Healthy spleen	<p>The spleen is the organ responsible for hematopoiesis, immunity, and blood supply. Being in a constant state of stress, irritability leads to problems with the functioning of this organ, and usually leads to a decrease in immunity. The program supports spleen function by reducing the adverse effects of negative psychological factors .</p> <p><b>Application:</b> daily for a month, in conjunction with DeVita Ritm programs (see "Stop stress" complex). In order to maintain immunity under constant stress, the scheme is used in a cycle of 2-3 times a year.</p>
102	Healing	<p>Healing of the body is impossible without the process of regeneration, which occurs at different levels of organization - system, organ, tissue, cellular, intracellular. The program activates the regeneration process in various inflammations (infectious, non-infectious) or damage to the integrity of organs and tissues (trauma), thereby enhancing the healing process.</p> <p><b>Application:</b> daily for one month. If necessary, the program is used until full recovery.</p>
103	Stop Gardnerella	<p>Gardnerella is a bacterium that is part of the normal biocenosis of women, actively multiplying when local immunity is weakened, causing inflammation of the genitourinary organs. The program induces the death of bacteria, and reduces inflammation.</p> <p><b>Application:</b> daily for 2-3 weeks, in conjunction with the "Without toxins" program on DeVita AP and the corresponding programs on DeVita Ritm (see "Women's Health" complex). To prevent re-infection, sexual partners should use the program at the same time.</p>
104	Male power	<p>Men have faced the problem of erectile dysfunction at all times and regardless of age. One of the causes of this problem is chronic prostatitis. The program reduces the inflammatory process, contributes to the restoration of normal functioning of the reproductive system, reduces the chance of complications (impotence, infertility).</p> <p><b>Application:</b> daily for a month, in conjunction with programs on DeVita Ritm (see "Men's Health" complex)"</p>
105	Without actinomycetes	<p>Actinomycetes are soil microorganisms. They form colonies on skin and mucous membranes without causing disease. They are capable of penetrating the human body only through damaged skin, mucous membranes. They become active under the influence of various factors, which include damage to the mucosa, the presence of infectious diseases, weakened immunity, diabetes mellitus, cancer, etc. The program prevents the spread of infection, and reduces the risks of complications.</p> <p><b>Application:</b> daily, in conjunction with the "Without toxins", "Restoration 2- Automatic" programs on DeVita AP .</p>



106	Stop gonococcus	<p>A bacterium that causes inflammation of the genitourinary tract. It's transmitted sexually. The program inhibits microorganisms, and reduces inflammation.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without toxins" program on DeVita AP and the corresponding programs on DeVita Ritm (see "Women's Health", "Men's Health" complexes ). To prevent re-infection, sexual partners should use the above scheme at the same time.</p>
107	Stop parvovirus	<p>A virus that affects the joints and causes chronic anemia. The virus is dangerous for people with weakened immune systems and blood diseases. The program induces the death of the virus, and reduces the risk of complications.</p> <p><b>Application:</b> daily month, in conjunction with the "Without toxins" program on DeVita AP.</p>
108	Without tuberculosis	<p>A chronic infectious disease caused by bacteria. It affects respiratory organs, bones and joints, genitourinary organs, eyes, peripheral lymph nodes. The following factors can provoke the disease: reduced immune defense, stress, chronic diseases of the lungs and GI tract, bad habits, living conditions. The program fights the root cause of the disease - the pathogenic flora.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without toxins" program on the DeVita AP and the "Healthy day" and "Health" programs on the DeVita Ritm.</p>
109	Without brucella	<p>Brucella is the bacteria that causes brucellosis. The disease is transmitted by ingestion of contaminated food (thermally unprocessed meat and dairy products), contact with an infected animal, or by aerosol transmission. The disease affects almost all organs, occurs in both acute and chronic forms. The program suppresses bacteria and reduces the toxic load on the body.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without toxins" program on DeVita AP and the corresponding programs on DeVita Ritm.</p>
110	Skin without acne	<p>Acne is a disease of the hair follicles and sebaceous glands that is more common in adolescents and young adults. There is a distinction between non-inflammatory and inflammatory acne. The program reduces inflammation, promotes the regulation of normal skin function, and improves the skin's condition.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Stop acne bacteria" program on DeVita AP and the corresponding programs on DeVita Ritm (see Clear Skin Complex).</p> <p>In case of skin predisposition to acne formation, in order to prevent inflammation, it is desirable to use the program every quarter according to the above scheme.</p>
111	Fungus-free skin	<p>Dermatomycosis is an infectious disease of the skin and its appendages, which is caused by various pathogenic, conditionally pathogenic fungi. The following factors contribute to the development of the disease: dermatitis, chronic diseases, poor personal hygiene, and professional factors. The program reduces toxic load and inflammation.</p> <p><b>Application:</b> daily for a month, combined with anti-fungal programs.</p>
112	Straight spine	<p>Spinal curvature is a disruption in the proper alignment of the spinal column. Any curvature of the spine can lead to not only back pain, but also internal organ problems. Various factors influence the curvature of the spine: weakness of the muscular corset of the back; sedentary behavior; asymmetrical types of physical activity; constant tilting of the torso forward when sitting; carrying a bag or backpack on one side; diseases of the musculoskeletal system. The program helps reduce muscle spasm, reduces pain, promotes muscle recovery, and improves the health of all tissues of the spine.</p> <p><b>Application:</b> daily, in conjunction with the "Healthy spine-automatic", "Healthy spine-automatic" programs on DeVita Ritm.</p>
113	Stop pain (low back)	<p>Discomfort in the lower back is common to almost every person in their 30s and older. There are many reasons for this: circulatory disorders; degenerative changes; trauma; infections, etc. The program reduces pain and improves blood flow in the lower back area.</p> <p><b>Application:</b> daily for a month or as needed in conjunction with the "Low back without inflammation" program on the DeVita AP and the "Joints without pain-automatic" and "Healthy spine-automatic" programs on the DeVita Ritm.</p>
114	Without allergies new	<p>The occurrence of allergies in humans is influenced by a number of factors. Among them: household allergens, working conditions, genetic predisposition and other factors. These conditions can not only cause allergic reactions, but can contribute to the development of asthma. The program reduces allergies, promotes better breathing by relieving spasm, and it improves well-being.</p> <p><b>Application:</b> daily for a month or as needed, in conjunction with the corresponding programs on program on DeVita AP and the DeVita Ritm (see Stop allergy complex).</p>

115	Without fluke	<p>The Siberian or cat liver fluke a flatworm which belongs to the class of trematodes. It causes opisthorchiasis, affecting organs such as the gallbladder, pancreatic ducts, and bile ducts of the liver. Infection occurs by eating low-salted, undercooked food, raw fish; by using household appliances that have come into contact with such fish; by failing to maintain personal hygiene after visiting a restroom or cleaning animal toilets and interacting with them. The program deals with the parasite at all stages of its development.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Without toxins" program on DeVita AP and the corresponding programs on DeVita Ritm (see Stop parasites complex). It is mandatory to maintain personal hygiene. To consolidate the result, it is desirable to repeat the cycle in a month.</p>
116	Stop dementia	<p>A number of reasons can lead to impaired thinking abilities in a person: vascular diseases, Alzheimer's disease, trauma, brain neoplasms, alcoholism, drug addiction, CNS infections and some other diseases. The program helps maintain brain function.</p> <p><b>Application:</b> daily for a month, repeat the cycle at regular intervals every quarter. Use in conjunction with the " Blood circulation in the head" program and the corresponding programs on the DeVita Ritm.</p>
117	Stop Parkinson	<p>The disease is characterized by damage to nerve cells in the brain stem and is expressed by slowness of movement, tremors at rest and impaired reflexes. It is common in the elderly. The program helps to slow the progression of the disease.</p> <p><b>Application:</b> daily for a month, repeat the cycle at regular intervals every quarter. Use in conjunction with the " Blood circulation in the head" program and the corresponding programs on the DeVita Ritm.</p>
118	Stop pain (injury)	<p>Nobody is safe from injury. Injuries are accompanied by pain, swelling, violation of tissue integrity and other signs. The program reduces pain, swelling, and speeds up the regeneration process.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Stop Injuries" and "Blood circulation regulation" program on DeVita Ritm.</p>
119	Strong nerves	<p>Amid a prolonged nervous overstress due to prolonged illness, life troubles, loss of loved ones and other reasons can arise nervous disorders characterized by irritability, general fatigue, apathy, mood swings. The program helps to normalize the psychological state, helps to improve mood and general condition.</p> <p><b>Application:</b> daily for a month, in conjunction with DeVita Ritm programs (see "Stop Stress" complex).</p>
120	Good sleep	<p>Healthy sleep is the key to good mood and well-being. But stressful situations at home and work, lifestyle changes and other reasons can trigger insomnia. The program reduces the irritability of the nervous system and promotes restful sleep.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Strong nerves" program and the corresponding programs on DeVita Ritm (see Stop insomnia).</p>
121	Calm intestines	<p>Stomach heaviness, heartburn, bloating, rumbling, meteorism, stool problems can all be signs of irritable bowel. Such a condition is typical for people in developed countries. The program reduces the severity of pain, relieves mucosal inflammation, normalizes GI motility.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Intestines without disorders" program on DeVita AP and the corresponding programs on DeVita Ritm (see Digestion without problems complex).</p>
122	Healthy Stomach	<p>Stomach pain is familiar to almost every other person. The causes of such pains are: functional disorders, intestinal infections, gastritis, etc. Stress, improper nutrition, taking some medications can provoke pain. The program relieves pain, reduces inflammation, helps normalize the stomach.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Healthy stomach" program on DeVita Ritm and the "Stop stress" or "Digestion without problems" complex programs. In order to prevent exacerbations, the cycle is carried out 2-3 times a year.</p>
123	Duodenitis, stop!	<p>Eating disorders, consumption of foods that irritate the mucosa of the GI tract, diseases of the stomach, liver and pancreas can cause inflammation of the mucous membrane of the duodenum and pancreas, which is often accompanied by severe pain. The program relieves pain and helps mucous membrane recovery.</p> <p><b>Application:</b> daily for a month, in conjunction with the corresponding programs on the DeVita AP and the "Stomach+ Automatic" program on the DeVita Ritm.</p>

124	Healthy pancreas	<p>The pancreas plays an important role in the digestive and endocrine systems, it produces enzymes that aid the digestive process and the hormone insulin that regulates blood sugar levels. Therefore, it is very important to keep it functioning properly. The program reduces discomfort and maintains the function of the pancreas.</p> <p><b>Application:</b> daily for a month, in conjunction with the corresponding programs on DeVita AP and DeVita Ritm.</p>
125	Blood sugar under control	<p>Stress, high-carbohydrate diets, medications, and endocrine and internal organ diseases lead to high blood sugar. It is possible to balance sugar levels by combining diet, exercise and the "Blood sugar under control" program.</p> <p><b>Application:</b> daily for one month.</p>
126	Stop diabetes	<p>Diabetes is a disease characterized by chronically elevated blood glucose levels. There is a risk of cardiovascular disease, loss of vision, decreased immunity, and impaired blood supply to the lower extremities. The program helps to minimize possible risks.</p> <p><b>Application:</b> daily for a month, in conjunction with the program "Blood sugar under control" on DeVita AP and basic programs on DeVita Ritm complex "Digestion without problems" and special programs "Normal blood sugar", "Healthy Pancreas". The full cycle is performed 2 times a year, the rest of the time the DeVita Ritm programs are used.</p>
127	Small intestine	<p>The small intestine performs several important functions: digestion of food; absorption of nutrients; movement of remaining substances to the GI tract; hormone production; and immunological defense. Most often people are faced with functional disorders and inflammatory processes in this part of the intestine. The program helps to reduce pain and inflammatory manifestations, helps to normalize bowel function.</p> <p><b>Application:</b> daily for a month, in conjunction with the corresponding programs on program on DeVita AP and the DeVita Ritm (see Digestion without problems complex).</p>
128	Intestines without disorders	<p>Disorders of the intestines can appear with functional impairments in its work, infectious diseases or stress. They are accompanied by abdominal pain, bloating, and constipation and/or diarrhea. The program reduces pain, helps normalize the work of the large intestine.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Calm intestines" program on DeVita AP and the corresponding programs on DeVita Ritm (see Digestion without problems complex).</p>
129	Without hemorrhoids	<p>Hemorrhoids are the result of impaired blood circulation and inflammation in the blood vessels of the rectum. The disease is equally common in men and women, and brings great discomfort. The program relieves pain, reduces inflammation, and improves vascular health.</p> <p><b>Application:</b> daily for one month. To minimize the occurrence of the disease, it is necessary to control the influence of risk factors; take a cycle using the schemes of complexes "Digestion without problems", "Heart and blood vessels" or others, depending on the cause of the problem.</p>
130	Stop pelvic inflammation	<p>Inflammatory processes in the pelvic area can be both independent diseases and complications of previously experienced diseases. One of the most common causes can be considered bacterial and viral infections. The program reduces inflammation and relieves pain.</p> <p><b>Application:</b> daily for a month, in conjunction with the DeVita AP and DeVita Ritm programs for relevant infectious agents (see "Women's Health" and "Men's Health").</p>
131	Stop chickenpox	<p>Chickenpox is an infectious disease caused by the varicella-zoster virus that affects the skin and nervous system. It is transmitted by airborne and contact. It is more severe in adults. The program fights the pathogen, reducing the manifestation of intoxication and the occurrence of complications.</p> <p><b>Application:</b> daily for 2 weeks, in conjunction with the "Without toxins" program on DeVita AP.</p>
132	Stop rubella	<p>Rubella is an acute viral infection that is very "contagious". It's transmitted by airborne droplets. The program inhibits the development of the virus and its spread in the body.</p> <p><b>Application:</b> daily month, in conjunction with the "Without toxins" program on DeVita AP.</p>
133	Without mumps	<p>Mumps is an acute infectious disease caused by a virus of the paramyxovirus genus. The main signs of the disease are intoxication and inflammation of the parotid salivary glands. The virus is transmitted by airborne and household contact. The program fights the virus, reducing the manifestation of intoxication and the occurrence of complications.</p> <p><b>Application:</b> daily for 2 weeks, in conjunction with the "Without toxins" program on DeVita AP.</p>

134	Stop scarlet fever	Scarlatina is an acute bacterial infection that is transmitted by airborne and household contact. Characterized by disease of the upper respiratory tract, tonsils. The program destroys pathogenic bacterial flora, reducing intoxication, reducing the risks of complications. <b>Application:</b> daily month, in conjunction with the "Without toxins" program on DeVita AP.
135	Stop measles	Rubella is an acute viral infection that is very "contagious". It's transmitted by airborne droplets. The program suppresses the development of the virus, relieves signs of intoxication, reduces the chances of complications. <b>Application:</b> daily month, in conjunction with the "Without toxins" program on DeVita AP.
136	No diphtheria	Diphtheria is an acute bacterial infection that is transmitted by airborne droplets and airborne dust. It is possible to be infected through the contact and domestic means (utensils, dirty hands, etc.), as well as alimentary (food). The disease is characterized by lesions of the upper respiratory tract, and oropharyngeal mucosa. The program fights the pathogen, reduces signs of intoxication and prevents the spread of inflammation to other organs. <b>Application:</b> daily month, in conjunction with the "Without toxins" program on DeVita AP. Use the appropriate programs on the DeVita Ritm and DeVita AP if necessary.
137	Stop cholera	Cholera is an acute bacterial infection, transmitted by fecal-oral route through water or contaminated food. The main symptoms are vomiting, diarrhea, and abdominal pain. The program inhibits the development of the bacteria, which helps to reduce intoxication and improve the person's condition. <b>Application:</b> daily month, in conjunction with the "Without toxins" program on DeVita AP. It is mandatory to maintain personal hygiene. In order to prevent the disease, for persons visiting countries in Asia (India), Africa, South America, each time during and/or after travel - daily for 2 weeks.
138	Stop pertussis	Pertussis is an acute bacterial infection with damage to the respiratory, cardiovascular and nervous systems. The path of transmission of the pathogen is airborne. The program fights the pathogen by reducing the effects of toxins on cells and preventing inflammation from spreading to other organs. <b>Application:</b> daily for a month, in conjunction with the "Without toxins" program on DeVita AP and the "Cellular respiration" and "Free Breath" programs on DeVita Ritm.
139	Without constipation	Improper nutrition (lack of dietary fiber, liquids), sedentary work, stress, change of climatic area during relocation, impaired motility of the large intestine and other factors can cause difficult or systematically obstructed defecation. Prolonged chronic constipation can contribute to the development of complications in the digestive tract. The program helps to stimulate bowel function and strengthen the intestinal muscle tone. <b>Application:</b> daily for a month, in conjunction with the corresponding programs on DeVita AP and DeVita Ritm (see "Digestion without problems", "Stop stress" complexes. To prevent further constipation, introduce fiber-rich foods into your diet and eliminate foods that aggravate constipation. Increase fluid intake. With a predisposition to this problem, frequent relocations, stress, the cycle is carried out 2-3 times a year.
140	Blood circulation in the head	Insufficient blood supply to the organs of the head can lead to undesirable consequences: headache, dizziness, tinnitus, decreased visual and hearing acuity, impaired perception, attention, memory. The program helps to improve blood flow to the tissues of the head organs, thereby increase the flow of oxygen and nutrition to them. <b>Application:</b> daily for a month, in conjunction with the corresponding programs on program on DeVita AP and the programs DeVita Ritm (see Heart and Blood Vessels complex).
141	For women	Abundant and painful menstruation, irregularity of the cycle, may be the cause of endometriosis. One of the complications of this disease is infertility. There are many causes of the disease: hormonal disorders, decreased immunity, heredity, etc. The program helps to restore hormonal balance, normal functioning of the reproductive system. <b>Application:</b> daily for a month, in conjunction with the "Stop pelvic inflammation" program on DeVita AP and the programs DeVita Ritm (see Women's health complex).

142	Stop spasm (muscles)	<p>The superficial and deep back muscles are like a corset, supporting the spine in the correct position. Lifting heavy weights, spending long periods of time in an uncomfortable position, poor posture, problems with the musculoskeletal system cause spasms of the back muscles, accompanied by pain. The program helps to relieve spasm, pain, reduces inflammation, and improves a person's condition.</p> <p><b>Application:</b> daily for a month, in conjunction with the DeVita AP "Muscles without tension" program and the DeVita Ritm "Healthy spine-automatic" and "Joints without pain-automatic" programs. To strengthen the muscular corset of the back, regular training of all back and abdominal muscles using special DeVita Ritm programs is essential (see "For Athletes" complex).</p>
143	Without motion sickness	<p>Motion sickness is caused by the vestibular apparatus – the organ of balance and orientation in space, which is located in the inner ear. When motion sickness occurs, feelings of nausea and dizziness are experienced. The program relieves cramps, nausea, and improves health.</p> <p><b>Application:</b> every time you travel in transport, in conjunction with the "Blood circulation regulation" program on the DeVita Ritm.</p>
144	Stop spasm (Intestines)	<p>The pain associated with intestinal spasm can occur in connection with violations of various mechanisms of intestinal motility, with intoxication, against the background of digestive disorders, stressful situations, helminthiasis, etc. The program helps to relax the intestinal wall, reduce pain.</p> <p><b>Application:</b> daily for a for one month, in conjunction with the corresponding programs on DeVita AP and DeVita Ritm (see "Digestion without problems", "Stop stress" complexes).</p>
145	Stop conjunctivitis	<p>Conjunctivitis is an inflammation of the mucous membrane, occurs most often with viral or bacterial infection, allergies. It is characterized by redness, swelling and itching. The program helps to reduce inflammation, swelling and itching.</p> <p><b>Application:</b> for infections daily for a month, in conjunction with the appropriate programs on the DeVita AP . For allergies, daily for the entire period of exposure to allergens, in conjunction with the corresponding programs on DeVita AP and DeVita Ritm (see "Stop Allergy" complex).</p>
146	Muscles without tension	<p>Muscle tension can be caused by injuries, uncomfortable posture, poor posture, osteochondrosis, stress, or disease. The program reduces pain, increases muscle elasticity, and reduces muscle tone.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Stop spasm (muscles) program on DeVita AP and "Muscles without pain" on DeVita Ritm.</p>
147	Stop fatigue	<p>Fatigue is characterized as a state of exhaustion, manifested by a feeling of weakness, lethargy, feelings of physiological discomfort, disturbances in the flow of mental processes. Fatigue can be the result of overwork, poor sleep, worries, lack of exercise, anxiety, and illness. It is most frequent in residents of large cities. The program reduces the irritability of the nervous system and helps to increase performance.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Strong nerves", "Blood circulation in the head" programs on DeVita AP and the programs DeVita Ritm (see Stop stress complex).</p>

# DEVITA AP AUTOMATIC PROGRAMS DESCRIPTION

Name	Effect of the program	Programs	Working time
Mode 1 Automatic	It inhibits the activity of viruses, inhibits the vital activity of protozoan microorganisms. Optimal disease prevention when changing climatic zones, visiting Asian countries.	Antivirus, No protozoans, Without toxins	6:47
Mode 2 Auto-matic	Relieves toxic load caused by streptococcal and staphylococci. It is effective both for colds and for the prevention of coccal infections.	No staphylococcus, No streptococcus, Without toxins	7:53
Mode 3 auto-matic	It contributes to the removal of bacterial aggravation of the digestive tract and the genitourinary system.	Healthy gastrointestinal, Intimate health, Without toxins	8:26
Mode 4 Automatic	It removes mycotoxins formed by mold and yeast. Reduces the risk of accumulation of carcinogenic factors	Stop candida, No mould & fungi, Without toxins	7:47
Cleansing-automatic	It promotes the elimination of toxins, the normalization of lymph outflow, and the strengthening of lymphatic vessels. Helps the removal of excess fluid from the intercellular space, has an antiseptic effect.	Healthy lymph, Drainage - General	7:10
No parasites Automatic	Reduces the stress parasitic burden of the liver, pancreas, intestines, as well as the toxic load caused by the vital activity of these forms.	No parasites, Digestive tract, No parasites, General, Without toxins	7:38
No flu - automatic	It suppresses the activity of various strains of the influenza virus, helps relieve symptoms of intoxication and alleviates the human condition during the illness.	Anti-Flu, Against flu - General, Without toxins	6:38
Stop herpes S-automatic	It prevents the development of herpes simplex and its toxic effect on the body: the formation of vesicles on the mucous membranes, skin, reduces itching.	No herpes, No herpes simplex, Without toxins	7:35
No herpes Z-automatic	It reduces the activity of zoster herpes viruses, reduces the toxic signs of its manifestation on the surface of the skin and pain.	No herpes zoster, No herpes, Without toxins	5:47
Without cold automatic	Prevents the development of respiratory viruses, reduces symptoms of intoxication, leads to a decrease in inflammatory reactions and the likelihood of complications	Antivirus, Anti-Cold, Without toxins	8:09
Support immune system-Automatic	It blocks the activity of various viruses, their distribution in the human body, promotes the activation of immunity. It is effective as prevention of colds infections during hypothermia and during the period of mass diseases.	Antivirus, Without toxins	5:03
Weight Reduction-Automatic	Reduces parasitic load and destroys viruses that are the cause of excess weight. It normalizes the functioning of the autonomic nervous system, maintains the body tone and improves the psychoemotional state.	Weight Control, Without toxins	4:28
Reduce Points-automatic	It cleanses the body of toxins, normalizes lymphatic drainage, supports the immune system and increases vital energy.	Healthy lymph, Without toxins	4:31
Cell Protection 1 -Automatic	The program helps to relieve viral, bacterial and toxic loads, increasing the body's defenses.	Antivirus + No streptococcus + No herpes simplex + Without toxins	12:09
Cell Protection 2 -Automatic	Prevents the effects of carcinogenic factors on the human body, thereby reducing the risk of cell degeneration.	Healthy Cell 1 + Healthy Cell 2 + Antivirus Gastrointestinal + No Trichomonas	13:54

Cell Protection 3 –Automatic	The program helps to resist viral, bacterial and protozoa infections, reducing the degree of intoxication.	No herpes zoster + No staphylococcus + No chlamydia + Intimate health + Stop Amoeba	12:42
Cell Protection 4 –Automatic	The program reduces the fungal toxic load, increasing your own protection and reduces the likelihood of tumor processes.	Healthy Cell 1 + Healthy Cell 2 + Stop candida + Stop Penicillum	12:48
Lyme 1 - Automatic	The program helps to remove viral, bacterial and toxic loads, increasing the body's defenses.	Antivirus + No streptococcus + No herpes simplex + Without toxins	12:09
Lyme 2 - Automatic	The program works to destroy pathogens that have been ingested by a tick bite and to reduce the total toxic load.	Against Lyme Disease 1 + Against Lyme Disease 2 + Antivirus Gastrointestinal + No Trichomonas	11:39
Lyme 3 - Automatic	The program helps to resist viral, bacterial and protozoa infections, reducing the degree of intoxication.	No herpes zoster + No staphylococcus + No chlamydia + Intimate health + Stop Amoeba	12:42
Lyme 4 - Automatic	The program reduces the fungal toxic load, increasing the body's defenses while fighting borrelia.	Against Lyme Disease 1 + Against Lyme Disease 2 + Stop candida + Stop Penicillum + Without toxins	11:51
Restoration 1 - Automatic	The program is aimed to relieve pain after dental interventions, repair damaged mucosa, reduce swelling, and accelerate tissue regeneration.	Stop Injury (teeth) + Healthy mouth	6:33
Restoration 2 - Automatic	The program reduces the influence of pathogenic factors, reducing the degree of the inflammatory process. Promotes speedy recovery.	Antiseptic+Anti-inflammation - General	6:09
Healthy Gums - Automatic	The program helps to reduce the inflammatory process of the oral cavity, improves trophic tissue, strengthens local immunity.	Healthy Gums+Mucous membrane Regeneration	6:48
Nerve Regeneration - Automatic	The program is aimed at improving the functioning of the nervous system in case of damage by pathogens that have entered the body due to a tick bite. It reduces swelling, helps restore the integrity of the nerve fiber and nerve connections of various parts of the nervous system.	Healthy Nerves+CNS Without Problems	4:51
Healthy Breathing	The program aims to suppress the activity of viruses and bacteria that cause respiratory tract diseases, helps to reduce the likelihood of exacerbations, reduce intoxication, increase the body's defenses.	Stop coronavirus NEW + Bacterial pneumonia + Without toxins	10:37
Without motion sickness -automatic	The program helps relax the intestinal wall, reduces pain, nausea, and promotes better health. <b>Application:</b> use preferably 30 min before travel.	Without motion sickness Stop spasm (Intestines)	4:35
Sos! Allergy - Automatic	The program reduces allergic manifestations, promotes better breathing by relieving spasm, reducing inflammation, relieving swelling and itching. <b>Application:</b> daily for the entire period of exposure to allergens or as needed, in conjunction with the corresponding programs on DeVita Ritm (see "Stop Allergy" complex).	Without allergies new Stop conjunctivitis	6:14

Restoration 3 -automatic	<p>The program improves blood supply to organs and tissues, activates the regeneration process especially after infectious diseases or injuries.</p> <p><b>Application:</b> daily for a month, if necessary use the program until full recovery, in conjunction with DeVita Ritm programs (see Heart and blood vessels complex).</p>	Restoration of blood circulation Healing	2:53
Deep sleep- automatic	<p>The program normalizes the psychological state, reduces the excitability of the nervous system, improves mood and promotes restful sleep.</p> <p><b>Application:</b> daily for a month, in conjunction with DeVita Ritm programs (see "Stop Stress" complex) or "Stop insomnia"</p>	Strong nerves Good sleep	6:20
Stop brain aging - automatic	<p>The program helps maintain brain function, slows down the process of neuron destruction, helps improve blood supply to the tissues of the head organs, thereby increasing the flow of oxygen and nutrition to them.</p> <p><b>Application:</b> daily for a month, in conjunction with the corresponding programs on program on DeVita Ritm (see Heart and blood vessels complex). Repeat the cycle regularly every quarter.</p>	Stop dementia Stop Parkinson Blood circulation in the head	9:36
Stop fatigue- automatic	<p>The program reduces the agitation of the nervous system, helps to normalize the psychological state, improves blood supply to the tissues of the organs of the head, helps to increase performance.</p> <p><b>Application:</b> daily, in conjunction with DeVita Ritm programs (see "Stop stress" complex).</p>	Strong nerves Blood circulation in the head Stop fatigue	9:09
Anti-inflammation (female) -automatic	<p>The program reduces the inflammatory process, relieves pain, helps to restore hormonal balance, normal functioning of the reproductive system.</p> <p><b>Application:</b> daily for a month, in conjunction with the DeVita AP and DeVita Ritm programs for relevant infectious agents (see "Women's Health" complex).</p>	Stop pelvic inflammation For women	6:35
Anti-inflammation (male) -automatic	<p>The program helps to reduce the inflammatory process, restore the normal functioning of the reproductive system, reduces the chance of complications (impotence, infertility).</p> <p><b>Application:</b> daily for a month, in conjunction with the DeVita AP and DeVita Ritm programs for relevant infectious agents (see "Women's Health" and "Men's Health" complexes).</p>	Stop pelvic inflammation Male power	5:56
Active Gastro -automatic	<p>The program helps to relieve pain, reduce inflammation, restore mucosa, normalize the stomach, gallbladder, and pancreas.</p> <p><b>Application:</b> daily for a month, in conjunction with the DeVita Ritm program "Stomach+ automatic" and the programs of the complex "Stop stress" or "Digestion without problems". In order to prevent exacerbations, the cycle is carried out 2-3 times a year.</p>	Healthy Stomach Duodenitis, stop! Healthy pancreas	9:21
Bowel without irritation	<p>The program reduces pain, meteorism, bloating, relieves inflammation of the mucosa, normalizes the motility of the GI tract.</p> <p><b>Application:</b> daily for a month, in conjunction with the corresponding programs on DeVita Ritm (see "Digestion without problems" complex).</p>	Calm intestines Intestines without disorders	7:56



Without constipation -automatic	<p>The program helps to relieve spasms, reduce pain, normalizes the motility of the large intestine.</p> <p><b>Application:</b> daily for a for one month, in conjunction with the corresponding programs on DeVita AP and DeVita Ritm (see "Digestion without problems", "Stop stress" complexes).</p>	Without constipation Stop spasm (Intestines)	3:44
Skin without pimples -automatic	<p>The program fights pathogenic microflora, reduces inflammation, removes toxins from the body, regulates the barrier function of the skin, improves its condition.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Mode 2-automatic" program on DeVita AP and the corresponding programs on DeVita Ritm (see Clear Skin complex).</p> <p>In case of skin predisposition to acne formation, in order to prevent inflammation, it is desirable to use the program every quarter according to the above scheme.</p>	Stop acne bacteria Without toxins Skin without acne Stop pimples	5:01
Healthy nails - automatic	<p>The program causes the death of pathogenic fungi, reduces inflammation, removes toxins, helps restore the integrity of the skin of the feet and nails, and improves their condition.</p> <p><b>Application:</b> daily for a month, in conjunction with the corresponding programs on DeVita Ritm (see "Clean skin" complex). It is mandatory to maintain personal hygiene. To consolidate the result, it is desirable to repeat the cycle in a month. In order to prevent the occurrence of the disease - once or twice a year (fall-spring).</p>	No onychomycosis Without trichophyton Nails without fungus Without microsporium Without toxins	10:41
Stop skin fungus -automatic	<p>The program inhibits the development of fungal flora on the scalp and body, reduces inflammation, reduces toxic load, and promotes a better quality of life.</p> <p><b>Application:</b> daily for a month, in conjunction with the corresponding programs on DeVita Ritm (see "Clean skin" complex). It is mandatory to maintain personal hygiene. To consolidate the result, it is desirable to repeat the cycle in a month.</p>	Without Epidermophyton Without trichophyton Without microsporium Fungus-free skin Without toxins	8:23
Low back without pain -automatic	<p>The program relieves pain and improves blood flow in the lower back area.</p> <p><b>Application:</b> daily for 2-4 weeks or as needed, in conjunction with the programs "Joints without pain-automatic", "Healthy spine-automatic" on DeVita Ritm .</p>	Low back without inflammation Stop pain (low back)	6:14
Healthy back -automatic	<p>The program helps to reduce muscle spasm, reduce pain, promotes muscle fiber repair, increases muscle elasticity, reduces muscle tone, and improves the condition of all tissues of the spine.</p> <p><b>Application:</b> daily for a month, in conjunction with the programs "Joints without pain-automatic", "Healthy spine - automatic" on DeVita Ritm .</p>	Straight spine Stop spasm (muscles) Muscles without tension	8:24

# LIST OF DEVITA RITM MANUAL PROGRAMS

Nº	The name of the program	Program description
1	Active protection	<p>The human body functions in a constantly changing environment to which it must always adapt. It is negatively affected by temperature and air pressure fluctuations, electromagnetic radiation from household appliances and office technology and frequent stress. The body is also affected by physical and mental stress and ever-changing working hours. The program regulates the functioning of the nervous and endocrine systems, increases performance and alertness, eliminates fatigue and irritability, and protects against harmful radiation.</p> <p><b>Application:</b> The program is to be used daily in the morning or afternoon. It is part of the automatic programs «Healthy Day-automatic» and «Complete cleansing-automatic».</p> <p><b>Attention!</b> We do not recommend using the program in the evening or at night.</p>
2	Blood circulation regulation	<p>Lack of exercise, bad habits, an unhealthy diet and stress have a negative impact on blood circulation. The program restores local and systemic circulation and optimal oxygen uptake, improving the blood supply to all organs.</p> <p><b>Application:</b> The program can be used at any time of day (and night) if required.</p> <p><b>Please note!</b> It is not recommended to use the program within 2 months after a heart attack.</p>
3	Digestion	<p>Factors such as eating disorders or transition to new foods put increased strain on the gastrointestinal tract. This can cause indigestion, flatulence, pain and discomfort. The program has a beneficial effect on the gastrointestinal tract, promotes digestion, reduces the risk of eating disorders, reduces the feeling of fullness and discomfort, relieves acid reflux and intestinal spasms.</p> <p><b>Application:</b> The program can be used at any time of day and night, optimally between 7:00 pm and 9:00 pm. If necessary, use it again after 0.5 hours.</p>
4	Detoxification	<p>In the course of life, toxins accumulate in the body - both in the cell and in the intercellular space. Toxins disrupt the metabolism and create a favorable environment for pathogens to thrive. The program activates the function of the excretory organs, contributes to the elimination of toxins and regulates the liver and gallbladder function.</p> <p><b>Application:</b> The program can be used at any time of the day (and night). Be sure to drink 1-1.5 liters of clean drinking water daily.</p>
5	Anti-stress	<p>Almost every day we experience stress: at work, on the road or at home. Normally, the human psyche neutralizes the effects of stress, but in constant negative situations or with a sudden change in the rhythm of life, the body needs help. The program helps in difficult situations, with physical and mental stress. It normalizes the nervous system, promotes fast falling asleep and healthy sleep. For severe nervous agitation, use it together with the "Healthy Sleep" program.</p> <p><b>Application:</b> The program is to be used before going to bed. The device is placed under the pillow. If you have not fallen asleep by the end of the program, repeat the session.</p>
6	7 Chakras	<p>A human being is a complex energetic structure, which is based on seven chakras located along the spine. Each chakra is associated with certain organs, and if the chakras are damaged, the functioning of their corresponding organs is disturbed. The chakras accumulate energy. Negative energy clogs the chakras and has a detrimental effect on the body as a whole. The program restores lost energy, distributes it correctly and thus restores the functioning of organs and systems.</p> <p><b>Application:</b> The program is used in the daytime.</p>
7	Hormones normalization	<p>The endocrine system regulates organ functions through hormones and maintains the body's constant internal environment (homeostasis). The endocrine system is involved in the body's development, growth and reproductive system. It also regulates the mental state and emotional reactions. The programme normalizes the function of the system "hypothalamus - pituitary gland - adrenal gland - sex glands" and regulates the excretory functions of the body.</p> <p><b>Application:</b> The program can be used at any time of day or night, every day for at least one month. It is part of the automatic programs «Healthy Day-automatic» and «Complete cleansing-automatic».</p> <p><b>Caution!</b> Not recommended for children under 13 years of age.</p>

8	Clean Kidneys	<p>The kidneys fulfill many functions in the human body. The most important one is controlling the water and electrolyte balance, which keeps the volume and pressure of the blood and other body fluids constant. The kidneys also regulate the acid-base balance. They excrete the end products of nitrogen metabolism, i.e. urea, creatinine, uric acid, etc. Metabolic and nutritional disorders (especially excessive meat consumption), vitamin A deficiency, hypothermia and infections lead to contamination of the kidneys and the formation of kidney stones. The program promotes better elimination of toxins from the body, regulates kidney function and prevents the formation of kidney stones.</p> <p><b>Application:</b> The program can be used at any time of day and night, preferably between 5:00 and 7:00 pm. In case of toxic or infectious damage to the body, repeat the program after 0.5 hours. For normal kidney function, make sure to drink pure water - 1-1.5 liters a day.</p> <p><b>Please note:</b> When using this program, this program, increased urination may be observed.</p>
9	Clean Lymph	<p>The lymphatic system works around the clock to keep the body clean inside. It plays an important role in the removal of metabolic waste products and in maintaining immunity. Lymph is naturally filtered and purified as it passes through the lymph nodes, but sometimes this is not enough. When the function of the lymphatic system is disturbed, up to 83% of harmful substances accumulate in the interstitial space. This increases the burden on all excretory organs: liver, intestines, kidneys. The program helps dilute the lymph, removes toxins and excess water from the body, improves cell metabolism, normalizes blood pressure and increases immunity.</p> <p><b>Application:</b> The program can be used at any time. It is recommended to perform a light massage with strokes (without much pressure) during the program.</p> <p><b>Please note:</b> For the most effective elimination of toxins, drink clean water - 1-1.5 liters per day.</p>
10	Healthy Eyesight	<p>Constant strain on the eyes (reading, working at the computer), as well as changes in conditions and rhythm of life cause eye fatigue, eye irritation and reduced vision. The program relieves eye strain, fatigue and pain and contributes to the restoration of vision.</p> <p><b>Application:</b> Use before going to bed daily for 7 days. Then rest for one day. Repeat the sessions until you get the desired result. After finishing the program do not strain the eyes. For best results, combine it with ciliary-muscle training.</p>
11	Female Sexuality	<p>Sexual desire indicates different conditions of a woman. If a woman is no longer interested in her sex life, it can be assumed that she has a medical condition. In most cases, a decreased sex drive develops from a combination of several factors - somatic, mental and social. The program will help when a woman's libido is low or absent, it will ease emotional tension and mental discomfort, increasing sexual desire.</p> <p><b>Application:</b> Daily once a day.</p>
12	Male Sexuality	<p>A man's sexual desire depends on his age, physical exertion, mental state, external stimuli, etc. Stress, problems at work, fears and anxiety disorders cause psychological erectile dysfunction in men. As a result, they become anxious about the appearance and "quality" of their erection, comparing and analyzing it. And the more a man prepares for the "event", the lower his chances of success. The program is aimed at psychological improvement, regulation of the 6th chakra, it reduces anxiety and strengthens a man's self-confidence in bed.</p> <p><b>Application:</b> Daily, once a day.</p>
13	Recovery From Injuries	<p>Contusions, dislocations, and fractures can happen to both athletes and everyday people. In combination with medical treatment, the program accelerates the healing process, reduces pain, relieves swelling, stops bleeding and restores muscle elasticity.</p> <p><b>Application:</b> For injuries, daily 2-3 times a day until desired results are obtained. If necessary, combine it with the programs "Joints without inflammation" and/or "Joints without pain".</p>
14	Warm-up	<p>Gym-goers often lack energy and motivation, especially after a long day at work. The "Warm-up" program energizes you, increases the intensity of training, helps you overcome fatigue, increases metabolic processes in muscles and makes training more effective.</p> <p><b>Application:</b> Before every workout.</p>

15	Workout	<p>Successful and vigorous training often leads to fatigue and muscle soreness, which prevent you from achieving greater results. The program increases muscle strength, relieves muscle fatigue, strengthens muscles, improves blood and oxygen supply to the muscles, increases endurance, allows you to achieve maximum results with minimal stress on the body.</p> <p><b>Application:</b> The program is to be used during training.</p> <p><b>Attention!</b> It is not doping.</p>
16	Healthy Heart	<p>Heart function determines the supply of oxygen and nutrients to the organs and thus the efficient functioning of the entire body. Heart malfunctions are caused by improper lifestyle, stress, physical overload. The program regulates blood circulation and blood pressure, ensures normal blood supply to the heart muscle, promotes self-regulation of the heart and its stable work.</p> <p><b>Application:</b> Use at any time of day or night once, if necessary, twice (optimally - from 11:00 a.m. to 1:00 p.m.).</p> <p><b>Attention!</b> It is not recommended to use the program for 2 months after myocardial infarction.</p>
17	Cardiac Rhythm Normalization	<p>When the heart beats as it should, we don't notice it working. But when it gets out of rhythm, it causes us great discomfort. Cardiac arrhythmias can also occur in healthy people for various reasons - as a reaction to tight clothing, insect bites, certain medications, stress, heavy meals, stool problems. The program restores the heart rhythm, improves the blood supply to the heart and normalizes the nervous system.</p> <p><b>Application:</b> The program is used once a day for cardiac arrhythmias.</p>
18	Anti-Cold	<p>The first signs of a cold are a stuffy nose and sore throat. These symptoms are often ignored and people work despite being ill. This can lead to complications, the most common of which are upper respiratory tract infections: maxillary sinusitis, sinusitis, bronchitis, otitis media. The program activates immunity, increases resistance, reduces the likelihood of complications, relieves swelling, promotes better drainage of mucus from the sinuses and helps to recover in a short time.</p> <p><b>Application:</b> The program can be used at any time of day and night, in the cold season or during epidemics of colds it can be used daily. It provides relief from allergies.</p>
19	Free Breath	<p>Bronchospasm and cough are the body's defense reactions to an irritant or infection of the respiratory tract, and can also be a result of the effect of certain medications on the body. Bronchospasms are characteristic of people who are prone to allergies and people who suffer from bronchial asthma. The program restores the functions of the bronchopulmonary system, normalizes lung function, regulates blood supply and oxygen uptake, relieves spasms.</p> <p><b>Application:</b> Daily as required - at any time of day and night. In case of spasms and severe cough - up to 2-3 times in a row.</p>
20	Women's Health	<p>Poor environmental conditions, unhealthy habits and the sun have a negative effect on the female body, accelerating the ageing process and causing problems in the body of a woman. The program addresses a wide range of female ailments and is suitable for all ages. It relieves premenstrual syndrome and pain, and corrects menopausal disorders.</p> <p><b>Application:</b> For girls in puberty, it is recommended to use this program 1-2 times a week. For women suffering from menstrual pain and cramps: Use daily during your period. To increase the regulatory effect in case of sexual problems or disorders of the female reproductive system: Use daily together with the "Female Sexuality" program.</p> <p><b>Attention!</b> The program is not recommended during pregnancy or when trying to become pregnant.</p>
21	Men's Health	<p>Men's health does not only mean healthy potency, but also proper functioning of other internal organs. It is necessary to take care of men's health from a young age. The program is designed to harmonize the male genitourinary system. It regulates urination and the functioning of male sex glands, including the prostate gland. The older a man is, the more useful this program is for him.</p> <p><b>Application:</b> The program is to be used once a week, optimally from 9:00 to 11:00 pm. For problems with potency and/or prostate use it daily in combination with the programs "Male sexuality", "Stop prostatitis".</p>

22	Healthy Sleep	<p>Healthy and sound sleep is an important part of any person's life. During sleep, the body restores the energy used up during the day. However, any emotional upheaval can cause sleep disorders. This in turn leads to the development of physical and mental disorders. Irritability, tiredness, lack of concentration - these are the first signs of sleep disorders. The program helps reduce the excitability of the nervous system, helps you fall asleep quickly and normalize your sleep cycles.</p> <p><b>Application:</b> Daily before going to bed, the program can be repeated up to 2-3 times.</p>
23	Good Mood	<p>Depression, anxiety, fear - these concepts are interrelated. Anxiety and fear are signs of depression. Very often this condition leads to psychosomatic disorders, emotional exhaustion and dissatisfaction with oneself. Fears are strong factors preventing a person from discovering their abilities and achieving success. The program "Good mood" helps relieve depressive states and reduce anxiety, it stabilizes sleep, improves your mood and mental activity.</p> <p><b>Application:</b> Daily, 1 time in the daytime.</p>
24	Joints Without Pain	<p>Pain in the joints or spine can have several causes. Most often it is pinched or inflamed nerves, which occurs due to sedentary lifestyle, bruises, falls, overweight, herniated discs of the spine, hypothermia, vigorous physical activity, destruction of cartilage tissue or inflammation of the joint capsule. The program helps relieve pain, improve blood microcirculation, normalize the tone of muscles, ligamentous apparatus, cartilage tissue, and activates metabolic processes in damaged tissues.</p> <p><b>Application:</b> In case of an exacerbation, use the program as a supplement to the main therapy and in combination with the program "Joints without inflammation". Daily, at any time of day and night, 2 weeks, one day break, then repeat the course until you get the desired results.</p>
25	Joints Without Inflammation	<p>Every third inhabitant of our planet experiences pain or crunching in joints and spine. This problem can be caused by infections, injuries, allergic reactions, overweight, professional sports, improper diet, etc. The program improves blood supply to the damaged area, relieves inflammation, swelling and pain. It has proven itself as a good supportive and restorative remedy during the "active farm and garden work" season, also in case of sprains and severe muscle fatigue.</p> <p><b>Application:</b> daily for a month, starting 2 weeks before garden activities or sports events. Then 1-2 times a week until the end of the season. For people suffering from recurring back and joint problems, in order to maintain health, it is recommended to use this program twice a year, daily for 1 month. In case of exacerbation combine it with the "Joints without pain" program.</p>
26	Intestine Without Problems	<p>Symptoms such as bloating, nausea, intestinal discomfort, loss of appetite or abdominal pain are signs of gastrointestinal disorders. These symptoms do not appear at the initial stage, but when the disease has become chronic. To relieve pain, eliminate nausea and flatulence and normalize bowel movements, use the program "Intestine without problems".</p> <p><b>Application:</b> Daily as needed and depending on the condition 1-2 times a day.</p>
27	Digestion Normalization / Metabolism Optimization	<p>Metabolism is a series of chemical reactions that serve to maintain bodily functions. The main factors for metabolic disorders are: an unhealthy diet, taking certain medications, stress and lack of exercise. Metabolic disorders can lead to diseases such as diabetes, gout, obesity, allergies, cardiovascular disorders, as well as growth and developmental delays in children. The program aims at balancing the function of the thyroid and pancreas, the proper absorption of potassium, calcium, magnesium and sodium, and maintaining the acid-base balance.</p> <p><b>Application:</b> 1-2 times a day as required.</p>
28	Beautiful Skin	<p>Skin problems are caused by functional disorders of the internal organs. In this case, the skin takes over the body's detoxification function. Acne, pigmentation spots, rosacea, etc. appear on the skin. This affects both adolescents who are going through hormonal changes and people suffering from endocrine and digestive system disorders. The program normalizes the function of the skin glands, strengthens the skin's barrier function and improves cell metabolism.</p> <p><b>Application:</b> 1-2 times a day until the desired result is achieved.</p>

29	Beautiful Hair and Nails	<p>Hair loss, dandruff and brittle nails are signs of disorders in the body or negative environmental influences. The program helps to improve blood microcirculation, activate cell metabolism and cell regeneration. It strengthens the hair follicles and nails, stimulates hair growth and increases the skin's defense functions.</p> <p><b>Application:</b> Daily until the desired result is achieved.</p>
30	Without Allergens	<p>Allergens are substances to which the body is hypersensitive due to an abnormal immune response. They surround us everywhere: on the street, at home, at work. The most common allergens are plant pollen, dust, animal hair, medicines, household chemicals etc. The body may react to these substances with such symptoms as tears, swelling of mucous membranes, skin rashes, etc. The program alleviates the symptoms of allergies and normalizes the functions of the immune system.</p> <p><b>Application:</b> Daily as required. People who suffer from a pollen allergy should start using the program one month before the beginning of the pollen season.</p>
31	Immune Protection	<p>The immune system is a very complex, self-regulating, self-learning and multilevel system. Its coordinated functioning ensures biological defense. The program reduces the load on the liver, activates the thymus, spleen and lymphatic system, promotes rejuvenation and longevity.</p> <p><b>Application:</b> Daily or every other day, at any time of day and night. It is part of the automatic program «Health-automatic».</p>
32	Healthy Kidneys	<p>The program is used to eliminate excess fluid from the body, as well as salts that have the ability to retain fluid in the tissues. The program has a pronounced diuretic effect, it also helps restore the physiological function of the kidneys and reduce swelling.</p> <p><b>Application:</b> Daily as required for one month.</p>
33	Vegetative Regulation	<p>The autonomic nervous system controls many of the body's vital functions. Problems in its performance can be manifested in heart rhythm irregularities, blood pressure changes, excessive sweating, fatigue, weakness, brittle hair and nails, etc. The program is aimed at restoring the functions of the autonomic nervous system, stimulating the brain, normalizing metabolism, heat exchange and appetite, and promoting cell regeneration.</p> <p><b>Application:</b> Daily or every other day, at any time of day and night. It is part of the automatic program «Health-automatic».</p>
34	Cellular Respiration	<p>Hypoxia is an oxygen deprivation state, which can occur either because of a lack of oxygen entering the body or because the body tissues have stopped absorbing it. The program helps improve blood supply to the head, heart and lungs, active absorption of oxygen by the blood, relieve vascular spasms, reduce acidification and eliminate toxins.</p> <p><b>Application:</b> Daily or every other day, at any time of day and night. It is part of the automatic program «Health-automatic».</p>
35	Healthy Joints	<p>Joint pain can occur not only in older people, but also in young individuals. The pain is accompanied by inflammatory reactions, which are followed by swelling, redness of the skin and restricted movement. This problem can be caused by infections, metabolic disorders, injuries, obesity, heavy mechanical and functional loads.</p> <p>The program helps improve the condition of damaged joint tissue and relieve pain.</p> <p><b>Application:</b> Daily for one month. It is recommended to use the program in combination with the programs "Joints without pain" and "Joints without inflammation". The program is most effective when you perform light exercises on the affected joint.</p>
36	Stop Prostatitis	<p>The prostate gland is an organ that performs a number of essential functions: reproductive, secretory, barrier functions, it controls the outflow of urine and influences the psycho-emotional state. The gland is very vulnerable to various types of pathological processes. The factors that disrupt the work of the prostate may be both external (poor ecology, ionizing radiation, stress, passive lifestyle, exposure to infectious agents) and internal (genetic predisposition, comorbid chronic and acute diseases, bad habits). The program is aimed at reducing inflammatory processes and normalizing the work of the gland.</p> <p><b>Application:</b> Daily for 1-3 months. It is recommended to combine it with the program "Men's health".</p>

37	Balanced nervous system	Overexcitement of the nervous system can be caused by stress that has accumulated over a period of time. It can manifest itself as: insomnia, depression, muscle aches or headaches, irritability, mood instability etc. The program relieves nervous excitement and harmonizes the nervous system. <b>Application:</b> In case of severe mental stress, daily, at the end of the working day and/or before going to bed.
38	Recovery After Training	Exercise stimulates your metabolism, which leads to the accumulation of metabolites in the body. The strain on the autonomic nervous system increases and there may be micro-injuries to the muscles, muscle soreness, etc. Not everyone can tolerate this discomfort, which is why many people stop exercising. The program relieves muscle fatigue, promotes rapid elimination of metabolic products and supports the nervous system. <b>Application:</b> After any type of training.
39	Healthy Venous System	Problems with veins become apparent once blue, bumpy "snakes" have formed on the legs, which initially cause aesthetic discomfort, and later lead to other functional disorders. The most common cause of this problem is a sedentary lifestyle, lack of exercise or regular walking, sitting for hours at the wheel or in front of a computer screen, and constant wearing of high heels. The program helps to relieve pain and swelling, strengthen the vein wall and increase its tone, reduce inflammatory reactions, improve blood circulation in the capillaries. <b>Application:</b> As required (for heaviness, tired legs) or daily (for varicose veins). The program should be used at night or before going to bed. In the second case, place a rolled-up blanket or pillow under the legs so that the legs are in an elevated position. Do not stand up for 30 minutes after the session.
40	Without Headache	Headaches can bring a lot of problems. The causes of headaches are different: blood pressure fluctuations, intoxication, vegetative dysfunctions, osteochondrosis etc. The program is aimed at eliminating pain sensations, relieving vascular spasms, improving blood supply and nutrition of brain cells. <b>Application:</b> As needed at any time of day and night. If the pain is not gone at the end of the program, repeat the session after 20-30 minutes.
41	Normal Blood Sugar	Elevated blood sugar can occur in a person under the influence of many factors: heavy physical or mental labor, excessive consumption of carbohydrate-rich foods, severe stress, panic attacks, certain medications and functional circulatory disorders. Monitor your blood sugar levels regularly. To prevent the development of a more serious disease, use the program "Normal blood sugar", which will help normalize the function of the pancreas. <b>Application:</b> Daily for a month, and further as needed. Can be combined with the programs "Digestion Normalization (Metabolism Optimization)" and "Healthy Pancreas".
42	Functional Memory +	The human brain contains a large amount of knowledge. With age and due to great mental strain, there can be a decrease in memorization skills and the memory is not always capable of producing the necessary information. The program stimulates brain activity and the memorization process, and, combined with a healthy lifestyle, it has a positive effect on the psycho-emotional condition of a person. <b>Application:</b> It is recommended to use the program daily in the morning and afternoon, preferably in the process of studying the necessary material.
43	Nervous System - General	A tick bite can very often lead to functional disorders of the nervous system, manifested by headaches, fatigue, insomnia, depression, movement disorders of facial muscles, optic nerve dysfunction etc. The program supports and harmonizes the nervous system, ensures good coordination between its divisions, balances the functional and mental state of the person, stabilizes mental activity and mood. <b>Application:</b> For 1-6 months. See the complex «Lyme Disease».
44	Healthy Ears	Ear infections and pain can be caused by injuries or pathogens such as viruses, bacteria or fungi. The problem can occur at any age, but mostly in children and older people. The program is designed to restore the mucous membrane of various parts of the hearing organ, to reduce unpleasant sensations caused by ear diseases and to improve hearing. <b>Application:</b> Daily for 1-2 weeks. It is desirable to alternate with the program «Ears Without Inflammation» on the DeVita AP.

45	Throat Without Problems	<p>Sore throat is a problem that occurs in both adults and children at different times of the year. It can be caused by decreased immunity, infectious diseases of the upper respiratory tract, decayed teeth and gum disease. The program helps to maintain local immunity, restore the nasopharyngeal mucosa, reduce pain, swelling and hoarseness.</p> <p><b>Application:</b> Daily for 1-2 weeks. It is advisable to alternate with the DeVita AP program «Throat Without Pain».</p>
46	Strong Teeth	<p>Gum inflammation is a major problem for people, not only causing discomfort but also potentially leading to tooth loss. The most common causes are: reduced immunity, impaired blood circulation and bad habits. The program helps to optimize blood microcirculation, strengthen the gums, restore the damaged mucous membrane and thus improve the condition of the mouth and teeth.</p> <p><b>Application:</b> Daily for 1 month, see the «Healthy teeth» complex.</p>
47	Healthy Pancreas	<p>The pancreas fulfills two important functions in the body: the exocrine and endocrine functions. A malfunction of the gland can cause serious problems in the digestive tract and metabolic disorders. The program is designed to regulate pancreatic function in the case of various disorders.</p> <p><b>Application:</b> Daily for 1-6 months. It can be combined with the "Digestion Normalization (Metabolism Optimization)" and «Normal blood sugar» programs.</p>
48	Healthy Thyroid	<p>The thyroid gland is part of the human endocrine system (hormone system). It plays a very important role in the body. Its dysfunction can lead to an imbalance in the cardiovascular, respiratory, urinary, reproductive and metabolic systems. The program supports the thyroid gland under various stresses, normalizes its function and the metabolic processes in the body.</p> <p><b>Application:</b> Daily for 1-6 months. It can be combined with the "Digestion Normalization (Metabolism Optimization)" program.</p>
49	Nerve restoration	<p>Certain infectious diseases, injuries, use of medications, alcohol consumption, and working in hazardous industries can cause problems with the peripheral nervous system. To support the nervous system and restore nerve endings use the Nerve Restoration program.</p> <p><b>Application:</b> daily for one month. If necessary, the cycle should be repeated regularly 2-3 times a year.</p>
50	Recovery after infection	<p>After suffering from infectious diseases, including coronavirus infection, various problems may occur in the cardiovascular, and nervous systems. The program helps to restore the functions of the nervous system, reducing such manifestations as insomnia, anxiety, irritability.</p> <p><b>Application:</b> daily for a month after infections, in conjunction with the programs "Restoration of blood supply" on DeVita Ritm and "Deep sleep-automatic" on DeVita AP.</p>
51	Restoration of blood supply	<p>A large number of people face circulatory problems especially after suffering from infectious diseases including coronavirus infection. This can cause not only risks of developing cardiovascular problems, but also deterioration of the general condition of the body. The program helps to improve blood circulation, restore hematopoietic function, and strengthen blood vessels.</p> <p><b>Application:</b> daily for a month after infections, as regular restoration of blood supply constitutes an essential process to support own health. Can be combined with "Recovery after infection" program.</p>
52	Healthy Stomach	<p>Signs of gastric disease manifest as pain, nausea, feeling of rapid satiety, heaviness in the stomach, and belching, etc. Stomach problems can be triggered by stress, poor diet, GI problems, and taking certain medications. The program helps to relieve pain, restore mucous membranes, and normalize stomach function.</p> <p><b>Application:</b> daily for a month, in conjunction with the "Healthy stomach" program on DeVita AP and the "Stop stress" or "Digestion without problems" complex programs on DeVita Ritm. In order to prevent exacerbations, the cycle is carried out 2-3 times a year.</p>



53	Healthy spine	<p>"A healthy spine means a healthy body." This expression emphasizes the importance of the spine, as it performs many very important functions related to maintaining the health of both individual organs and the entire body. The program helps relieve pain, improves the ligamentous apparatus, increases blood supply, and restores metabolic processes in the tissues.</p> <p><b>Application:</b> daily until recovery, combining with the "Muscles without pain", "Joints Without Inflammation" program on DeVita Ritm and "Healthy back - automatic" on DeVita AP. It is advisable to apply appropriate exercises to strengthen the skeletal muscles.</p>
54	Muscles without pain	<p>Muscle pain can be caused by several factors: increased physical activity, injuries, infectious diseases, sedentary lifestyle, chronic fatigue. The program helps to relieve muscle spasm, reduce pain, improve blood supply, and increase the elasticity of muscle fibers.</p> <p><b>Application:</b> daily until recovery, in conjunction with the corresponding programs on DeVita Ritm and DeVita AP (see "Healthy back" and "Movement without pain").</p>
55	Stop asthma	<p>Asthma is a chronic inflammation of the airways that has its own symptoms. It can develop at any age, especially if there are risk factors: allergies, heredity, smoking, harmful working or living conditions. The program helps relieve allergic symptoms, bronchial spasm, and improve well-being.</p> <p><b>Application:</b> daily for a month or as needed, in conjunction with the "Without toxins" program on DeVita Ritm and "Without allergy" on DeVita AP.</p>
56	Stomach and intestines	<p>Statistically, one in two persons have problems with their stomach, or intestines, or both at the same time. Proper functioning of the GI tract is essential to a person's overall health. This is why it is essential to keep your digestive system working properly. And that's where the "Stomach and Intestines" program can help.</p> <p><b>Application:</b> daily for a month, in conjunction with the corresponding programs on DeVita Ritm and DeVita AP (see "Digestion without problems", "Gastro+" complexes).</p>
57	Normal intestines	<p>Civilized lifestyle, unfortunately, is not favorable to facilitate the work of the stomach and intestines - office work, fast food and limited mobility lead to intestinal problems in the form of pain, disruption of microflora, peristalsis and others. The program reduces spasms and pain, improves peristalsis, and stimulates recovery of mucous membranes.</p> <p><b>Application:</b> daily for 1-6 months or as needed, in conjunction with the corresponding programs on DeVita Ritm and DeVita AP (see "Digestion without problems", "Gastro+" complexes).</p>
58	Without constipation	<p>The disturbance of the process of defecation is a fairly common phenomenon and can be manifested by a feeling of incomplete emptying of the bowel, abdominal pain. Constipation causes the development of intestinal dysbiosis, chronic colitis and intoxication of the body. The program improves intestinal motility, stimulates mucosal recovery, improves blood supply to the given area, and reduces pain.</p> <p><b>Application:</b> daily for 1-6 months or as needed, in conjunction with the corresponding programs on DeVita Ritm and DeVita AP (see "Digestion without problems", "Gastro+" complexes).</p>
59	Stop pain (head)	<p>Everyone has experienced headache attacks at least once in their lives. It has a different nature and origin. Whatever the origin of the headache, the program will help you deal with it and make you feel better.</p> <p><b>Application:</b> as needed. In the presence of triggering diseases, in conjunction with the corresponding programs on DeVita Ritm and DeVita AP.</p>
60	Healthy skin	<p>Inflammation of the skin is a very unpleasant sensation manifested by burning, itching, and tissue damage. It provokes the appearance of dermatitis can cause skin infections, allergic reactions, physical or mechanical damage. The program relieves inflammatory conditions, itching, accelerates the regeneration process and helps to restore the immune barrier.</p> <p><b>Application:</b> daily for a month or as needed, in conjunction with the corresponding programs on DeVita Ritm and DeVita AP (see "Healthy skin", "Stop allergies" complexes).</p>

61	Eye rest	<p>Visual discomfort or fatigue that occurs quickly during visual work - does this sound familiar to you? This feeling is especially common when the eye is working at close distances from the object of attention. Prolonged work on computer; prolonged driving in poor visibility; professional activities requiring significant eye muscle strain; prolonged reading, especially in poor lighting; watching TV programs for hours. All of this leads to eye fatigue. The "Eye rest" program will help relieve fatigue.</p> <p><b>Application:</b> at any time of day, as needed. In addition to using the program, it is advisable to take breaks and do eye exercises</p>
62	Attention	<p>Attention determines successful orientation of a person in the surrounding world and more complete assimilation of the studied material. The "Attention" program helps you concentrate at the right time and achieve great results in various activities.</p> <p><b>Application:</b> as needed at any time of the day.</p>
63	Accurate calculation	<p>The "Accurate calculation" program is designed to help you develop your counting skills. Confident counting, like any other skill, is acquired through practice. Using the program during such training (for example, learning the multiplication table, solving mathematical problems) activates the function of specific areas of the brain, removes the syndrome of uncertainty and thus contributes to a good perception and memorization of the material.</p> <p><b>Application:</b> as needed, in the process of learning the relevant material.</p>
64	Fast Reading	<p>Today's world offers us a gigantic amount of information and knowledge from different fields that must be able to be absorbed quickly and efficiently. But it happens that psychological barriers cause not only difficulties in reading for both children and adults, but can become the reasons for the formation of emotional stress in the assimilation of material. The "Fast Reading" program, in conjunction with exercises, helps to completely eliminate the fear of reading and assimilate the material completely.</p> <p><b>Application:</b> as needed, in the process of learning the relevant material.</p>
65	Smoothing of wrinkles	<p>We perceive the appearance of wrinkles on the face as signs of aging. To eliminate this problem, different methods are used: massage, masks, creams, etc. To enhance the effect of these products and prevent the appearance of new wrinkles, use the "Smoothing of wrinkles" program. It helps to address the weakened collagen matrix, the main cause of wrinkles.</p> <p><b>Application:</b> daily for a month, then 2-3 times a week. Preferably in conjunction with the "Facial Tonus" program.</p>
66	Facial Tonus	<p>Daily facial skin care requires products with active action to improve skin complexion, deep moisturizing and toning. Regular use of the "Face Tone" program helps restore the natural processes of cell regeneration, prevents further aging of the skin, as well as eliminates fine wrinkles and smoothes deeper ones. In addition, the program gives energy to the whole body that will enable it to resist negative environmental influences.</p> <p><b>Application:</b> daily in the morning or afternoon for a month, then 2-3 times a week. Preferably in conjunction with the "Smoothing of wrinkles" program.</p>

# DEVITA RITM AUTOMATIC PROGRAMS DESCRIPTION

Name	Effect of the program	Programs	Working time
Healthy Day automatic	It contributes improving the blood supply of organs and tissues, activates the excretory systems, improves the metabolism, correctly distributes the energy balance, relieves nervous tension, strengthens the protective functions of the body. Acts against toxic loads collected in accordance with the daily biological activity of the organs.	Active protection, Blood circulation regulation, Anti-stress, Detoxification, Digestion, Hormones normalization, 7 chakras	5:21
Health automatic	It helps correcting the state of the immune system for colds viral diseases, improves cellular metabolism, which facilitates faster clearance of toxins. Also helps to restore the work of systems and organs disrupted in the course of illness.	Immune protection, Cellular respiration, Vegetative regulation, Intestine without problems	8:03
Complete cleansing automatic	Helps activating the work of excretory systems, removes toxins, increases the adaptive reserves and defenses of the body. Acts against toxic loads collected in accordance with the daily biological activity of the organs.	Active protection, Blood circulation regulation, Intestine without problems, Clean kidneys, Clean lymph, Hormones normalization	14:10
High tonus automatic	Promotes the removal of toxic loads from the liver, kidneys and lymphatic system, activates their work, improves blood circulation. Starts the natural processes of rejuvenation.	Immune protection, Cellular respiration, Clean kidneys	6:20
Healthy metabolism automatic	Promotes harmonizing the work of the body, regulates the activity of the autonomic nervous system, stabilizes the appetite, activates the work of the digestive organs and secretions. Promotes the normalization of metabolic processes and psycho-emotional state, reduction of body weight and improvement of quality of life.	Vegetative regulation, Intestine without problems, Clean kidneys	6:00
No poisoning - automatic	Helps on the reduction of the harmful effects on the body of toxins (food, drug, alcohol). Contributes to neutralize and more likely to withdraw them from the intestine and liver, reducing the burden on these organs.	Detoxification, Intestine without problems	2:26
Normal blood pressure - automatic	Relieves spasm of blood vessels, normalizes blood pressure. It helps restoring the heart and kidneys, normalizing cellular metabolism, reducing hypoxia and relieving from toxins. Improves blood circulation of organs, helps reduce swelling.	Healthy heart, Cellular respiration, Clean kidneys	6:07
Healthy sleep - automatic	It contributes removing psychological stress, relaxes the central nervous system, improves the quality and duration of sleep.	Balanced nervous system, Healthy sleep	3:06
Joints without pain - automatic	Relaxes the tensed muscles, improves the condition of the ligamentous apparatus, activates the blood supply, reduces pain and inflammation in the joints and spine. Effective with an active load on the musculoskeletal system, trauma.	Joints without pain, Joints without inflammation	3:06
Healthy joints - automatic	Relieves from inflammation and pain in the joints and spine. Improves their blood supply and nutrition. Promotes the restoration of damaged joints tissues, including articular cartilage and synovial membrane, as well as surrounding tissues: ligaments, tendons and muscles. Effective for any problems with joints: arthritis, artorosa, stiffness.	Joints without pain, Joints without inflammation, Healthy joints	5:22

Healthy skin and hair - automatic	It helps improving the blood supply of the skin, corrects the work of sebaceous and sweat glands, which increases skin immunity and reduces the manifestation of bacterial and fungal lesions. Also improves blood circulation in the capillaries of the hair follicles and the nail plates, promoting increased nutrition of hair and nails.	Beautiful skin, Beautiful hair and nails	3:36
Anti cold - automatic	It is used during colds and allergic reactions (pollinosis). Relieves inflammation and reduces the respiratory manifestation of colds and allergies: puffiness, stuffy nose, cough.	Anti-cold, Free breath	3:31
Men's Health - automatic	Promotes normalization of the prostate gland, general improvement of the urogenital system in men, increases activity and sexual endurance. It also helps to prevent the emergence of problems associated with the disruption of the functions of the urogenital system.	Men's health, Stop prostatitis	3:31
Without gout - automatic	Promotes normalizing the metabolism in the body, regulates the balance of K, Ca, Na, P. It removes inflammation in the joints, reduces swelling, promotes the removal of salts and toxins. Improves the functional state of the kidneys.	Healthy joints, Digestion normalization, Clean kidneys	6:07
Active recovery - automatic	The program promotes the recovery from infectious diseases. It normalizes nervous system function, improves blood circulation, helps restore hematopoietic function, strengthen blood vessels, reduce manifestations such as insomnia, and anxiety. <b>Application:</b> daily for a month after infections, in conjunction with the programs "Restoration of blood supply" on DeVita Ritm and "Deep sleep-automatic" on DeVita AP.	Recovery after infection Restoration of blood supply	6:22
Stomach + -automatic	The program helps to reduce pain, restore mucosa, normalize the stomach and intestines, including pH balance. <b>Application:</b> daily for a month, combined with the "Healthy Stomach" program on DeVita AP and the "Stop stress" or "Digestion without problems" programs on DeVita Ritm. In order to prevent exacerbations, the cycle is carried out 2-3 times a year.	Healthy Stomach Stomach and intestines	4:52
Without disbacteriosis - automatic	The program reduces spasms and painful symptoms, improves peristalsis, and promotes recovery of the intestinal mucosa. <b>Application:</b> daily for 1-6 months or as needed, in conjunction with the corresponding programs on DeVita Ritm and DeVita AP (see "Digestion without problems", "Gastro+" complexes).	Stomach and intestines Normal intestines	4:27
Healthy Intestine- automatic	The program improves intestinal motility, which contributes to better evacuation of the bowel contents, helps in the restoration of mucous membranes, improves blood supply to the given area, and reduces pain. <b>Application:</b> daily for 1-6 months or as needed, in conjunction with the corresponding programs on DeVita Ritm and DeVita AP (see "Digestion without problems", "Gastro+" complexes).	Normal intestines Without constipation	4:12

Healthy spine- automatic	<p>The program helps relieve pain and muscle spasm, improves the condition of the ligamentous apparatus, restores metabolic processes in tissues, increases the elasticity of muscle fibers.</p> <p><b>Application:</b> daily until recovery, in conjunction with the "Joints without inflammation" program on DeVita Ritm and "Healthy back- automatic" on DeVita AP. It is advisable to apply appropriate exercises to strengthen the skeletal muscles.</p>	Healthy spine Muscles without pain	4:07
Day (child) - automatic	<p>The program strengthens the protective functions of the body, relieves nervous tension, improves blood supply to organs and tissues, activates the work of excretory systems, improves metabolism, and proper distribution of energy balance.</p> <p><b>Application:</b> daily in the morning or afternoon.</p>	Active protection, Blood circulation regulation, Anti-stress, Detoxication, and Digestion.	3:40

# EXPLANATIONS FOR THE USE OF PROGRAMS DEVITA RITM

DeVita Ritm (manual programs)				
Name of the program	In which cases the program can be used	Effect	First results	Lasting effect
Active protection	With physical and psychological stress, decreased performance, attentiveness and tone, with chronic fatigue and irritability and daily work with household and computer equipment. To remove the bioenergetic load when working with people.	Immunostabilizing Gives energy Anti-stress	Within an hour	1 month
Blood circulation regulation	Violation of blood microcirculation (blood supply to organs), lack of oxygen in the blood, headaches, impaired blood pressure and/or vascular wall tone.	Relieves from spasms Supports tone Expands blood vessels Normalizes blood pressure	Within an hour	1 month
Digestion	Overeating, heaviness in the stomach, heartburn, bloating.	Relieves from spasms Anti-inflammatory Normalizes the acid balance	Within an hour	1 month
Detoxification	High physical activity, malnutrition, food poisoning, pain in the right hypochondrium, alcohol consumption.	Antitoxic Protects the liver Improves the outflow of bile Relieves from spasms	Within an hour	1 month
Anti-stress	Stress, high physical and psychological stress (study, work, sport)	Anti-stress Calming effect	Within an hour	1 month
7 chakras	Chronic fatigue, headaches, decreased immunity, high physical and psychological stress.	Gives tonus Harmonizing effect Anesthetic	Within an hour	1 month
Hormones normalization	In case of puberty disorders in adolescents from 12 years old, after surgical interventions, with hormonal failures, with reproductive system dysfunction in men and women.	Normalizes the hormonal background	Within 1 month	6 months
Clean kidneys	Inflammatory processes in the kidneys, bladder, swelling.	Improves urine output Decongestant Anti-inflammatory	Within 1 day	1-2 months
Clean lymph	Impaired lymphatic system, swelling, intoxication, impaired cellular metabolism, impaired blood pressure, decreased immunity, skin diseases.	Lymphatic drainage effect Decongestant Antitoxic Immunostabilizing	Within 1 day	1 month
Healthy eyesight	Decreased muscle tone of the eyes, visual impairment, visual fatigue. For prevention of visual impairment during prolonged work with computer equipment or texts on paper.	Correction effect Gives tonus Strengthening	Within 1 day	1-3 months
Female sexuality	Decreased libido, increased emotional excitability, a tendency to nervous excitement.	Antidepressant Activating	2-3 days	1 - 2 months
Male sexuality	Increased stress, excitability of the nervous system, decreased potency, fear of "failure".	Antidepressant Activating	2-3 days	1 - 2 months

Recovery from injuries	Any injuries of bones, cartilage and soft tissues, ligamentous apparatus: fractures, bruises, dislocations, cuts. After surgery.	Anesthetic Decongestant Anti-inflammatory Regenerating	Within an hour	1 month
Warm -up	Used before physical activity.	Gives tonus Gives energy	Within an hour	1 month
Workout	Used during exercise.	Gives tonus Gives energy Regulatory	Within an hour	1 month
Healthy heart	Disorder of the functional activity of the heart, impaired heart rate, impaired blood pressure, high stress loads.	Relieves from spasms Normalizes blood pressure Gives tonus Anesthetic	Within an hour	1-2 months
Cardiac rhythm normalization	Heart rhythm disturbances: tachycardia, arrhythmia, including respiratory and atrial fibrillation.	Relieves from spasms Normalizes the cardiac rhythm	Within an hour	1-2 months
Anti-cold	Diseases of the upper respiratory tract. Allergy, accompanied by swelling of the nasopharynx, lacrimation.	Anti-inflammatory Decongestant Immunostabilizing	Within an hour	Within 1 month
Free breath	Diseases of the lower respiratory system. Allergy. Spasms of the bronchopulmonary system of various origins.	Anti-inflammatory Relieves from spasms Expands blood vessels	Within an hour	1 month
Women's health	Dysfunction and inflammatory processes of the reproductive system.	Anti-inflammatory Expands blood vessels Hormone-normalizing Anti-aging	2-3 days	6-12 months
Men's health	Reproductive system dysfunctions.	Anti-inflammatory Expands blood vessels Hormone-normalizing Anti-aging	2-3 days	6-12 months
Healthy sleep	Insomnia, sleep phase disturbance, increased nervous irritability.	Antidepressant Calming Normalizes sleep	Within 1 day	1 month
Good mood	Depression, irritability, psycho-emotional stress, fatigue.	Antidepressant Regulating mood	Within 1 day	1 month
Joints without pain	Inflammation and pain of the musculoskeletal system (joints, muscles, spine), high physical exertion, including sports.	Anesthetic Regulatory	Within 1 day	1 month
Joints without inflammation	Inflammation of the musculoskeletal system, swelling, trauma, with high physical exertion, including sports.	Anti-inflammatory Decongestant Regulatory	Within 1 day	1 month
Intestine without problems	Pain, nausea, bloating, dysbiosis, decreased tone of the intestines, gall bladder, violation of the gastric mucosa and duodenum, liver function.	Relieves from spasms Anesthetic Improves the outflow of bile Anti-inflammatory Normalizes metabolism Decongestant	Within 1 day	6-12 months
Digestion normalization	Thyroid, pancreas and gastrointestinal dysfunction, swelling, impaired absorption of K, Ca, Mg, Na, impaired acid-base balance.	Regulatory Harmonizing	5-7 days	6-12 months

Beautiful skin	Decreased barrier function of the skin and cellular metabolism. Acne, dryness or excessive oily skin, relieves from psoriasis and vitiligo.	Regulates the immune and barrier functions of the skin Normalizes cell nutrition	5-7 days	1 month
Beautiful hair and nails	Hair loss and fragility, brittle nails. Dry and oily seborrhea.	Regulating the function of the skin glands Normalizes cell nutrition	5-7 days	1 month
Without allergens	Various types of allergies (respiratory - odors, dust, pets, hay fever, skin, food), swelling of the upper respiratory tract.	Antitoxic Antiallergic Decongestant Anti-inflammatory Immunostabilizing	Within 1 day	1-2 months
Immune protection	Decrease in immune defense, increased physical and mental stress.	Immunostabilizing Anti-inflammatory	Within 1 month	1 month
Healthy kidneys	Violation of metabolic processes associated with the withdrawal of certain acids and salts, impaired functional activity of the kidneys (formation of kidney stones, deposition of salts in the joints, high blood pressure).	Improves urine output Decongestant	Within 1 day	1-2 months
Cellular respiration	Headaches, problems with the vessels of the head, lethargy, drowsiness, lack of oxygen in the blood, a sedentary lifestyle. For residents of large cities.	Relieves from spasms Expands blood vessels Normalizing Fills with oxygen	Within 1 day	1-2 months
Vegetative regulation	When doing hard work for prolonged periods of time, psychological stress, imbalance of the ANS.	Regulatory	Within 1 day	1-2 months
Healthy venous system	Decreased tone of veins, stagnation in the veins of the lower extremities.	Anti-inflammatory Anesthetic	Within an hour	1-2 months
Without headache	Spasms of cerebral vessels, intense mental work, instability of blood pressure, impaired blood microcirculation in the head.	Anesthetic Relieves from spasms	Within an hour	1-2 months
Normal blood sugar	Pancreatic dysfunction associated with an increase in blood sugar.	Regulatory Normalizing	Within 1 day	6-12 months
Healthy joints	Inflammatory changes in the shell of joints and periarticular tissues, hypothermia. Physical overload of joints of a professional, domestic or sports nature. Joint injuries.	Anti-inflammatory Anesthetic	Within an hour	3-6 months
Stop prostatitis	Dysfunction of the prostate gland, decreased functional activity of the testicles, impaired urination. With irregular sexual relations.	Anti-inflammatory Regulatory	2-3 days	6-12 months
Recovery after training	Post-training period. To relieve general discomfort, headaches, lethargy, drowsiness, irritability, reduce the concentration of lactic acid, normalize respiratory processes.	Regulatory Restores energy Relieves from spasms Antitoxic	Within an hour	1 month
Functional memory +	Decrease in the speed of memorization, concentration of attention, the duration of information storage.	Stimulating	5-7 days	1 month



Balanced nervous system	Stress, irritability, emotional instability, high physical and psychological stress (study, work, sport)	Anti-stress Calming	Within an hour	1 month
Nervous System - General	Functional problems of the nervous system, a violation of the integrity of the tissue and the connection between the departments of the central nervous system, leading to a violation of psychological stability.	Anti-inflammatory Regenerating Regulatory Anti-stress	Within 1 month	3-4 months
Healthy Ears	Pain and discomfort in the ear area arising from contact with infectious agents, water, pressure drop. Also in case of hearing loss.	Regulatory Relieves from spasms Anti-inflammatory Anesthetic	Within 1 day	1 month
Throat Without Problems	Pain and swelling in the throat, nasopharynx, feeling of stuffiness, hoarseness or loss of voice arising from the inflammatory process.	Anti-inflammatory Anesthetic Decongestant	1-2 days	Within 1 month
Strong Teeth	Bleeding, friability of the gums, violation of the integrity of the mucous membranes (microtrauma, dental procedures).	Anti-inflammatory Regenerating Normalizes microcirculation	5-7 days	Within 1 month
Healthy Pancreas	Violation of the externally secretory (digestive) or within the secretory (endocrine) functions of the gland.	Normalizing Regulating Anti-inflammatory	1 month	6 - 12 months
Healthy Thyroid	Violation of metabolic processes associated with the functional imbalance of the gland.	Normalizing Regulating Anti-inflammatory	1 month	6 - 12 months
Nerve restoration	The program is used for skin sensitivity disorders and reduced limb mobility after injuries.	Anti-inflammatory, Regenerating, Regulating	1 month	3-4 months
Recovery after infection	After recovering from infectious diseases.	Normalizing, Regulating, Anti-inflammatory	1 month	3-4 months
Restoration of blood supply	After recovering from infectious diseases, injury, etc.	Anti-inflammatory, Regenerating, Normalizing microcirculation	1 month	1-2 months
Healthy Stomach	When experiencing discomfort in the stomach area: pain, heartburn, acid reflux, nausea, etc.	Pain-relieving, Anti-inflammatory, Regenerating, Normalizing	Within a day	1-2 months
Healthy spine	When experiencing back pain and fatigue associated with functional disorders, overcooling, heavy lifting, or prolonged work in an uncomfortable position.	Pain-relieving, Anti-inflammatory, Regulating	Within a day	3-4 months
Muscles without pain	When experiencing active physical exertion, injuries, infectious diseases, or prolonged periods in a static position.	Pain-relieving, Relieving spasms, Anti-inflammatory effect	Within a day	3-4 months
Stop asthma	Various types of allergies accompanied by coughing.	Anti-inflammatory, Relieving spasms, Expands blood vessels	Within an hour	1 month

Stomach and intestines	In case of improper diet, issues with the stomach and duodenum, presence of pain, heaviness, sour burps.	Pain-relieving, Anti-inflammatory, Regenerating, Normalizing	Within a day	1-2 months
Normal intestines	In case of an improper lifestyle, disruption of intestinal microflora, and peristalsis.	Relieving spasms Pain-relieving Anti-inflammatory effect Normalizing motility Anti-edematous	1 month	6-12 months
Without constipation	In case of defecation disorders, due to malnutrition, prolonged sitting, drug therapy, etc.	Pain-relieving Anti-inflammatory effect Normalizing motility	1 weeks	6-12 months
Stop pain (head)	Headaches of various origins.	Pain-relieving Relieving spasms	Within an hour	1-2 months
Healthy skin	Redness, itching, irritation, and damage to the integrity of the skin.	Anti-inflammatory, Regulating, Normalizing cellular nutrition, Regenerating	1 weeks	1 month
Eye rest	When performing tasks that require active visual strain, prolonged work on a computer, reading, and driving a car.	Tonifying, Strengthening	Within an hour	1 month
Attention	When unable to concentrate on completing a task.	Stimulating	Within an hour	1 month
Accurate calculation	When experiencing difficulties in performing arithmetic calculations.	Stimulating	Within a month	3-4 months
Fast Reading	When experiencing difficulties related to reading, memorizing poetry, etc.	Stimulating	Within a month	3-4 months
Rejuvenation	When signs of skin aging appear, and the regeneration processes decrease.	Regenerating, Normalizing cellular nutrition, Regulating	Within a month	3-6 months
Wrinkle-free skin	When fine lines start to appear on the face and to enhance the effects of facial skincare products.	Regenerating, Strengthening	1 weeks	1 month
Facial Tonus	When skin elasticity decreases and wrinkles begin to appear.	Tonifying, Anti-edematous, Strengthening	1 weeks	1 month

DeVita Ritm (automatic programs)				
Name of the program	In which cases the automatic program can be used	Effect	First results	Lasting effect
Healthy Day - automatic	Improves blood supply to organs and tissues, activates the work of excretory systems, improves metabolism, the correct distribution of energy balance, relieves nervous tension, strengthens the protective functions of the body.	Immunostabilizing Gives energy Anti-stress Antitoxic Normalizes the hormonal background	5-7 days	1 month

Health - automatic	Corrects the state of the immune system for colds and viral diseases, improves cellular metabolism, which contributes to faster cleansing of toxins. It also restores the work of systems and organs that are affected during the illness.	Anti-inflammatory Immunostabilizing Relieves from spasms Decongestant	5-7 days	1 month
Complete cleansing - automatic	Activates the work of excretory systems, removes toxins, increases adaptive reserves and the defense mechanism of the organism. The daily biological activity of organs was taken into account.	Gives energy Antitoxic Decongestant Normalizes the hormonal background	5-7 days	1 month
High tonus - automatic	Relieves the toxic load from the liver, kidneys and lymphatic system, activates their work, improves blood circulation. Starts natural rejuvenation processes.	Anti-inflammatory Improves urine output Decongestant Immunostabilizing	Within 1 month	6-12 months
Healthy metabolism - automatic	Harmonizes the work of the organism, regulates the activity of the autonomic nervous system, stabilizes the appetite, activates the digestive and excretory organs. It contributes to the normalization of metabolic processes and psycho-emotional state, weight loss and improves the quality of life.	Anti-inflammatory Normalizes metabolism Decongestant Regulates appetite	Within 1 month	6-12 months
No poisoning - automatic	Reduces the harmful effects of toxins on the organism (food, drug, alcohol). Neutralize and removes them from the intestines and liver, reducing the burden on these organs.	Antitoxic Protects the liver Improves the outflow of bile Relieves from spasms	Within 1 day	5-7 days
Normal blood pressure - automatic	Relieves from spasm of blood vessels, normalizes blood pressure. Restores the work of heart and kidneys, normalizes cellular metabolism, reduces hypoxia and removes toxins. Improves blood supply to organs, reduces swelling.	Normalizes blood pressure Relieves from spasms Improves urine output	Within 1 day	1 month
Healthy sleep - automatic	Relieves from psychological stress, relaxes the central nervous system, improves the quality and duration of sleep.	Anti-stress Calming	Within 1-2 hours	1 month
Healthy joints - automatic	Relieves from inflammation and pain in the joints and spine. Improves their blood supply and nutrition. Restores damaged joint tissues, including articular cartilage and synovial membrane, as well as surrounding tissues: ligaments, tendons and muscles. Effective for any problems with joints: arthritis, arthrosis, stiffness	Anti-inflammatory Anesthetic Decongestant Regulatory Supports the function of joints	Within 1 day	3-6 months
Joints without pain - automatic	Relaxes tense muscles, improves the condition of the ligamentous apparatus, activates blood supply, and reduces pain and inflammation in the joints and spine. Effective with active load on the musculoskeletal system, injuries.	Anti-inflammatory Decongestant Anesthetic	Within 2-3 hours	1 month

Healthy skin and hair – automatic	Improves blood circulation to the skin, corrects the functioning of the sebaceous and sweat glands, which enhances skin immunity and reduces the manifestation of bacterial and fungal infections. It also improves blood circulation in the capillaries of the hair follicles and the bed of nail plates, contributing to enhanced nutrition of hair and nails.	Regulates the immune and barrier functions of the skin, the function of the glands Normalizes cell nutrition	5-7 days	1 month
Anti cold - automatic	This program is used during colds and allergic reactions (hay fever). Relieves from inflammation and reduces the respiratory manifestation of colds and allergies: swelling, nasal congestion, cough.	Anti-inflammatory Decongestant Immunostabilizing Relieves from spasms	1-2 hours	1 month
Men's Health - automatic	Contributes to the normalization of the prostate gland, the overall improvement of the genitourinary system in men, increases activity and sexual stamina. It also prevents problems associated with the impaired function of the organs of the genitourinary system.	Anti-inflammatory Expands blood vessels Hormone-normalizing Anti-aging	2-3 days	6-12 months
Without gout - automatic	Normalizes the metabolism in the organism, regulates the balance of K, Ca, Na, P. Relieves from inflammation in the joints, reduces swelling, and removes salts and toxins. Improves the functional state of the kidneys.	Anti-inflammatory Anesthetic Decongestant Improves urine output Regulates salt metabolism	1-2 days	3-6 months
Active recovery -automatic	Helps to recover faster after having experienced infectious diseases.	Normalizing, Regulating, Anti-inflammatory	1 month	3-4 months
Stomach + -automatic	Restores the function of the stomach and duodenum, reduces inflammation, and activates the regeneration process of the mucous membrane.	Pain-relieving, Anti-inflammatory, Regenerating, Normalizing	Within a day	2-3 months
Without dysbacteriosis -automatic	Improves intestinal function, reduces pain, and spasms.	Relieving spasms Pain-relieving Anti-inflammatory effect Normalizing motility Anti-edematous	1 month	6-12 months
Healthy Intestine- automatic	Normalizes intestinal motility, promotes better evacuation of contents, and activates the process of mucosal regeneration.	Pain-relieving Anti-inflammatory effect Normalizing motility	1 weeks	6-12 months
Healthy spine- automatic	Improves the condition of the ligamentous apparatus, increases the elasticity of muscle fibers, facilitates movement, and relieves fatigue.	Pain-relieving, Relieving spasms, Anti-inflammatory effect	Within a day	3-4 months
Day (child) - automatic	Enhances the body's protective functions, relieves nervous tension, and promotes proper distribution of energy balance.	Immuno-stabilizing, Energizing Anti-stress Detoxifying	1 weeks	1 month

# EXPLANATIONS FOR THE USE OF PROGRAMS DEVITA AP

DeVita AP (manual programs)				
Name of the program	In which cases the program can be used	Effect	First results	Lasting effect
No staphylococcus	Purulent wound infections, infection, purulent inflammatory diseases of various organs (upper respiratory tract, throat, skin, musculoskeletal system, external genital organs, urinary system, meninges, foodborne toxicosis). Prevention of fetal infection during pregnancy.	Antibacterial Anti-inflammatory	2-3 hours	10-14 days
No streptococcus	Diseases of the ear, throat, nose, inflammatory processes of the musculoskeletal system, heart, kidneys, skin. Prevention of fetal infection during pregnancy.	Antibacterial Anti-inflammatory	2-3 hours	10-14 days
No chlamydia	Infections of the genitourinary system, joints, eyes, upper and lower respiratory tracts, as well as the cardiovascular system. Prevention of occupational diseases among people working with birds. As well as prevention of infection through water during bathing and transmission of the infection to the fetus during pregnancy.	Antibacterial Anti-inflammatory	1-2 days	10-14 days
No helicobacter	Inflammation of the stomach, duodenum, damage to the mucous membranes of these organs.	Antibacterial Anti-inflammatory	5-7 days	3-5 weeks
No bacteria	Problems caused by Klebsiella: inflammation of the respiratory system, urinary tract, meninges, eyes, intestinal infections, blood poisoning. Problems caused by listeria: inflammation of the digestive tract, meninges, joints, respiratory organs. In order to prevent occupational diseases of people involved in animal husbandry, poultry farming, primary processing of meat and skins. Problems caused by pneumococcus: inflammation of the respiratory tract, ears, joints, meninges, blood poisoning.	Antibacterial Anti-inflammatory	Within 1 day	14-21 days
No hepatitis A	Prevention of the virus infection when visiting countries of Asia and Africa, reduces the risk of complications.	Antiviral Anti-inflammatory	1-2 weeks	6 months
Without toxins	Mandatory program for any toxic load (viral, bacterial, protozoal, mycotic, parasitic)	Antitoxic Energy saving	1-2 hours	1 month
No herpes zoster	Shingles - skin rashes, accompanied by severe pain along the sensory nerve.	Antiviral Anti-inflammatory Anesthetic	Within 1 day	1 month
No herpes simplex	Genital herpes and herpes (cold sore) on the lips.	Antiviral Anti-inflammatory	Within 1 day	1 month

No herpes	As an addition to any of the antiherpetic programs to reduce the time of recovery and reduce inflammatory and pain manifestations.	Antiviral Anti-inflammatory Antipruritic	Within 1 day	1 month
Healthy lymph	Disorders in the lymphatic system manifested by the following symptoms: swollen fingers and a tendency to swelling, chronic fatigue, constipation, discolored skin on the legs, chest pain and swelling, dry irritated skin, inability to concentrate, cold hands and feet, cellulite.	Anti-inflammatory Decongestant Antitoxic	Within 1 day	1 month
Anti-Cold	Diseases of the nose, throat, bronchi. In conjunction with the programs: "No Staphylococcus", "No Streptococcus", "No Bacteria".	Antiviral Anti-inflammatory Decongestant	2-3 hours	5-7 days
No mould & fungi	Lesions of the skin, subcutaneous fat, mucous membranes, visual and hearing apparatus, bones, joints, as well as internal organs (bronchopulmonary system, liver, spleen, heart, blood vessels, gastrointestinal tract, kidneys, central nervous system), allergic reactions in the form of cough, nasal congestion, skin rashes.  Effective in prevention of fungal infections after visiting the sauna, pool, beach, with prolonged exposure to damp rooms.	Antifungal Anti-inflammatory	1-2 days	2 months
No candida	Diseases of the mucous membranes of the internal organs, skin, accompanied by rashes, dysbiosis, burning sensation in the genitals.	Antifungal Anti-inflammatory	2-3 days	1 month
No fungi	Damage to the skin, hair, nails, mucous membranes of the internal organs by fungi of the genus: trichophyton, microsporum, mucor. Effective in prevention of fungal infections after visiting the sauna, pool, beach, getting in contact with an infected person.	Antifungal Anti-inflammatory Antipruritic	2-3 days	3-6 months
No protozoans	Inflammatory processes of the digestive tract caused by a number of protozoa: giardia, amoeba, etc. Prevention of infection through water when swimming.	Anti-protozoan Anti-inflammatory	2-3 days	1 month
Intimate health	Diseases of the mucous membrane of the respiratory tract, genitourinary tract caused by pathogens such as: ureaplasma, mycoplasma, toxoplasma. Prevention of transmission of infection to the fetus, sexual partner.	Antibacterial Anti-protozoan Anti-inflammatory	2-3 days	1 month
No schistosomes (blood-flukes) & trematods	Damage to the liver, pancreas, lungs, circulatory and genitourinary systems, toxic-allergic reactions caused by a number of flukes: feline, hepatic, pulmonary, lanceolate, as well as schistosomes.	Antiparasitic	1-5 days	1 - 2 months

No nematodes and ascarides	Infection with roundworms, such as ascarides, pinworms, etc. Causing weight loss, peeling of the skin, rash, hair loss, bowel obstruction, inflammation of the respiratory tract, liver and pancreas. The program fights parasites at different stages of their life cycle.	Antiparasitic	2-3 days	1-2 months
No parasites, General	Suppression of the vital activity of various types of parasites. Used in combination with other helminth extermination programs.	Antiparasitic	2-3 days	1 - 2 months
No parasites, Digestive tract	Parasites of the liver and pancreas, duodenum ulcer. Used in combination with other helminth extermination programs.	Antiparasitic	2-3 days	1 -2 months
Anti-Flu	The program aims to destroy various strains of influenza viruses (A, B, C).	Antiviral	2-3 hours	5-10 days
Against flu - General	Mutated forms of flu, fever, pain.	Antiviral Antipyretic Anesthetic Anti-inflammatory	1-2 hours	5-10 days
Healthy gastrointestinal	Normalizes the condition and stops dehydration in case of dysbacteriosis caused by rotavirus infection, E. coli.	Antibacterial Antiviral Anti-inflammatory Protects cells Normalizes metabolism	2-3 hours	5-7 days
Help against poisoning	Normalizes the condition with Enterovirus infection and Salmonella infection.	Antibacterial Antiviral Anti-inflammatory Protects cells Normalizes metabolism Antitoxic	2-3 hours	10 days
Antivirus	Reduces the total viral load, cytomegalovirus infection, proceeding in the form of a persistent cold with headaches and joint pain, inflammation of the salivary glands.	Antiviral Anti-inflammatory	2-3 hours	1 month
Drainage, General	Infections of the upper respiratory tract, lymphoid tissue, liver, spleen caused by the Epstein-Barr virus.	Antiseptic Antitoxic Lymphatic drainage	Within 1 day	1 month
Drainage Antitox	With chemical toxic loads of the tissue sector with salts of heavy metals, radionuclides, food toxins.	Antitoxic Lymphatic drainage	1-2 days	14-21 days
Mexican detoxification	After surgical interventions, with chronic and drug intoxication.	Lymphatic drainage Decongestant Antitoxic Stabilizing Regulates excretory functions	6-7 days	1 month
No allergies	Various types of allergies (respiratory - odors, dust, pets, hay fever, skin, food), stabilizes the immune system.	Antitoxic Antiallergic Decongestant Anti-inflammatory Immunostabilizing	1-2 days	4-6 months

Weight Control	Presence of excess weight, impaired appetite, destroys one of the causes of obesity - Ad-36 adenovirus.	Antiviral Antitoxic	2-4 weeks	4-6 months
Against Lyme Disease 1	Tick bites. Those programs are used in pairs. If the bite happened earlier, you must take the full application scheme of the complex "Lyme Disease".	Antibacterial	1-2 days	14 – 28 days
Against Lyme Disease 2		Antiviral		
Against Cold - General	Diseases of the nose, throat, bronchi. Use together with the programs: "No staphylococcus", "No streptococcus", "No bacteria".	Decongestant Anti-inflammatory Anesthetic	1-2 days	1 month
Anti-inflammation - General	With any inflammatory processes in the body.	Decongestant Anti-inflammatory Anesthetic	1-2 days	1 – 2 months
Antiseptic	After surgical interventions, with any damage to the integrity of the skin and mucous membranes, with inflammatory processes.	Antiseptic Decongestant Anti-inflammatory Anesthetic Accelerates regeneration	2-3 hours	7 – 21 days
Anti-virus Gastrointestinal	Prevention of Coxsackie virus disease when visiting countries of Asia and Africa, reducing the risk of complications. The following symptoms can indicate the presence of the virus: dehydration, sore throat, the presence of a vesicle rash on the skin.	Antiviral Protects cells Normalizes metabolism	1-2 days	14 – 21 days
Clean Sinuses	Acute inflammation of the nasal mucosa and sinuses (sinusitis, sinusitis, etc.)	Decongestant Anti-inflammatory Anesthetic	1-2 days	14 - 21 days
Clean Sinuses - General	Together with antiviral and antibacterial programs.	Decongestant Anti-inflammatory Anesthetic Stabilizing	1-2 days	14 days
CNS Without Problems	For chronic problems with the nasal mucosa and sinuses. As one of the means of preventing an acute process in the upper respiratory tract for colds.	Decongestant Anesthetic Protects cells Revitalizing	2-3 days	6-12 months
Ears Without Inflammation	Pain and discomfort in the ear area arising from contact with infectious agents, water, pressure drop.	Decongestant Anti-inflammatory Anesthetic	1-2 days	14 - 21 days
Healthy Bile	Violation of the timely and uniform separation of bile into the intestine, which may be a consequence of irregular nutrition, consumption of harmful products, excessive consumption of fats, medicines, and the presence of parasites.	Anti-inflammatory Anesthetic Cholagogic	2-3 hours	1-4 months
Healthy Cell - 1	Prevents the development of cell mutations, reduces the negative effects of the environment and toxic load. Particular attention should be paid to those who live in adverse ecological areas, experience great emotional stress and the influence of negative production factors. Programs should be used in the "Healthy Cell" complex.	Antiviral Anti-inflammatory Protects cells	3-4 days	6-12 months
Healthy Cell - 2				



Healthy Gums	Presence of dentures, braces, crowns, as well as infections of the oral cavity that cause gum disease.	Anti-inflammatory Anesthetic Antitoxic Improves nutrition	2-3 hours	1 month
Healthy Mouth	Inflammation of the oral mucosa (sores, mucosal damage) that occurs when there are infectious agents, damage to the oral cavity, during dental interventions, an unbalanced diet.	Antiviral Anti-inflammatory Anesthetic Accelerates regeneration	1-2 days	1 month
Stop Injury (teeth)	The best effect will be acquired if used in combination with the program "Mucous Membrane Regeneration".	Anesthetic Decongestant Accelerates regeneration	2-3 hours	7-14 days
Mucous Membrane Regeneration	Inflammation of the oral mucosa (sores, mucosal damage) that occurs when there are infectious agents, damage to the oral cavity, during dental interventions, an unbalanced diet.	Anti-inflammatory Accelerates regeneration Immunostabilizing	1-2 days	1 month
Without Fistula	The best effect will be acquired if used in combination with the program "Healthy Mouth".	Anti-inflammatory Decongestant Antitoxic	1-2 days	7-14 days
Stop Adenovirus - General	After dental interventions (removal, treatment), in the presence of swelling of the gums, soreness.	Antiviral Anti-inflammatory	1-2 days	14-21 days
Healthy Nerves	Inflammation of the mucous membranes of the nose, throat, eyes, intestines.	Anti-inflammatory Protects cells	2-3 days	6-12 months
Stop Amoeba	Normalizes the condition with dysentery, gastroenteritis, intestinal ulcers, caused by the presence of protozoa. It is recommended to use to prevent infection when visiting tropical and subtropical countries of Asia, Africa, South America.	Anti-protozoan Anti-inflammatory Normalizes metabolism	2-3 days	1 month
Stop Coronavirus	Inflammation of the mucous membranes of the nose, throat, lungs.	Antiviral Anti-inflammatory	1-2 days	14-21 days
Without Klebsiella	Inflammation of the respiratory system, urinary tract, meninges, eyes, intestinal infections, blood poisoning. It is recommended as an addition for viral load, malnutrition, also for people with reduced immunity, people suffering from alcoholism, diabetes, chronic lung diseases, liver diseases.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Stop Migraine	Headaches of various origins.	Anesthetic Antispasmodic Normalizes blood circulation	Within 1 day	1 month
Stop Penicillium	Reduces toxic infections when using products affected by the fungus, in case of chronic diseases of the respiratory tract, kidneys, allergic reaction to penicillium.	Antifungal Anti-inflammatory	2-3 days	1 month

Throat Without Pain	Pain and swelling in the throat, nasopharynx, feeling of stuffiness, hoarseness or loss of voice arising from the inflammatory process.	Anti-inflammatory Decongestant Anesthetic Antitoxic	Within 1 day	14-21 days
Without Campylobacter	As an additional method to prevent infection and the development of diarrhea, for people working in livestock and poultry farms, as well as tourists visiting developing countries, people who use contaminated meat products, milk, water.	Antibacterial Anti-inflammatory	1-2 days	14-28 days
Without Clostridia	In case of damage to the integrity of the skin with contamination (cuts, bites, other injuries). When using expired or doubtful foods, canned foods, etc. Alternate with the "Without toxins" program throughout the day. Make sure to contact/see a doctor!	Antibacterial Anti-inflammatory	1-2 days	14-28 days
Without Cryptococcus	As one of the methods for preventing the disease, for people working in poultry farms engaged in breeding birds in an apartment or on their own farm. Particular attention is given to people with reduced immunity and chronic lung diseases located in places where birds gather.	Antifungal Anti-inflammatory	2-3 days	1 month
Without Enterococcus	Frequent and painful urination, pain during bowel movements, decreased potency, increased fatigue, depression, gastritis, dysbiosis, food poisoning.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Without Echinococcus	Prevention of helminthic invasion and damage to internal organs to people raising animals, livestock, collecting berries, fruits, vegetables, hunting. Especially living in South America and Europe, North Africa, New Zealand, Australia, in Russia in areas: Volga region, Western Siberia, Krasnodar Territory, in the Far East.	Antiparasitic	2-3 days	1-2 months
Without Leptospira	Jaundice, hemorrhagic fever, occupational diseases of pastoralists and agricultural workers.	Antibacterial Anti-inflammatory	1-2 days	14-28 days
Without nematodes and ascarididae - New	Infection with roundworms, such as ascariides, pinworms, etc. Causing weight loss, peeling of the skin, rash, hair loss, bowel obstruction, inflammation of the respiratory tract, liver and pancreas. An enhanced, expanded program fights parasites at different stages of their life cycle.	Antiparasitic Antitoxic	2-3 days	1-2 months
Without Kidney Papilloma	HPV damage to the mucous surface of the urinary system (renal pelvis and bladder).	Antiviral Protects cells	2-3 days	1 month
Without Papilloma - General	Damage to the skin and mucous membranes caused by the Human Papilloma Virus. It is recommended to use with the "Clean Skin" complex.	Antiviral Protects cells	2-3 days	1 month

Without Proteus	Gastritis, diarrhea, abdominal pain, vomiting, exacerbation of dermatitis, allergies, problems in the genitourinary system, long non-healing wounds.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Without Pseudomonas aeruginosa	Local purulent-inflammatory processes in the soft tissues, organs of hearing, lungs, joints, etc. Also for people who are in hospital for a long time.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Without Rickettsia	Tick bites, especially in Asia and Australia, typhus, tick-borne fevers, etc.	Antibacterial Anti-inflammatory Protects cells	1-2 days	14-28 days
Without Rotavirus	Normalizes the condition and stops dehydration in case of dysbacteriosis caused by rotavirus infection.	Antiviral Protects cells	1-2 days	14-21 days
Without Salmonella	Normalizes the condition during toxic infection, which appeared as a result of eating foods with pathogenic microorganisms.	Antibacterial Anti-inflammatory Antitoxic	1-2 days	14-28 days
Without Shigella	Normalizes the condition and stops dehydration with diarrhea caused by the presence of this bacterium.	Antibacterial Anti-inflammatory Protects cells	1-2 days	14-28 days
Without Toxoplasma	Parasitic invasion caused by intracellular protozoa and accompanied by the development of the inflammatory process in the lymph nodes, liver, brain, lungs, heart, muscles, etc.	Anti-protozoan Anti-inflammatory	2-3 days	1 month
Without Trematodes - General 1	Damage to the liver, pancreas, lungs, circulatory and genitourinary systems, toxic-allergic reactions caused by a number of flukes: feline, hepatic, pulmonary, lanceolate, etc. Programs are used together with each other.	Antiparasitic Protects cells Antitoxic	2-3 days	1-2 months
Without Trematodes - General 2				
Without Trichomonas	Inflammatory processes of the genitourinary sphere, pain during intercourse, itching and burning, discharge and unpleasant odor.	Antiprotozoal Anti-inflammatory	2-3 days	1 month
Without Ureaplasma	Inflammatory processes of the genitourinary sphere, impaired reproductive function, pain in the lower abdomen, rapid and painful urination.	Anti-protozoan Anti-inflammatory	2-3 days	1 month
Without Yersinia	Toxic-allergic reaction that occurs when hygiene measures are not followed and accompanied by skin rashes, a violation of the digestive tract, inflammation of the joints.	Antibacterial Anti-inflammatory	1-2 days	14-28 days
Stop Coronavirus - New	Inflammation of the mucous membranes of the nose, throat, lungs, weakness, decreased sense of smell and hearing. Improved and advanced program. It is advisable to use according to the methodology of the "Stop Coronavirus" complex.	Antiviral Protects cells Anti-inflammatory	1-2 days	1 month

Bacterial pneumonia	Inflammation of the respiratory system caused by Streptococcus pneumoniae, Chlamydomphila pneumoniae, Mycoplasma pneumoniae, Klebsiella pneumoniae, etc., accompanied by fever, weakness, headache, cough, shortness of breath, pain in the chest, muscles and joints.	Antibacterial Anti-inflammatory	1-2 days	1 month
Without cough	In the presence of both dry and wet cough with various infections. To get the best effect, it is recommended to use in conjunction with the "Anti-cold" and "Free Breath" programs on the DeVita Ritm device.	Anti-inflammatory Decongestant Relieves from spasms	Within 1 day	1 month
Without pneumonia	In the presence of respiratory tract infection, cough, general weakness, fatigue, drowsiness, shortness of breath, chest pain. It is recommended to use with other anti-inflammatory programs aimed at normalizing the function of the respiratory system.	Anti-inflammatory Decongestant	2-3 days	1 month
Lower back without inflammation	For lower back and sacral pain.	Anti-inflammatory Anesthetic Improves trophic	2-3 hours	1 month
Stop acne bacteria	In the presence of acne	Antibacterial Anti-inflammatory	5-7 days	1 month
Without microsporium	When ringworm (trychophytosis) appears on the scalp areas of the skin. And in order to prevent infection through contact with sick people, their personal hygiene items, clothing. Or when interacting with sick animals.	Antimycotic Anti-inflammatory	3-5 days	1-2 months
No onychomycosis	When the nail plates are affected by fungus. And in order to prevent infection when visiting public baths, saunas, swimming pools, beaches.	Antimycotic Anti-inflammatory	5-7 days	3-6 months
Zygomycetes, stop!	To reduce the likelihood of infection by fungal flora in certain chronic diseases.	Antimycotic Anti-inflammatory	2-3 days	1-2 months
Nails without fungus	When fungus affects the skin of the feet and nails. And in order to prevent infection when visiting public wet places, violation of personal hygiene.	Antimycotic Anti-inflammatory	5-7 days	3-6 months
Without trichophyton	When the fungus infects smooth skin. And in order to prevent infection by contact with infected people with their belongings or contact with sick animals.	Antimycotic Anti-inflammatory	5-7 days	1 month
Mycosis fungoides	To reduce the likelihood of infection by fungal flora in certain chronic diseases.	Antimycotic Anti-inflammatory	2-3 days	1-2 months
Solar keratosis	To prevent the spread or appearance of flaky age spots, especially in individuals who have been exposed to sun or photo radiation for long periods of time.	Anti-inflammatory Protecting cells	5-7 days	2-3 months
Stop warts	To prevent the spread of the virus in the body and reduce the formation of skin growths.	Antiviral Anti-inflammatory	5-7 days	1 month

Stop pimples	In the presence of inflammatory foci on the skin of the face and body.	Antibacterial Anti-inflammatory	2-3 days	14-28 days
Fingers without inflammation	If there is an inflammatory process in the finger area, caused by injuries sustained while performing a manicure, or injection, burr removal, ingrown nail.	Anti-inflammatory Anti-edema Accelerating regeneration	2-3 days	14-28 days
No scleroderma	To reduce the likelihood of developing the disease if there is a hereditary predisposition.	Anti-inflammatory Stabilising	5-7 days	6-12 months
Erysipelas, stop	To reduce the area of spread and prevent the development of infection when the integrity of the skin and mucous membranes is compromised.	Antibacterial Anti-inflammatory	2-3 days	1 month
Relieve psoriasis	To improve quality of life and reduce complications when the disease is present.	Anti-inflammatory Stabilising	2-3 days	1-2 months
Leprosy	In the presence of infection and to prevent contamination by contact with sick people and monkeys.	Antibacterial Anti-inflammatory	2-3 days	14-28 days
Without Epidermophyton	When a fungal infection is present on any parts of the human body and to prevent infection through contact with a sick person.	Antimycotic Anti-inflammatory	5-7 days	1 month
Restoration of blood circulation	In circulatory disorders associated with chronic or infectious diseases.	Anti-inflammatory Antitoxic Normalizing blood circulation	1-2 days	1 month
Healthy spleen	When immunity and organ problems are lowered especially when under prolonged stress.	Anti-inflammatory Immunostabilizing Antitoxic	2-3 days	1 month
Healing	For recovery of the body after illnesses or injuries.	Anti-inflammatory Antitoxic Accelerating regeneration	2-3 days	1-2 months
Gardnerella, stop!	To prevent the development of inflammatory processes in the genitourinary sphere, with hormonal and physiological changes, promiscuous sexual relations, improper use of personal hygiene products, etc.	Antibacterial Anti-inflammatory	2-3 days	14 days
Men's power	With erectile dysfunction, violation of the reproductive system in general and the prostate gland in particular.	Anti-inflammatory Improves trophic	2-3 days	6-12 months
Without actinomycetes	To prevent the spread of infection, as a result of infection through damaged skin when working with earth. Agricultural workers to use to prevent infection.	Antimycotic Anti-inflammatory	2-3 days	14 – 28 days
Stop gonococcus	A program for sexually active people to prevent infection and inflammatory process of the genitourinary sphere, as well as to improve the condition of the disease.	Antibacterial Anti-inflammatory	2-3 days	14 – 28 days
Stop parvovirus	For relief from joint pain, chronic anemia, weakened immunity.	Antiviral Anti-inflammatory	2-3 days	14 days

Without tuberculosis	As an additional method to prevent infection and inflammation, especially in persons at risk (health care workers, shelters and correctional facilities).	Antibacterial Anti-inflammatory	2-3 days	14-28 days
Without brucella	Eating insufficiently heat-treated dairy and meat products; close contact with farm animals. People keeping cows, goats, sheep are at risk	Antibacterial Anti-inflammatory	5-7 days	3-5 weeks
Skin without acne	In the presence of signs of inflammation on the skin (acne), formed as a result of changes in hormonal background (teenagers), stress, excessive ultraviolet light, use of inappropriate cosmetics.	Antibacterial Anti-inflammatory	5-7 days	1-2 months
Fungus-free skin	To prevent fungal infection in case of poor hygiene, chronic skin diseases. It is desirable to use in conjunction with the programs "Without Epidermophyton", "Without Trichophyton".	Antimycotic Anti-inflammatory Anti-itching	5-7 days	1 month
Straight spine	When the muscle corset of the back is weak; sedentary behavior; uneven physical activity, diseases of the musculoskeletal system.	Anti-inflammatory Anesthetic Improves trophic	2-3 hours	1 month
Stop pain (lower back)	In pain of various origins in the lumbar region.	Anti-inflammatory Anesthetic	2-3 hours	1 month
Without allergy new	In the presence of allergic manifestations in persons with hereditary or occupational predisposition.	Antitoxic Anti-allergic Anti-edema Anti-inflammatory Immunostabilizing	1-2 days	1-2 months
Without salmonella	When using insufficiently heat-treated fish for food, in contact with household appliances for cutting fish, non-compliance with personal hygiene.	Antiparasitic	5-7 days	1 month
Stop dementia	To prevent a decrease in mental abilities, especially in the elderly.	Anti-inflammatory Antitoxic Regulating	5-7 days	6-12 months
Stop Parkinson	As an additional means of slowing down the development of the disease in the elderly.	Anti-inflammatory Antitoxic Stabilising	5-7 days	6-12 months
Stop pain (injury)	For all types of injuries. For the best effect, the "Antiseptic" program can be added.	Anti-edema Anti-inflammatory Anesthetic Accelerating regeneration	During first 24 hours	2-3 weeks
Strong nerves	Persons suffering from mood swings, irritability, a state of prolonged fatigue.	Anti-inflammatory Accelerating regeneration Anti-stress	2-3 days	1 month
Good sleep	When the quality of sleep is disturbed, short duration of sleep.	Antitoxic Regulating	2-3 days	1 month
Calm bowel	For abdominal bloating, rumbling, flatulence, stool problems.	Anti-inflammatory Anesthetic Antitoxic	1-2 days	6-12 months

Healthy stomach	With pain in the stomach, heartburn, functional disorders of the stomach.	Anti-inflammatory Anesthetic Accelerating regeneration	1-2 days	1-4 months
Duodenitis, stop!	If you have pain in the upper abdomen, the presence of problems with the initial part of the intestine. Preferably used with the "Healthy Stomach" program.	Anti-inflammatory Anesthetic	1-2 days	1-4 months
Healthy pancreas	When the pancreas is dysfunctional.	Normalizing Anti-inflammatory	3-4 days	1-2 months
Blood sugar control	When the pancreas is malfunctioning and sugar levels are elevated.	Normalizing Anti-inflammatory	3-4 days	6-12 months
Stop diabetes	As a supplemental remedy in order to maintain good health.	Stabilising	5-7 days	6-12 months
Small intestine	In functional disorders of the intestine.	Anti-inflammatory Protecting cells	1-2 days	6-12 months
Intestines without disturbances	With functional disorders of the intestine, with transmitted infectious diseases of the intestine.	Anti-inflammatory Anesthetic Normalizing motility	1-2 days	6-12 months
Without hemorrhoids	If there is discomfort, accompanied by a feeling of heaviness and burning in the anus, bleeding at the end of defecation.	Anti-inflammatory Anti-edema Anesthetic	5-7 days	1-2 months
Stop pelvic inflammation	In inflammatory processes of the genital organs in men and women.	Anti-inflammatory Anesthetic	2-3 days	1 month
Stop chickenpox	To speed up the healing process and to prevent infection.	Antiviral Anti-edema	2-3 days	14-28 days
Stop rubella	To speed up the healing process and to prevent infection.	Antiviral Anti-edema	2-3 days	14-28 days
Without mumps	To speed up the healing process and to prevent infection.	Antiviral Anti-edema	2-3 days	14-28 days
Stop scarlet fever	To speed up the healing process and to prevent infection.	Antibacterial Anti-inflammatory	2-3 days	14-28 days
Stop measles	To speed up the healing process and to prevent infection.	Antiviral Anti-edema	2-3 days	14-28 days
No diphtheria	To speed up the healing process and to prevent infection.	Antibacterial Anti-inflammatory	2-3 days	14-28 days
Stop cholera	To speed up the recovery process and to prevent infection, persons visiting countries in Africa, South America, and India.	Antibacterial Anti-inflammatory Protecting cells	2-3 days	14-28 days
Stop whooping cough	To speed up the healing process and to prevent infection.	Antibacterial Anti-inflammatory	2-3 days	14-28 days
Without constipation	To improve the process of defecation, reduce discomfort.	Anti-inflammatory Normalizing motility	1-2 days	6-12 months
Head blood circulation	For memory impairment, presence of headaches, dizziness, decreased attention and visual acuity.	Anesthetic Normalizing blood circulation	1-2 days	1 month
For women	To improve the condition in the presence of inflammation of the reproductive organs.	Anti-inflammatory Anti-edema	2-3 days	1-2 months

Stop spasm (muscles)	When the muscle corset is weak, being in an uncomfortable position, sudden lifting of heavy weights, problems with the musculoskeletal system.	Anti-inflammatory Anesthetic Anti-edema	In 2-3 hours	1 month
No motion sickness	To reduce discomfort when driving in transport.	Antispasmodic Soothing Stabilising	20-30 min.	1 month
Stop spasm (intestine)	For intestinal disorders, impaired motility.	Antispasmodic Normalizing motility	1-2 days	6-12 months
Stop conjunctivitis	For inflammation of the mucous membrane of the eye, which is accompanied by redness, irritation or lacrimation.	Antibacterial Antiviral Anti-inflammatory Anti-edema Antitoxic	1-2 days	14-21 days
Muscles without clamps	When the muscle corset is weak, being in an uncomfortable position, sudden lifting of heavy weights, problems with the musculoskeletal system.	Antispasmodic Anti-inflammatory Antitoxic	20-30 min.	1 month
Stop fatigue	With a feeling of weakness, fatigue, lethargy, arising after prolonged stress, fatigue, insomnia.	Anti-stress Restoring energy	2-3 days	1 month

DeVita AP (automatic programs)				
Name of the program	In which cases the program can be used	Effect	First results	Lasting effect
Mode 1 - automatic	Viral and protozoan loads. Optimal disease prevention when changing climatic zones, visiting Asian countries.	Antiviral Anti-inflammatory Anti-protozoan Antitoxic	1-3 days	2-4 weeks
Mode 2 - automatic	Toxic stress caused by streptococcus and staphylococcus. Effective both for colds and for the prevention of coccal infections.	Antibacterial Anti-inflammatory Antitoxic	2-3 hours	10-14 days
Mode 3 - automatic	Bacterial burden in the digestive tract and the genitourinary system.	Antiviral Anti-inflammatory Antibacterial Antitoxic	1-3 days	2-4 weeks
Mode 4 - automatic	Presence of mycotoxins in the organism formed by mold and yeast. Reduces the risk of accumulation of carcinogenic factors.	Antifungal Anti-inflammatory Antitoxic	2-3 days	1-2 months
Cleansing - automatic	Violation of the lymphatic system, swelling, decreased immunity.	Antiseptic Antitoxic Lymphatic drainage Decongestant	Within 1 day	1 month
No parasites - automatic	Parasitic burden in the liver, pancreas, intestines, as well as toxic load caused by the vital activity of these microorganisms.	Antiparasitic Antitoxic	2-3 days	1 - 2 months



No flu - automatic	Inflammatory processes caused by various strains of the influenza virus.	Antiviral Antipyretic Anesthetic Anti-inflammatory Antitoxic	2-3 hours	5-10 days
Stop herpes S - automatic	In case of vesicles forms on the mucous membranes of the lips, genitals, eyes, mouth, skin, which are formed due to activation of the Herpes Simplex virus.	Antiviral Anti-inflammatory Antipruritic Antitoxic	Within 1 day	1 month
Stop herpes Z - automatic	Shingles - skin rashes, accompanied by severe pain along the sensory nerve.	Antiviral Anti-inflammatory Antipruritic Anesthetic Antitoxic	Within 1 day	1 month
Without cold - automatic	Inflammation of the nose, throat, bronchi.	Antiviral Anti-inflammatory Antitoxic	Within 1 day	5-7 days
Support immune system - automatic	Effective for preventing colds infections during hypothermia and during mass illnesses.	Antiviral Antitoxic	2-3 hours	5-7 days
Without helminths - automatic	Effective for the prevention of helminthiasis in children and adults.	Antiparasitic Antitoxic	2-3 days	1 - 2 months
Weight Reduction - automatic	Excess weight, resulting from viral complications. Puffiness.	Antiviral Antitoxic	1-2 days	4-6 months
Reduce Points - automatic	Violation of the lymphatic system, fluid retention, decreased immunity due to increased intoxication.	Antitoxic Lymphatic drainage Decongestant	1-2 days	4-6 months
Cell Protection 1	For removal of the general viral, herpetic and bacterial toxic loads. Reduces the activity of viruses and the negative effect of streptococcal toxins.	Antiviral Antibacterial Antitoxic	3-4 days	4-6 months
Cell Protection 2	Wrong lifestyle, great stress and physical exertion, living in adverse ecological areas. It is recommended to use the program within the application scheme of the "Healthy Cell" complex.	Antiviral Anti-protozoan Anti-inflammatory Protects cells	Within 1 week	6-12 months
Cell Protection 3	The program is effective against staphylococcal and chlamydial infections, and also reduces the negative impact of protozoan lesions of the genitourinary tract.	Antiviral Antibacterial Anti-protozoan Anti-inflammatory	3-4 days	4-6 months
Cell Protection 4	Fungal load and in order to prevent mutations and negative cell degeneration. It is recommended to use the program within the application scheme of the "Healthy Cell" complex.	Antifungal Anti-inflammatory Antitoxic Protecting cells	Within 1 week	6-12 months
Lyme 1 - automatic	Removal of the general viral, herpetic and bacterial toxic load. Reduces the activity of viruses and the negative effect of streptococcal toxins.	Antiviral Antibacterial Antitoxic	3-4 days	4-6 months

Lyme 2 - automatic	Tick bites, the detection of tick-borne borreliosis or Lyme disease. It is recommended to use the program within the application scheme of the "Lyme disease" complex.	Antiviral Antibacterial Anti-protozoan Anti-inflammatory	3-4 days	6-12 months
Lyme 3 - automatic	The program is effective against staphylococcal and chlamydial infections, and also reduces the negative impact of protozoan lesions of the genitourinary tract.	Antiviral Antibacterial Anti-protozoan Anti-inflammatory	3-4 days	4-6 months
Lyme 4 - automatic	Tick bites, the detection of tick-borne borreliosis or Lyme disease. It is recommended to use the program within the application scheme of the "Lyme disease" complex.	Antibacterial Antifungal Anti-inflammatory Antitoxic	3-4 days	6-12 months
Restoration 1 - automatic	Inflammation and damage to the mucous cavity, with dental interventions, an unbalanced diet.	Anesthetic Decongestant Anti-inflammatory Accelerates regeneration	1-2 days	1 month
Restoration 2 - automatic	After surgical interventions, with any damage to the integrity of the skin and mucous membranes, with inflammatory processes.	Anesthetic Decongestant Anti-inflammatory Accelerates regeneration	Within 1 day	2-4 weeks
Healthy Gums - Automatic	Infections and injuries of the oral cavity. It is recommended to use according to the application scheme of the "Healthy Teeth" complex	Decongestant Anti-inflammatory Accelerates regeneration	1-2 days	4-6 months
Nerve Regeneration - automatic	Negative disorder of the nervous system, a violation of the integrity of the nerve fiber caused by tick bites or other factors.	Anti-inflammatory Decongestant Anesthetic Protects cells	Within 1 week	6-12 months
Healthy Breathing - automatic	Diseases of the respiratory tract: nose, throat, bronchi, lungs.	Antiviral Antibacterial Antitoxic Anti-inflammatory	1-2 days	1 month
No motion sickness automatic	To relieve nausea attacks and improve the condition when driving in any type of transport.	Antispasmodic Soothing Stabilising	20-30 min.	1 month
SOS! Allergy automatic	To relieve discomfort and improve respiratory allergies to dust, pollen, odors, etc.	Antibacterial Antiviral Antitoxic Anti-allergic Anti-edema Anti-inflammatory	1-2 days	1-2 months
Restoration 3 automatic	After injuries or infectious diseases to accelerate regeneration processes.	Anti-inflammatory Antitoxic Accelerating regeneration Normalizing blood circulation	2-3 days	1-2 months
Deep sleep automatic	To improve the quality and duration of sleep.	Anti-inflammatory Accelerating regeneration Anti-stress Antitoxic Regulating	2-3 days	1 month

Stop brain aging automatic	For older people with decreased cognitive functions.	Anti-inflammatory Antitoxic Regulating Normalizing blood circulation	5-7 days	6-12 months
Stop fatigue automatic	To maintain a normal physiological and psychological state during strenuous work, stress.	Anti-inflammatory Accelerating regeneration Anti-stress Normalizing blood circulation Restoring energy	2-3 days	1 month
Anti-inflammatory (female) automatic	As an additional method to relieve inflammation in diseases of the reproductive system.	Anti-inflammatory Anti-edema Anesthetic Normalizing	2-3 days	1-2 months
Anti-inflammatory (male) automatic	As an additional method to relieve inflammation in diseases of the reproductive system.	Anti-inflammatory Normalizing trophics	1-2 days	1 month
Gastro active automatic	To support the stomach, gallbladder, pancreas.	Anti-inflammatory Anesthetic Accelerating regeneration Normalizing	1-2 days	6-12 months
Bowel without irritation automatic	For flatulence, bloating, intestinal disorders.	Anti-inflammatory Anesthetic Normalizing motility Accelerating regeneration	1-2 days	6-12 months
Without constipations automatic	To relieve discomfort during defecation, improving the evacuation function of the intestine.	Anti-inflammatory Anesthetic Anti-edema Normalizing motility	2-3 days	6-12 months
Skin without pimples automatic	In the presence of acne.	Antibacterial Anti-inflammatory Antitoxic	5-7 days	1 month
Healthy nails automatic	With inflammation of the skin around the nails and between the fingers, itching, changes in the color and thickness of the nail plates.	Antimycotic Anti-inflammatory Antitoxic Anti-itching	7-14 days	3-6 months
Stop skin fungus automatic	For inflammation of the scalp, in the body folds and exposed areas, itching, change and thickening of the skin and hair texture.	Antimycotic Anti-inflammatory Antitoxic Anti-itching	5-7 days	1 month
Lower back without pain automatic	For musculoskeletal problems.	Anti-inflammatory Anti-edema Anesthetic	2-3 hours	1 month
Healthy back automatic	With problems of the musculoskeletal system, weakness of the muscular corset.	Anti-inflammatory Anti-edema Anesthetic Antitoxic Improves trophic	2-3 hours	1 month

# THE BIOLOGICAL ROLE OF DELIXIR

Organ system	DeLixir	Biological Role
Cardiovascular system	Multi-Energy	Regulates the functions of blood coagulation, hematopoiesis (for women during menstruation), myocardium, relieves from swelling, prevents thinning and fragility of blood vessels. Normalizes blood pressure, improves blood composition.
	pH Balance	Strengthens blood vessels, prevents the development of varicose veins, hemorrhoids, activates the process of hematopoiesis, regulates blood pressure.
	Collagen+	Regulates the process of hematopoiesis (blood formation), supports the work of the heart muscle, lowers cholesterol, and improves blood circulation.
	Detox	Prevents blood thickening, normalizes blood pressure, cleanses the walls of blood vessels and strengthens them.
Gastrointestinal tract	Multi-Energy	Relieves from inflammation of the mucous membrane, improves the process of digestion and assimilation of nutrients, regulates the secretory function.
	pH Balance	Normalizes the acidity of gastric juice, improves digestion, prevents the development of intestinal diseases, inhibits fermentation, and stabilizes the flora.
	Collagen+	Improves the functioning of the gastrointestinal tract, normalizes appetite.
	Detox	It facilitates the digestion of proteins, helps with disorders of the intestines (atony, constipation) and indigestion. Removes toxins, suppresses pathogenic flora.
Nervous system	Multi-Energy	Improves brain function, relieves from nervous tension. Favorably affects the processes of inhibition and excitation in the cerebral cortex. Improves mood.
	pH Balance	It contributes to the formation of stress resistance, helps with irritability, sleep disturbances.
	Collagen+	It improves the state of depression, relieves from fatigue, and has a positive effect on brain activity.
	Detox	Relieves from nervous tension, improves mood.
Endocrine system, metabolism	Multi-Energy	It activates tissue respiration reactions and prevents disturbances in carbohydrate metabolism. In particular, pyruvic acid, which reduces tissue acidosis and hypoxia. Protects cells from destruction. Stimulates cell metabolism. Increases vitality. Facilitates the start of a diet for weight loss. Lowers cholesterol. Improves protein and lipid metabolism. Regulates the hormonal background.
	pH Balance	Stimulates carbohydrate and protein metabolisms, improves the absorption of calcium, iron, activates the metabolism, promotes weight loss. Lowers cholesterol, stabilizes sugar levels.
	Collagen+	Regulates protein and fat metabolism, slows down aging.
	Detox	It accelerates protein biosynthesis, cell metabolism, carbohydrate metabolism, inhibits the aging process, protects against free radicals.
Reproductive system	Multi-Energy	It normalizes hormonal levels (infertility), enhances potency, increases libido.
	pH Balance	Soothes menstrual ailments.
	Collagen+	Normalizes hormonal balance.
	Detox	Enhances sexual activity.
Urinary system	Multi-Energy	Removes water, eliminates swelling.
	pH Balance	Mild diuretic effect, prevention of urolithiasis.
	Collagen+	Reduces swelling.
	Detox	It removes excess water, eliminates swelling. Soothes kidney and bladder problems.

Respiratory system	Multi-Energy	Soothes asthma attacks.
	pH Balance	Enhances local immunity.
	Collagen+	Reduces inflammation, restores vascular elasticity.
	Detox	Strengthens the mucous membranes.
Skin / mucous membranes	Multi-Energy	It speeds up the regeneration of the skin, participates in the formation of collagen in the skin and connective tissue, improves the quality of hair and skin, strengthens the mucous membrane, eliminates bleeding gums, nosebleeds. Slows down the aging process.
	pH Balance	It activates the renewal of connective tissue, the growth of hair and nails, reduces bleeding gums, protects the mucous membranes.
	Collagen+	Accelerates skin regeneration, supports skin immunity (dermatitis).
	Detox	Revitalizes the color of hair and skin, strengthens the mucous membrane, reduces bleeding gums.
Musculoskeletal system	Multi-Energy	It improves the absorption of calcium and phosphorus, prevents the deformation of the pelvic bones in adolescent girls, a decrease in bone strength in adults, and muscle atrophy. Activates their development.
	pH Balance	Participates in bone and teeth strengthening, in renewing the connective tissue, and activates muscle work.
	Collagen+	Maintains the structure and function of bones, joints and muscles.
	Detox	Relieves from muscle spasms, strengthens the connective tissue, reduces joint pain.
Immune system	Multi-Energy	Restorative effect.
	pH Balance	Prevents colds, acidosis with infections.
	Collagen+	It has an anti-inflammatory, general strengthening effect.
	Detox	Protects against colds, restores the barrier function of the mucosa, inhibits inflammatory processes.
Organs of vision	Multi-Energy	It is necessary for the formation of visual pigment, the formation of visual acuity, and the reduction of dryness of the cornea and mucous membranes.
	pH Balance	Improves vision.
	Collagen+	Strengthens the muscles of the eye, inhibits the process of vision reduction.
	Detox	Increases visual acuity, reduces night blindness. Relieves from swelling and inflammation of the eyelids.
Energy in the body / Rejuvenation	Multi-Energy	Gives strength and energy. Protects cells from destruction (antioxidants).
	pH Balance	Reduces fatigue, promotes rejuvenation (cell renewal).
	Collagen+	Relieves from fatigue, improves performance.
	Detox	Reduces toxic load, slows down aging (antioxidants).

**RECOMMENDED DEVITA PROGRAMS  
FOR USE WITH VARIOUS FUNCTIONAL PROBLEMS**

Organ system	Disorders of the organs of the system	Programs on DeVita Ritm	Programs on DeVita AP
Cardiovascular system	Cardiac / coronary insufficiency (coronary heart disease, angina pectoris, post-infarction state, toxic damage to the heart muscle, functional overload of the heart, etc.)	Healthy Day - automatic Health - automatic Complete cleansing - automatic Healthy sleep Healthy heart	Antivirus Mode 2 – automatic No chlamydia No schistosomes (blood-flukes) & trematods Cleansing - automatic Intimate health Anti-virus Gastrointestinal Restoration 2 - automatic
	Heart rhythm disorders (arrhythmia, tachycardia)	Healthy Day - automatic Health - automatic Complete cleansing - automatic Healthy heart Cardiac rhythm normalization Healthy sleep - automatic	Antivirus Mode 2 – automatic No chlamydia No schistosomes (blood-flukes) & trematods Cleansing - automatic Anti-virus Gastrointestinal Intimate health
	Hypertension, hypotension, varicose veins, atherosclerosis and other diseases of the vascular system	Healthy Day - automatic Health - automatic Complete cleansing - automatic Normal blood pressure - automatic Healthy venous system Without headache Digestion normalization Healthy sleep - automatic	Antivirus Mode 2 – automatic Mode 3 – automatic No chlamydia, No schistosomes (blood-flukes) & trematods Cleansing - automatic Anti-virus Gastrointestinal Without Yersinia Stop Migraine Head circulation
Digestive system	Gastritis, esophagitis, heartburn, ulcer	Healthy Day - automatic Health – automatic Healthy sleep - automatic Stomach + automatic	No helicobacter Mode 4 – automatic No parasites - automatic Cleansing - automatic Restoration 2 - automatic Active gastro automatic
	Duodenitis, colitis, enteritis, ulcer, irritable bowel syndrome, dysbiosis, atony.	Healthy Day - automatic Health - automatic Complete cleansing - automatic Without disbacteriosis automatic Healthy Intestine automatic	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic No parasites – automatic Without nematodes and ascariidae - New Cleansing - automatic Without Salmonella Without Proteus Without Klebsiella Stop Amoeba Without Clostridia Active Gastro automatic Intestines without irritation automatic Without constipation automatic

Digestive system	Hepatitis, biliary dyskinesia, cholecystitis, pancreatitis, diabetes, cholelithiasis.	Healthy Day - automatic Health - automatic Complete cleansing - automatic Digestion normalization Normal blood sugar Healthy Pancreas	Mode 1 - automatic Mode 3 - automatic No parasites – automatic Without nematodes and ascaridae - New No bacteria Cleansing - automatic Healthy Bile Without Trematodes – General 1,2 Without Salmonella Without Leptospira Restoration 2 - automatic Healthy pancreas Blood sugar under control Stop diabetes Without fluke
	Acute diseases (poisoning, allergies)	Healthy Day - automatic Digestion normalization Without allergens No poisoning - automatic	Healthy Gastrointestinal Help against poisoning Drainage Antitox No allergies Cleansing - automatic Stop Penicillium Mexican detoxification Without Clostridia Without Shigella Without Proteus Without Rotavirus
Nervous system	Neurosis	Healthy Day - automatic Health – automatic Good mood Healthy sleep - automatic Functional memory +	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic Stop Migraine No parasites – automatic No bacteria No fungi Deep sleep automatic automatic
	Autonomic dysfunction	Healthy Day - automatic Health – automatic Good mood Healthy sleep - automatic Healthy heart Healthy kidneys Without headache	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No parasites – automatic No bacteria No fungi No chlamydia Anti-virus Gastrointestinal Stop Migraine
	Stroke	Healthy Day - automatic Health – automatic Recovery from injuries Healthy sleep - automatic Healthy heart Without headache Nerve restoration	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No parasites – automatic No bacteria, No fungi No chlamydia Anti-virus Gastrointestinal Stop Migraine

Nervous system	Radiculitis, neuritis, shingles	Healthy Day - automatic Health – automatic Recovery from injuries Healthy sleep – automatic Joints without pain – automatic Nervous System - General	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No parasites – automatic No bacteria, No fungi Stop herpes Z - automatic Anti-virus Gastrointestinal Restoration 2 - automatic
Endocrine system	Hypothalamic-pituitary disorders	Healthy Day - automatic Health - automatic Complete cleansing - automatic Women's health, Men's health	Mode 1 - automatic Mode 4 - automatic Cleansing - automatic Drainage Antitox Anti-virus Gastrointestinal Without Cryptococcus
	Pancreatic disease	Healthy Day - automatic Health - automatic Complete cleansing - automatic Digestion normalization Normal blood sugar Healthy Pancreas	Mode 1 - automatic No parasites – automatic Cleansing - automatic Without Trematodes - General 1,2 Without Salmonella Restoration 2 - automatic
	Thyroid disease	Healthy Day - automatic Health - automatic Complete cleansing - automatic Women's health, Men's health Digestion normalization Healthy Thyroid	Mode 1 - automatic Mode 4 - automatic Cleansing - automatic Drainage Antitox Restoration 2 - automatic
	Adrenal disease	Healthy Day - automatic Health - automatic Complete cleansing - automatic Women's health, Men's health	Mode 1 - automatic Mode 2 - automatic Cleansing - automatic Drainage Antitox Restoration 2 - automatic.
	Genital Diseases	Healthy Day - automatic Health - automatic Complete cleansing - automatic Women's health, Men's health	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Stop herpes S - automatic Anti-virus Gastrointestinal, Without Papilloma - General Without Trichomonas Without Enterococcus Restoration 2 - automatic



Reproductive system	Female Genital Disorders (cycle disorders, painful menstruation, menopause, neoplasms)	Healthy Day - automatic Health – automatic Women's health Female sexuality Complete cleansing - automatic	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Stop herpes S - automatic Without Proteus Anti-virus Gastrointestinal Without Papilloma - General Without Trichomonas Without Enterococcus Restoration 2 - automatic. Anti-inflam. (female) automatic Stop gonococcus Stop gardnerella
	Male Genital Disorders (prostatitis, adenoma, problems with potency)	Healthy Day - automatic Health – automatic Men's Health - automatic Male sexuality Complete cleansing - automatic	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Stop herpes S - automatic Without Proteus Anti-virus Gastrointestinal Without Papilloma - General Without Trichomonas Without Enterococcus Restoration 2 - automatic Anti-inflam. (male) automatic Stop gonococcus Stop gardnerella
Urinary system	Urethritis, Cystitis, Urolithiasis, gout, pyelonephritis	Healthy Day - automatic Health – automatic Complete cleansing – automatic Healthy joints - automatic Without gout - automatic	Mode 1 - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No chlamydia Stop herpes S - automatic Without Proteus Anti-virus Gastrointestinal Without Enterococcus Stop Penicillium Without Kidney Papilloma Without Trichomonas Restoration 2 - automatic Stop gardnerella

Respiratory system	Problems of the upper respiratory tract (acute respiratory infections, rhinitis, tracheitis, laryngitis, pharyngitis, sinusitis)	Healthy Day - automatic Health – automatic Without allergens Anti-cold	Without cold - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No allergies No bacteria Clean Sinuses Against Cold - General, Restoration 2 - automatic Anti-virus Gastrointestinal Healthy Breathing - automatic
	Lower respiratory problems (bronchitis, bronchial asthma, pneumonia)	Healthy Day - automatic Health – automatic Without allergens Anti-cold - automatic Stop asthma	Without cold - automatic Mode 2 - automatic Mode 3 - automatic Mode 4 - automatic Cleansing - automatic No allergies No bacteria No parasites - automatic No nematodes and ascarides Healthy Breathing - automatic Without pneumonia Without cough Restoration 2 - automatic Without Pseudomonas aeruginosa
Skin / mucous membranes	Dermatitis (contact, seborrheic, dermatomycosis)	Healthy Day - automatic Health – automatic Complete cleansing – automatic Healthy sleep – automatic Healthy skin and hair - automatic Healthy skin	Mode 1 - automatic Mode 2 - automatic Mode 4 - automatic No fungi No allergies Drainage Antitox Cleansing - automatic Restoration 2 - automatic Healthy nails automatic Stop skin fungus automatic Fingers without inflammation
	Neurodermatitis, eczema	Healthy Day - automatic Health – automatic Complete cleansing – automatic Healthy sleep – automatic Without allergens Women's health, Men's health Beautiful skin Healthy skin	Mode 2 - automatic Mode 3 - automatic No allergies Drainage Antitox Cleansing - automatic Without Proteus Without Pseudomonas aeruginosa Restoration 2 - automatic Solar keratosis Relieve psoriasis
	Allergy	Healthy Day - automatic Health - automatic Complete cleansing - automatic Without allergens Anti-cold	Mode 4 - automatic No parasites - automatic Without nematodes and ascarididae - New No allergies Drainage Antitox Cleansing - automatic SOS! Allergy Automatic

Musculoskeletal system	Osteochondrosis, arthritis, arthrosis, gout, protrusion, trauma, bursitis, myositis	<p>Healthy Day - automatic  Health – automatic  Complete cleansing – automatic  Joints without pain - automatic  Healthy joints – automatic  Without gout - automatic  Recovery from injuries  Healthy spine automatic</p>	<p>Mode 1 - automatic  Mode 2 - automatic  Mode 3 - automatic  Mode 4 - automatic  Cleansing - automatic  No chlamydia  Drainage Antitox  Against Lyme Disease 1,2  Without Pseudomonas aeruginosa  Without Enterococcus  Low back without pain automatic  Healthy back automatic  Restoration 3 automatic  Stop pain (injury)</p>
Immune system	Allergy, problems of the thymus, spleen, lymph nodes	<p>Healthy Day - automatic  Health – automatic  Complete cleansing – automatic  Without allergens  Clean lymph  Without headache</p>	<p>Without cold - automatic  Mode 2 - automatic  Mode 4 - automatic  No parasites - automatic  Without nematodes and ascarididae - New  Drainage Antitox  Cleansing - automatic  Mexican detoxification  No allergies  Against Cold - General  Throat Without Pain  Restoration 2 - automatic  Stop Penicillium  Restoration 3 automatic  Healthy spleen</p>
Sensory organs	Eye problems (optic atrophy, blepharitis, keratitis, conjunctivitis, cataracts, glaucoma, myopia, hyperopia, astigmatism)	<p>Healthy Day - automatic  Health – automatic  Complete cleansing – automatic  Healthy eyesight  Without allergens  Nervous System - General  Eye rest</p>	<p>Support immune system - automatic  Mode 2 - automatic  Mode 3 - automatic  Cleansing - automatic  Stop herpes S - automatic  Restoration 2 - automatic  No parasites - automatic  Without nematodes and ascarididae - New  No chlamydia  No fungi  No allergies  Drainage Antitox  Stop conjunctivitis</p>
	Hearing problems (otitis media, hearing loss)	<p>Healthy Day - automatic  Health – automatic  Complete cleansing – automatic  Healthy Ears  Nervous System - General  Recovery from injuries  Anti-cold</p>	<p>Without cold - automatic  Mode 2 - automatic  Mode 4 - automatic,  Intimate health  No allergies  Anti-virus Gastrointestinal  Ears Without Inflammation  Clean Sinuses  Against Cold - General  Restoration 2 - automatic</p>

Sensory organs	Problems of the sense of smell (acute respiratory disease, sinusitis, rhinitis, trauma)	Healthy Day - automatic Health – automatic Complete cleansing – automatic Nervous System - General Anti-cold Without allergens Recovery from injuries	Without cold - automatic Mode 2 - automatic Intimate health No allergies No bacteria Anti-virus Gastrointestinal Clean Sinuses - General Clean Sinuses Against Cold - General Stop Coronavirus - New Restoration 2 - automatic
Dental problems	Damage and inflammation of the mucosa, bleeding gums, pain after dental interventions	Healthy Day - automatic Health – automatic Complete cleansing – automatic Strong Teeth Recovery from injuries	Mode 1 - automatic Mode 2 - automatic Mode 4 - automatic No fungi Healthy Gums - Automatic Restoration 1-2 - automatic Stop Injury (teeth) Without Fistula Drainage Antitox

Organ system	Complex	Pendent Combination	Directions for pendant combinations
Cardiovascular system	Heart and Vessels	Healthy heart	It normalizes the heart by improving the regulation of heart rhythm and blood flow to the cardiac muscle. It prevents vascular spasms, heart and vascular disease, normalize blood pressure in the heart and kidney hypertension, and have a diuretic effect.
		Healthy kidneys	
		Healthy Liver	
Gastrointestinal	Digestion without problems	Detox S (V)	It improves the function of the liver and biliary tract, stomach, and the primary part of the intestines. It prevents the appearance of gallstones. It supports the immune and nervous systems, reduces the effects of poisoning (including drugs and heavy metals), reduces the toxic load, and normalizes the psychological state. It has an antispastic and anti-inflammatory effect .
		Healthy Liver	
		Healthy Stomach	
	Stop parasites	Healthy Liver	It reduces the toxic load, support the work of the gastrointestinal organs. It has an antispastic and anti-inflammatory effect.
		Healthy Stomach	
	Stop allergy- for food allergies	Detox S (V)	It strengthens the body's resistance to pathogenic factors, supports the immune and nervous systems, reduces the toxic load, and normalizes the psychological state, and reduce the symptoms of inflammation.
Healthy Liver			
Nervous system	Stop stress	Rehabilitation from viral infections	It helps the body to restore its energy potential, reduce the toxic load, normalize the work of organs and systems, harmonizes the psychological condition, and activates the body's defense processes, relieve stress load, normalize sleep.
	Stop insomnia	Stop stress	
Endocrine system, Metabolism	Body figure	Detox S (V)	It normalizes metabolic processes, hormonal background, contributes to the removal of toxins, harmonizes the psychological state, and reduces the toxic load.
		Healthy Liver	
		Male Health/Female Health	
Reproductive system	Women's health	Female Health	It improves blood circulation in the small pelvis, inhibits reproduction of pathogenic microflora and inflammatory processes in the genitourinary system, reduces toxic load, regulates hormonal balance, reduces pain during menstruation, and supports the immune and nervous systems.
		Healthy Liver	
		Healthy kidneys	
	Men's health	Male Health	It normalizes the urinary system, reduces the possible urinary tract infection and inflammatory reactions, reduces the toxic load, maintain full male reproductive function, and increases the quality of intimate relationships
		Healthy Liver	
		Healthy kidneys	
Urinary system	Women's health	Female Health	It improves pelvic circulation, prevents the reproduction of pathogenic microflora and inflammatory processes, regulates hormonal balance, and reduces pain during menstruation.
		Healthy Liver	
		Healthy kidneys	
	Men's health	Male Health	It normalizes the urinary system, reduces the possible urinary tract infection and inflammatory reactions, reduces the toxic load, maintain full male reproductive function, and increases the quality of intimate relationships.
		Healthy Liver	
		Healthy kidneys	

Respiratory system	Active immune system	Coronavirus Prevention	It is designed to maintain human immunity, used to prevent infection and ease the condition during the disease as an additional method of recovery, has anti-inflammatory and antiviral effect, promotes the removal of toxins.  It helps the body to restore its energy potential, reduce the toxic load, normalize the work of organs and systems, and harmonize the psychological condition, and activate the body's defense processes.
		Rehabilitation from viral infections	
		Detox S (V)	
	Stop allergy- for respiratory and asthma	Rehabilitation from viral infections	It strengthens the body's resistance to environmental factors, supports the immune and nervous systems, reduces the toxic load, normalizes the psychological state, reduces the symptoms of inflammation, eases the condition.
		Detox S (V)	
	Stop Coronavirus	Coronavirus Prevention	It is designed to maintain human immunity, used to prevent infection and ease the condition during the disease as an additional method of recovery, has anti-inflammatory and antiviral effect, promotes the removal of toxins.  It helps the body to restore its energy potential, reduce the toxic load, normalize the work of organs and systems, harmonize the psychological condition, and activate the body's defense processes.
Rehabilitation from viral infections			
Detox S (V)			
Skin/mucous	Clean skin	Clean skin	It normalizes the work of the gastrointestinal tract, helps to eliminate toxins and reduce the symptoms of inflammation, has antifungal, antibacterial and antiviral action, improve the health of the skin.
		Healthy Liver	
		Healthy Stomach	
	Stop allergy- for skin manifestations	Detox S (V)	It strengthens the body's resistance to environmental factors, supports the immune system, reduces toxic load, reduces inflammatory symptoms.
		Healthy Liver	
Musculoskeletal system	Movement without pain	Healthy joints	It has analgesic and anti-inflammatory effect, reduces pain, strengthens and regenerates connective tissue, reduces pathogenic and toxic load on the musculoskeletal system.
		Healthy kidneys	
		Healthy Liver	
	Sport performers	Healthy joints	It contributes to the strengthening and regeneration of connective tissue, reduces the toxic load, normalizes the hormonal background, reduces the symptoms of inflammation.
		Healthy Liver	
		Male Health/ Female Health	
Immune system	Active immune system	Coronavirus Prevention	It is designed to maintain human immunity, used to prevent infection and ease the condition during the disease as an additional method of recovery, has anti-inflammatory and antiviral effect, promotes the removal of toxins.  It helps the body to restore its energy potential, reduce the toxic load, normalize the work of organs and systems, harmonize the psychological condition, and activate the body's defense processes.
		Rehabilitation from viral infections	
		Detox S (V)	
	Stop Coronavirus	Coronavirus Prevention	It is designed to maintain human immunity, used to prevent infection and ease the condition during the disease as an additional method of recovery, has anti-inflammatory and antiviral effect, promotes the removal of toxins.  It helps the body to restore its energy potential, reduce the toxic load, normalize the work of organs and systems, harmonize the psychological condition, and activate the body's defense processes.
		Rehabilitation from viral infections	
		Detox S (V)	

Immune system	Stop allergy	Rehabilitation from viral infections	It strengthens the body's resistance to environmental factors, supports the immune and nervous systems, reduces the toxic load, normalizes the psychological state, reduces the symptoms of inflammation, eases the condition.
		Detox S (V)	







The combined use of all three DeVita devices is the key to a healthy and active life. Only by using all devices will your organism be completely cleansed of toxins, your organs will begin to function harmoniously, and you will feel cheerfulness and a surge of strength.

Also, the necessary smart cellular nutrition will establish vital processes.

The combined use of the technologies and smart cell nutrition supplements made by DEHolding is the key to health and longevity.

[WWW.DEHOLDING.INFO](http://WWW.DEHOLDING.INFO)